



SCIEX Triple Quad™ 3500 Instrument

System User Guide



This document is provided to customers who have purchased SCIEX equipment to use in the operation of such SCIEX equipment. This document is copyright protected and any reproduction of this document or any part of this document is strictly prohibited, except as SCIEX may authorize in writing.

Software that may be described in this document is furnished under a license agreement. It is against the law to copy, modify, or distribute the software on any medium, except as specifically allowed in the license agreement. Furthermore, the license agreement may prohibit the software from being disassembled, reverse engineered, or decompiled for any purpose. Warranties are as stated therein.

Portions of this document may make reference to other manufacturers and/or their products, which may contain parts whose names are registered as trademarks and/or function as trademarks of their respective owners. Any such use is intended only to designate those manufacturers' products as supplied by SCIEX for incorporation into its equipment and does not imply any right and/or license to use or permit others to use such manufacturers' and/or their product names as trademarks.

SCIEX warranties are limited to those express warranties provided at the time of sale or license of its products and are SCIEX's sole and exclusive representations, warranties, and obligations. SCIEX makes no other warranty of any kind whatsoever, expressed or implied, including without limitation, warranties of merchantability or fitness for a particular purpose, whether arising from a statute or otherwise in law or from a course of dealing or usage of trade, all of which are expressly disclaimed, and assumes no responsibility or contingent liability, including indirect or consequential damages, for any use by the purchaser or for any adverse circumstances arising therefrom.

For research use only. Not for use in diagnostic procedures.

AB Sciex is doing business as SCIEX.

The trademarks mentioned herein are the property of AB Sciex Pte. Ltd. or their respective owners.

AB SCIEX™ is being used under license.

© 2017 AB Sciex



AB Sciex Pte. Ltd.
Blk 33, #04-06
Marsiling Ind Estate Road 3
Woodlands Central Indus. Estate.
SINGAPORE 739256

Contents

Chapter 1 Operational Precautions and Limitations.....	8
General Safety Information.....	8
Regulatory Compliance.....	8
Australia and New Zealand.....	8
Canada.....	9
Europe.....	9
United States.....	9
International.....	10
Electrical Precautions.....	10
AC Mains Supply.....	10
Protective Earth Conductor.....	11
Chemical Precautions.....	11
System Safe Fluids.....	12
Ventilation Precautions.....	13
Environmental Precautions.....	14
Electromagnetic Environment.....	14
Decommissioning and Disposal.....	15
Qualified Personnel.....	16
Laboratory Conditions.....	16
Operating Conditions.....	16
Performance Specifications.....	16
Equipment Use and Modification.....	17
Contact Us.....	17
Technical Support.....	18
Documentation Symbols and Conventions.....	18
Related Documentation.....	19
Chapter 2 Principles of Operation.....	20
System Overview.....	20
Hardware Overview.....	21
Theory of Operation—Hardware	24
Analyst® Software Overview	26
Theory of Operation—Analyst® Software.....	31
Chapter 3 Operating Instructions — Sample Workflows.....	34
Chapter 4 Operating Instructions—Hardware.....	43
Start Up the System.....	43
Reset the Mass Spectrometer.....	45
Shut Down and Vent the System.....	45

Contents

Adjust the Integrated Syringe Pump Position.....	46
Plumb the Diverter Valve.....	50
Plumb the Diverter Valve in Injector Mode.....	51
Plumb the Diverter Valve in Diverter Mode.....	52
Chapter 5 Operating Instructions — Hardware Profiles and Projects.....	54
Hardware Profiles.....	54
Create a Hardware Profile.....	54
Add Devices to a Hardware Profile.....	58
Troubleshoot Hardware Profile Activation.....	60
Projects and Subprojects.....	61
Create Projects and Subprojects.....	61
Create Subprojects.....	62
Copy Subprojects.....	62
Switch Between Projects and Subprojects.....	63
Installed Project Folders.....	63
Back Up the API Instrument Folder.....	64
Recover the API Instrument Folder.....	64
Chapter 6 Tune and Calibrate.....	65
Verify Instrument Performance.....	65
About the Verifying or Adjusting Performance Dialog.....	66
Results Summary.....	66
Chapter 7 Automatic Optimization.....	68
About Automatic Optimization.....	69
Types of Sample Introduction.....	69
Automatically Optimize for an Analyte Using Infusion.....	70
Confirm the Presence of Compounds.....	70
Perform Automatic MS and MS/MS Optimization Using Infusion with a Known Precursor Ion and an Unknown Product Ion.....	72
Review the Optimization Results.....	75
Automatically Optimize for an Analyte Using FIA.....	76
Chapter 8 Operating Instructions — Acquisition Methods.....	84
Create an Acquisition Method Using the Acquisition Method Editor.....	84
Configure the Integrated Syringe Pump.....	85
Add an Experiment.....	85
Add a Period.....	86
Copy an Experiment into a Period.....	86
Copy an Experiment within a Period.....	86
Scan Techniques.....	87
Quadrupole-Mode Scan Types.....	87
LIT-Mode Scan Types.....	88
About Spectral Data Acquisition.....	88
Chapter 9 Operating Instructions — Analyze and Process Quantitative Data.....	90
Quantitative Analysis.....	90
Quantitation Methods.....	90

About Results Tables.....	91
Quantitation Methods and Results Tables.....	91
Create a Method Using the Quantitation Method Editor.....	91
Create a Results Table Using the Quantitation Wizard.....	93
Create a Standard Query (Optional).....	94
Results Table Right-click Menu.....	97
Peak Review and Manual Integration of Peaks.....	98
Review Peaks.....	98
Manually Integrate Peaks.....	103
Peak Review Right-Click Menu.....	104
Calibration Curves.....	105
View Calibration Curves.....	105
Overlay Calibration Curves.....	107
Calibration Curve Right-Click Menu.....	107
Sample Statistics.....	108
View the Statistics for Standards and QCs.....	108
Compare Results Between Batches.....	109
Chapter 10 Operating Instructions — Analyze and Explore Data.....	110
Open Data Files.....	110
Navigate Between Samples in a Data File.....	110
View Experimental Conditions.....	111
Show Data in Tables.....	112
Show ADC Data.....	113
Show Basic Quantitative Data.....	114
Chromatograms.....	115
Show TICs from a Spectrum.....	116
Show a Spectrum from a TIC.....	116
About Generating XICs.....	117
Generate an XIC Using a Selected Range.....	118
Generate an XIC Using the Maximum Peak.....	118
Generate an XIC Using Base Peak Masses.....	119
Extract Ion by Selecting Masses.....	119
Generate BPCs.....	120
Generate XWCs.....	122
Generate DAD Data.....	122
Generate TWCs.....	123
Adjust the Threshold.....	123
Chromatogram Panes.....	124
Spectra Panes.....	125
Graphical Data Processing.....	125
Graphs.....	126
Manage Data.....	126
Zoom In on the Y-axis.....	128
Zoom In on the X-axis.....	128
Chapter 11 Reporter Software.....	129
Analyst Reporter User Interface.....	130

Contents

Generate Reports.....	132
Chapter 12 Operating Instructions — Batches.....	134
Set Queue Options.....	134
Create and Submit a Batch.....	135
Add Sets and Samples to a Batch.....	135
Equilibrate the System.....	138
Submit a Sample or Set of Samples.....	139
Change Sample Order.....	139
Acquire Data.....	139
Set Sample Locations in the Batch Editor.....	140
Select Vial Positions Using the Locations Tab (Optional).....	141
Set Quantitation Details in the Batch Editor (Optional).....	142
Stop Sample Acquisition.....	142
Import Batch Files.....	143
Build a Batch as a Text File.....	143
Import a Batch from a Text File.....	143
Batch Editor Right-click Menu.....	144
Queue States and Device Status.....	145
Queue States.....	145
View Instrument and Device Status Icons.....	146
Queue Right-click Menu.....	147
Chapter 13 Service and Maintenance Information.....	149
Recommended Maintenance Schedule.....	149
Clean the Surfaces.....	151
Clean the Front-End.....	151
Symptoms of Contamination.....	152
Required Materials	152
Cleaning Best Practices.....	153
Prepare the Mass Spectrometer.....	154
Clean the Curtain Plate.....	155
Clean the Front of the Orifice Plate.....	156
Put the Mass Spectrometer Back in Service.....	156
Inspect the Roughing Pump Oil Level.....	157
Empty the Source Exhaust Drain Bottle.....	157
Storage and Handling.....	159
Chapter 14 Mass Spectrometer Troubleshooting.....	160
Appendix A Operating Instructions — Manual Compound Optimization.....	163
About Manual Compound Optimization.....	164
About Scan Types.....	164
Manually Optimize an Analyte.....	164
Confirm the Presence of Compounds.....	165
Optimize MS-Specific Parameters.....	166
Determine the Product Ions for Optimization.....	168
Optimize Collision Cell Exit Potential for each Product Ion.....	169
Manually Optimize the Ion Source and Gas Parameters.....	170

Prepare the Ion Source.....	170
Optimize the Ion Source Parameters.....	171
Advanced Parameters.....	172
Optimize AF2.....	172
About Collision Energy Spread (CES)	172
Appendix B Parameters for SCIEX Triple Quad 3500 Instruments.....	173
Appendix C Calibration Ions and Solutions.....	175
Appendix D Toolbar Icons.....	177
Appendix E Glossary of Symbols.....	186
Appendix F Glossary of Warnings.....	191

Operational Precautions and Limitations

1

Note: Before operating the system, carefully read all of the sections of this guide.

This section contains general safety-related information and provides regulatory compliance information. It also describes potential hazards and associated warnings for the system and the precautions that should be taken to minimize the hazards.

In addition to this section, refer to [Glossary of Symbols on page 186](#) for information about the symbols and conventions used in the laboratory environment, on the system, and in this documentation. Refer to the *Site Planning Guide* for site requirements, including AC mains supply, source exhaust, ventilation, compressed air, nitrogen, and roughing pump requirements.

General Safety Information

To prevent personal injury or system damage, read, understand, and obey all of the safety precautions and warnings in this document, the manufacturer chemical safety data sheet (SDS), and product label information. These labels are shown with internationally recognized symbols. Failure to heed these warnings could result in serious injury.

This safety information is intended to supplement federal, state, provincial, and local environmental health and safety (EHS) regulations. The information provided covers system-related safety information applicable to the operation of the system. It does not cover every safety procedure that should be practised. Ultimately, the user and the organization are responsible for compliance with federal, state, provincial, and local EHS regulations and for maintaining a safe laboratory environment.

Refer to the appropriate laboratory reference material and standard operating procedures.

Regulatory Compliance

This system complies with the regulations and standards listed in this section. Refer to the Declaration of Conformity included with the system and the individual system components for dated references. Applicable labels have been affixed to the system.

Australia and New Zealand

- **Electromagnetic Compatibility (EMC):** Radio Communications Act 1992 as implemented in these standards:
 - Electromagnetic Interference—AS/NZS CISPR 11/ EN 55011/ CISPR 11 (Class A). Refer to [Electromagnetic Interference on page 15](#).

- **Safety:** AS/NZ 61010-1 and IEC 61010-2-061

Canada

- **Electromagnetic Interference (EMI):** CAN/CSA CISPR11. This ISM device complies with Canadian ICES-001. Refer to [Electromagnetic Interference on page 15](#).
- **Safety:**
 - CAN/CSA C22.2 No. 61010-1
 - CAN/CSA C22.2 No 61010-2-061

Europe

- **Electromagnetic Compatibility (EMC):** Electromagnetic Compatibility directive 2014/30/EU as implemented in these standards:
 - EN 61326-1
 - EN 55011 (Class A)Refer to [Electromagnetic Compatibility on page 14](#).
- **Safety:** Low Voltage Directives 2014/35/EU as implemented in these standards:
 - EN 61010-1
 - EN 61010-2-061
- **Waste Electrical and Electronic Equipment (WEEE):** Waste Electrical and Electronic Equipment 2012/96/EEC, as implemented in EN 40519. Refer to [Waste Electrical and Electronic Equipment on page 16](#).
- **Packaging and Packaging Waste (PPW):** Packaging and Packaging Waste Directive 94/62/EC
- **RoHS Restriction of Hazardous Substances:** RoHS Directive 2011/65/EU

United States

- **Radio Emissions Interference Regulations:** 47 CFR 15, as implemented in FCC Part 15 (Class A)
- **Safety:** Occupational Safety and Health Regulations, 29 CFR 1910, as implemented in these standards:
 - UL 61010-1
 - IEC 61010-2-061

International

- **Electromagnetic Compatibility (EMC):**

- IEC 61326-1
- IEC CISPR 11 (Class A)
- IEC 61000-3-2
- IEC 61000-3-3

Refer to [Electromagnetic Compatibility on page 14](#).

- **Safety:**

- IEC 61010-1
- IEC 61010-2-061

Electrical Precautions



WARNING! Electrical Shock Hazard. Do not remove the covers. Removing the covers might cause injury or malfunctioning of the system. The covers need not be removed for routine maintenance, inspection, or adjustment. Contact a SCIEX Field Service Employee (FSE) for repairs that require the covers to be removed.

- Follow required electrical safe work practices.
- Use cable management practices to control electrical cables. This will reduce the chance of a tripping hazard.

For information about system electrical specifications, refer to the *Site Planning Guide*.

AC Mains Supply

Connect the system to a compatible AC mains supply as instructed in this guide.



WARNING! Electrical Shock Hazard. Use only qualified personnel for the installation of all of the electrical supplies and fixtures, and make sure that all of the installations adhere to local regulations and safety standards.



WARNING! Electrical Shock Hazard. Make sure that the system can be disconnected from the mains supply outlet in an emergency. Do not block the mains supply outlet.



WARNING! Electrical Shock Hazard. Use only the power cables supplied with the system. Do not use power cables that are not properly rated for the operation of this system.

An external line transformer is not needed for the mass spectrometer, optional bench, or roughing pump.

Protective Earth Conductor

The mains supply must include a correctly installed protective earth conductor. The protective earth conductor must be installed or checked by a qualified electrician before the system is connected.



WARNING! Electrical Shock Hazard. Do not intentionally interrupt the protective earth conductor. Any interruption of the protective earth conductor creates an electrical shock hazard.



WARNING! Electrical Shock Hazard. Make sure that a protective earth conductor (grounding cable) is connected between the sample loop and an appropriate grounding point at the ion source. This supplementary grounding will reinforce the safety configuration specified by SCIEX.

Chemical Precautions



WARNING! Radiation Hazard, Biohazard, or Toxic Chemical Hazard. Determine whether decontamination is required prior to cleaning or maintenance. The customer must decontaminate the system prior to cleaning or maintenance if radioactive materials, biological agents, or toxic chemicals have been used with the system.



WARNING! Environmental Hazard. Do not dispose of system components in municipal waste. Follow local regulations when disposing of components.



WARNING! Biohazard, Toxic Chemical Hazard. Connect the drain tubing to the mass spectrometer and the source exhaust drain bottle properly, to prevent leaks.

Operational Precautions and Limitations

- Determine which chemicals have been used in the system prior to service and regular maintenance. Refer to the *Safety Data Sheets* for the health and safety precautions that must be followed with chemicals. SCIEX Safety Data Sheets can be found at sciex.com/tech-regulatory.
- Work in a well-ventilated area or fume hood.
- Always wear assigned personal protective equipment, including powder-free neoprene or nitrile gloves, safety glasses, and a laboratory coat.
- Avoid ignition sources when working with flammable materials, such as isopropanol, methanol, and other flammable solvents.
- Take care in the use and disposal of any chemicals. Potential risk of personal injury if proper procedures for handling and disposing of chemicals are not followed.
- Avoid skin contact with chemicals during cleaning and wash hands after use.
- Make sure that all exhaust hoses are connected properly and that all connections are functioning as designed.
- Collect all spent liquids and dispose of them as hazardous waste.
- Comply with all of the local regulations for the storage, handling, and disposal of biohazardous, toxic, or radioactive materials.
- (Recommended) Use secondary containment trays beneath the roughing pump, the solvent bottles, and the waste collection container to capture potential chemical spills.

System Safe Fluids

The following fluids can safely be used with the system. Refer to [Required Materials on page 152](#) for information about safe cleaning solutions.

CAUTION: Potential System Damage. Do not use any other fluid until confirmation is received from SCIEX that it does not present a hazard. This is not an exhaustive list.

- **Organic Solvents**
 - MS-grade acetonitrile; up to 100%
 - MS-grade methanol; up to 100%
 - Isopropanol; up to 100%
 - HPLC-grade or higher water; up to 100%
- **Buffers**
 - Ammonium acetate; less than 1%
 - Ammonium formate; less than 1%

- **Acids and Bases**

- Formic acid; less than 1%
- Acetic acid; less than 1%
- Trifluoroacetic acid (TFA); less than 1%
- Heptafluorobutyric acid (HFBA); less than 1%
- Ammonia/ammonium hydroxide; less than 1%

Ventilation Precautions

The venting of fumes and disposal of waste must comply with all of the federal, state, provincial, and local health and safety regulations. It is the responsibility of the customer to make sure that the air quality is maintained in compliance with local health and safety regulations.

The source exhaust system and roughing pump must be vented to a dedicated laboratory fume hood or an external exhaust system.



WARNING! Fire Hazard. Make sure that the source exhaust system is connected and functioning to prevent flammable vapor from accumulating in the ion source.



WARNING! Radiation Hazard, Biohazard, or Toxic Chemical Hazard. Take care to vent exhaust gases to a dedicated laboratory fume hood or exhaust system and make sure that the ventilation tubing is secured with clamps. Make sure that the laboratory has appropriate air exchange for the work performed.



WARNING! Radiation Hazard, Biohazard, or Toxic Chemical Hazard. Do not operate the mass spectrometer if the source exhaust drain and roughing pump exhaust hoses are not properly connected to the laboratory ventilation system. Perform a regular check of the exhaust tubing to make sure that there are no leaks. The use of mass spectrometers without proper system ventilation might constitute a health hazard and might result in serious injury.



WARNING! Radiation Hazard, Biohazard, or Toxic Chemical Hazard. Use the ion source only if you have knowledge of and training in the proper use, containment, and evacuation of toxic or injurious materials used with the ion source.

Operational Precautions and Limitations



WARNING! Puncture Hazard, Radiation Hazard, Biohazard, or Toxic Chemical Hazard. Discontinue use of the ion source if the ion source window is cracked or broken and then contact a SCIEX Field Service Employee (FSE). Any toxic or injurious materials introduced into the equipment will be present in the source exhaust output. Dispose of sharps following established laboratory safety procedures.

Environmental Precautions

Use qualified personnel for the installation of electrical mains, heating, ventilation, and plumbing supplies and fixtures. Make sure that all of the installations comply with local bylaws and biohazard regulations. For information about the required environmental conditions for the system, refer to the *Site Planning Guide*.

Allow access space around the equipment when setting up the system.



DANGER! Explosion Hazard. Do not operate the system in an environment containing explosive gases. The system is not designed for operation in an explosive environment.



WARNING! Biohazard. For biohazardous material use, always comply with local regulations for hazard assessment, control, and handling. This system or any part is not intended to act as a biological containment.

CAUTION: Potential Mass Shift. Maintain a stable ambient temperature. If the temperature changes by more than 2 °C per hour, then the resolution and mass calibration might be affected.

Electromagnetic Environment

Electromagnetic Compatibility

Basic Electromagnetic Environment: Environment existing at locations characterized by being supplied directly at low voltage from the public mains network.

Performance Criteria A (Criteria A): Equipment shall operate as intended with no degradation of performance and no loss of function during or after test.

Performance Criteria B (Criteria B): Equipment may experience loss of function (one or more) during test but shall operate as intended with some degradation of performance and functions self-recoverable after test.

Performance Criteria C (Criteria C): Equipment may experience loss of function (one or more) during test but shall operate as intended with some degradation of performance and functions recoverable by operator after test.

The equipment is intended for use in a basic electromagnetic environment.

The expected performance loss under the electromagnetic immunity conditions is less than 20% change in total ion count (TIC).

Make sure that a compatible electromagnetic environment for the equipment can be maintained so that the device will perform as intended. If the power supply line is subject to high electrical noise, then install a surge protector.

Electromagnetic Interference

Class A Equipment: Equipment which is suitable for use in all establishments other than domestic and those directly connected to a low voltage power supply network which supplies buildings used for domestic purposes. [Derived from CISPR 11:2009, 5.3] Class A equipment shall meet Class A limits.

This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to Part 15 of the FCC (Federal Communications Commission) Compliance Rules.

These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the operator's manual, can cause harmful interference to radio communications.

Operation of this equipment in a residential area is likely to cause harmful interference in which case you will be required to correct the interference, at your own expense. Changes or modifications not expressly approved by the manufacturer could void your authority to operate the equipment.

Decommissioning and Disposal



WARNING! Environmental Hazard. Do not dispose of system components in municipal waste. Follow local regulations when disposing of components.

Before decommissioning, decontaminate the entire system following local regulations.

When removing the system from service, separate and recycle different materials according to national and local environmental regulations. Refer to [Storage and Handling on page 159](#).

Note: SCIEX will not accept any system returns without a completed Decontamination Form. Contact an FSE to obtain a copy of the form.

Do not dispose of system components or subassemblies, including computer parts, as unsorted municipal waste.

Waste Electrical and Electronic Equipment

Follow local municipal waste ordinances for proper disposal provisions to reduce the environmental impact of waste, electrical, and electronic equipment (WEEE). To safely dispose of this equipment, contact a local Customer Service office for complimentary equipment pick-up and recycling.

Qualified Personnel

Only qualified SCIEX personnel shall install, inspect, and service the equipment. After installing the system, the Field Service Employee (FSE) uses the *Customer Familiarization Checklist* to orient the customer on system operation, cleaning, and basic maintenance.

Only personnel qualified by the manufacturer shall maintain the equipment. A laboratory designate can be familiarized with the Qualified Maintenance Person (QMP) procedures during the installation. A QMP is a person who is suitably aware of the electrical and chemical risks associated with servicing laboratory equipment.

Laboratory Conditions

Operating Conditions

The system is designed to operate safely under these conditions:

- Indoors
- Altitude: Up to 2 000 m (6 400 feet) above sea level
- Ambient temperature: 5 °C (41 °F) to 40 °C (104 °F)
- Relative humidity: 80% for temperatures up to 31 °C (88 °F), decreasing linearly to 50% at 40 °C (104 °F)
- Mains supply voltage fluctuations: $\pm 10\%$ of the nominal voltage
- Transient overvoltages: Up to the levels of Overvoltage Category II
- Temporary overvoltages on the mains supply
- Pollution degree: Pollution Degree 2

Performance Specifications

The system is designed to meet specifications under these conditions:

- Ambient temperature: 15 °C (59 °F) to 30 °C (86 °F)
Over time, the temperature must remain within a range of 4 °C (7.2 °F), with the rate of the change in temperature not exceeding 2 °C (3.6 °F) per hour. Ambient temperature fluctuations exceeding the limits might result in mass shifts in spectra.
- Relative humidity: 20% to 80%, non-condensing

Equipment Use and Modification



WARNING! Personal Injury Hazard. Contact the SCIEX representative if product installation, adjustment, or relocation is required.



WARNING! Electrical Shock Hazard. Do not remove the covers. Removing the covers might cause injury or malfunctioning of the system. The covers need not be removed for routine maintenance, inspection, or adjustment. Contact a SCIEX Field Service Employee (FSE) for repairs that require the covers to be removed.



WARNING! Personal Injury Hazard. Use SCIEX-recommended parts only. Use of parts not recommended by SCIEX or use of parts for any purpose other than their intended purpose can put the user at risk of harm or negatively impact system performance.

Use the mass spectrometer and ion source indoors in a laboratory that complies with the environmental conditions recommended in the *Site Planning Guide*.

If the mass spectrometer and ion source is used in an environment or in a manner not prescribed by the manufacturer, then the protection provided by the equipment might be impaired.

Unauthorized modification or operation of the mass spectrometer and ion source might cause personal injury and equipment damage, and might void the warranty. Erroneous data might be generated if the mass spectrometer and ion source is operated either above or below the recommended environmental conditions or operated with unauthorized modifications. Contact an FSE for information on servicing the system.

Contact Us

SCIEX Support

- sciex.com/contact-us
- sciex.com/request-support

Customer Training

- In North America: NA.CustomerTraining@sciex.com
- In Europe: Europe.CustomerTraining@sciex.com
- Outside the EU and North America, visit sciex.com/education for contact information.

Online Learning Center

- [SCIEXUniversity](#)

CyberSecurity

For the latest guidance on cybersecurity for SCIEX products, visit sciex.com/Documents/brochures/win7-SecurityGuidance.pdf.

Technical Support

SCIEX and its representatives maintain a staff of fully-trained service and technical specialists located throughout the world. They can answer questions about the system or any technical issues that might arise. For more information, visit the SCIEX website at sciex.com.

Documentation Symbols and Conventions

The following symbols and conventions are used throughout the guide.



DANGER! Danger signifies an action which leads to severe injury or death.



WARNING! Warning signifies an action that could cause personal injury if precautions are not followed.

CAUTION: Caution signifies an operation that could cause damage to the system or corruption or loss of data if precautions are not followed.

Note: Note emphasizes significant information in a procedure or description.

Tip! Tip provides useful information that helps apply the techniques and procedures in the text for a specific need and provides shortcuts, but is not essential to the completion of a procedure.

Related Documentation

To find software product documentation, refer to the release notes or software installation guide that comes with the software. Documentation for the hardware products can be found on the *Customer Reference* DVD that comes with the system or component.

For the latest versions of the documentation, visit the SCIEX website at sciex.com.

Principles of Operation

2

The SCIEX Triple Quad™ 3500 system is designed for the qualitative and quantitative analysis of chemical species. This section includes information about the mass spectrometer and the Analyst® software. Refer to the ion source *Operator Guide* for an overview of the ion source.

For information on the computer and software, refer to the *Software Installation Guide* for the Analyst® software.



WARNING! Lifting Hazard. Do not move the system. Risk of personal injury or system damage. If the system must be moved, then contact a Field Service Employee (FSE).

System Overview



WARNING! Lifting Hazard. Make sure that at least six people are available to lift the mass spectrometer. Follow established safe lifting procedures. Refer to the *Site Planning Guide* for the weights of system components.


This system is intended for the qualitative and quantitative analysis of chemical species.

The SCIEX Triple Quad™ 3500 system includes the following components:

- A SCIEX Triple Quad™ 3500 mass spectrometer with a roughing pump and a source of compressed air and nitrogen
- Turbo V™ ion source that uses either the TurbolonSpray® probe or the atmospheric pressure chemical ionization (APCI) probe. Refer to the *Turbo V™ Ion Source Operator Guide*.
- SCIEX-supplied computer and monitor with the Analyst® software for instrument optimization, acquisition method development, processing, and data acquisition. For computer specifications and requirements, refer to the *Analyst® Software Installation Guide*.

Hardware Overview

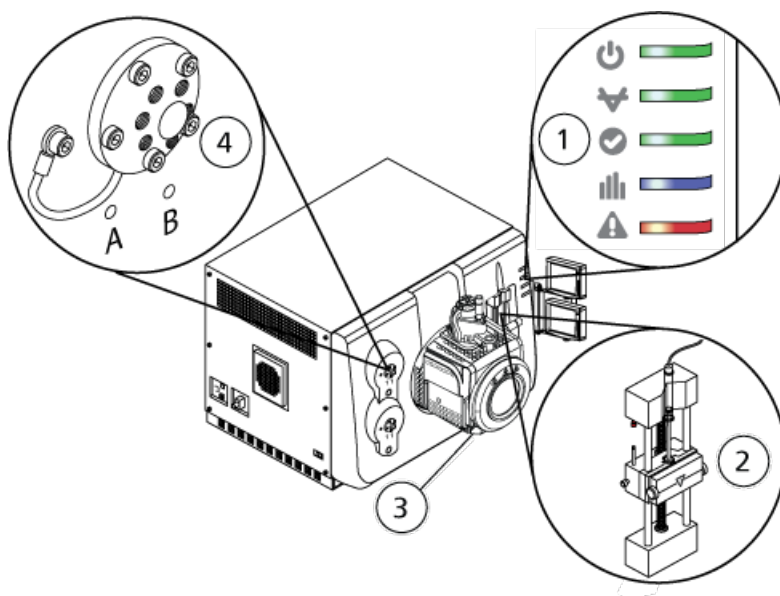Figure 2-1 Front View








Item	Description	Refer to...
1	Panel symbols	Panel Symbols on page 21.
2	Syringe pump	Adjust the Integrated Syringe Pump Position on page 46
3	Ion source	<i>Turbo VTM Ion Source Operator Guide</i> , available from the ion source documentation DVD, or from the SCIEX Web site, at sciex.com
4	Diverter valve	Plumb the Diverter Valve on page 50.

Panel Symbols

[Table 2-1](#) describes the mass spectrometer status LEDs.

Principles of Operation

Table 2-1 Panel Symbols


LED	Color	Name	Description
	Green	Power	Lit when the system is powered up.
	Green	Vacuum	Lit when the correct vacuum level has been achieved. Flashing if the vacuum is not at the correct level (during pump down and venting).
	Green	Ready	Lit when the system is in the Ready state. The system must be in the Ready state to operate.
	Blue	Scanning	Flashing when the system is acquiring data.
	Red	Fault	Lit when the system encounters a system fault.

After the system is turned on, all of the LEDs illuminate. The power LED remains lit. The other LEDs flash for two seconds and then turn off. The vacuum LED starts flashing. After the correct vacuum level is achieved this LED remains lit.

Connections

Figure 2-2 shows the location of the mass spectrometer connections, including the locations of the **RESET** and **VENT** buttons and the mass spectrometer convenience switch.

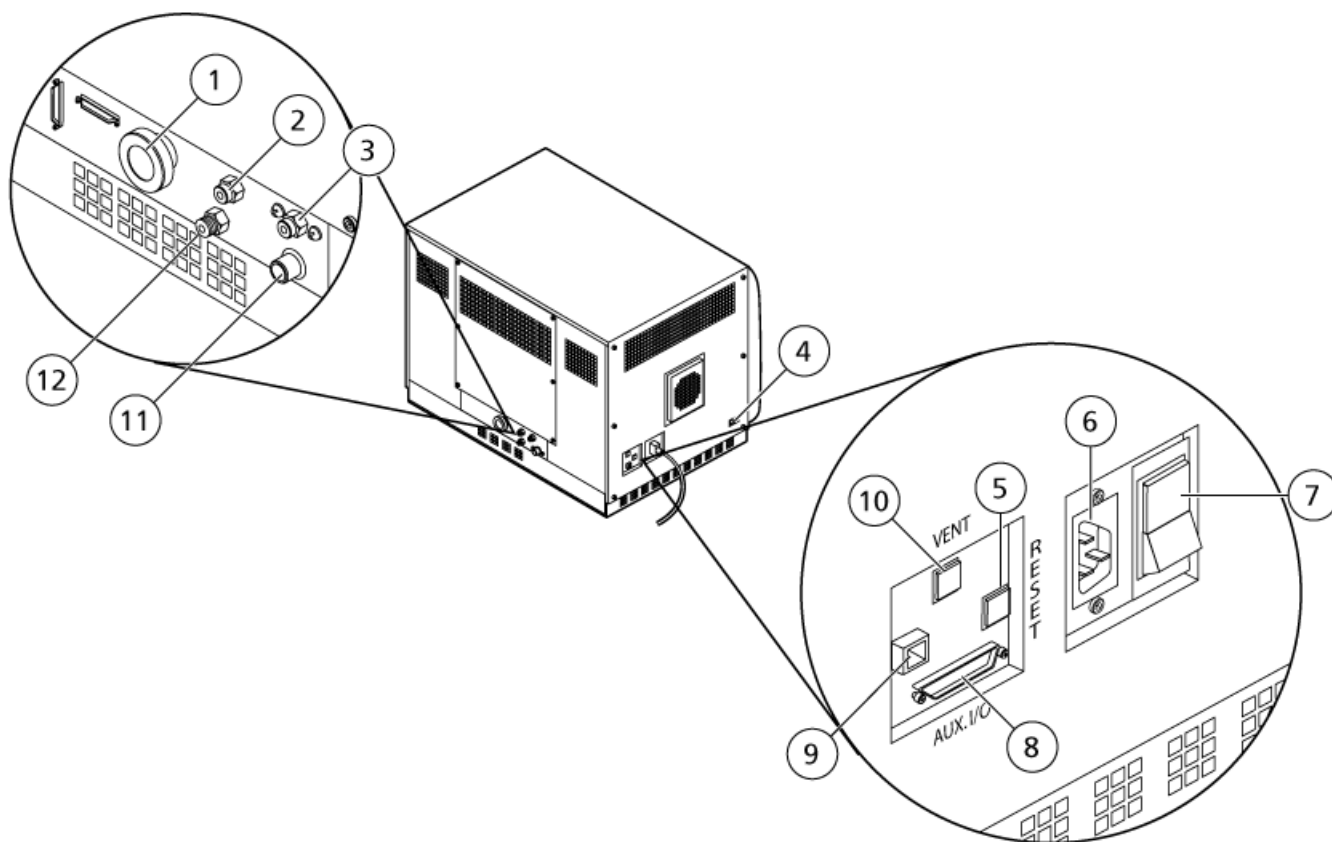Figure 2-2 Back and Side Views



Item	Description	Primary Materials	For more information...
1	Roughing pump vacuum connection	Aluminum (hose fitting), zinc-plated steel (hose clamp)	Contact an FSE.
2	Air supply (Gas 1/Gas 2)	Plastic	Refer to the <i>Site Planning Guide</i> . Gas tubing from the Gas 1/Gas 2 supply connects to this port.
3	Source exhaust supply	Stainless steel	Refer to the <i>Site Planning Guide</i> . Gas tubing from the source exhaust gas supply connects to this port.
4	Source communication connection	Aluminum	Contact an FSE.
5	RESET button	Plastic	Refer to Reset the Mass Spectrometer on page 45 .

Principles of Operation


Item	Description	Primary Materials	For more information...
6	Mains supply connection	Aluminum/plastic	Refer to Start Up the System on page 43 or Shut Down and Vent the System on page 45 .
7	Mass spectrometer convenience switch (Up = On; Down = Off)	Plastic	Refer to Start Up the System on page 43 or Shut Down and Vent the System on page 45 .
8	Aux I/O connection	Sheet metal (zinc-plated)	Refer to the <i>Peripheral Devices Setup Guide</i> .
9	Ethernet connection (connects the mass spectrometer and the computer)	Sheet metal (zinc-plated)	Contact an FSE.
10	VENT button	Plastic	Refer to Start Up the System on page 43 or Shut Down and Vent the System on page 45 .
11	Source exhaust waste (to waste bottle)	Stainless steel	Refer to the <i>Site Planning Guide</i> .
12	Nitrogen gas supply (Curtain Gas™ supply, CAD gas)	Stainless steel	Refer to the <i>Site Planning Guide</i> . Gas tubing from the nitrogen gas supply connects to this port.

Theory of Operation—Hardware

Mass spectrometry measures the mass-to-charge ratio of ions to identify unknown compounds, to quantify known compounds, and to provide information about the structural and chemical properties of molecules.

The mass spectrometer has a series of quadrupole filters that transmit ions according to their mass-to-charge (m/z) ratio. The first quadrupole in this series is the QJet® ion guide located between the orifice plate and the Q0 region. The QJet® ion guide does not filter ions, but focuses them before they enter the Q0 region. By prefocusing the larger ion flux created by the wider orifice, the QJet® ion guide increases system sensitivity and improves the signal-to-noise ratio. In the Q0 region, the ions are again focused before passing into the Q1 quadrupole.

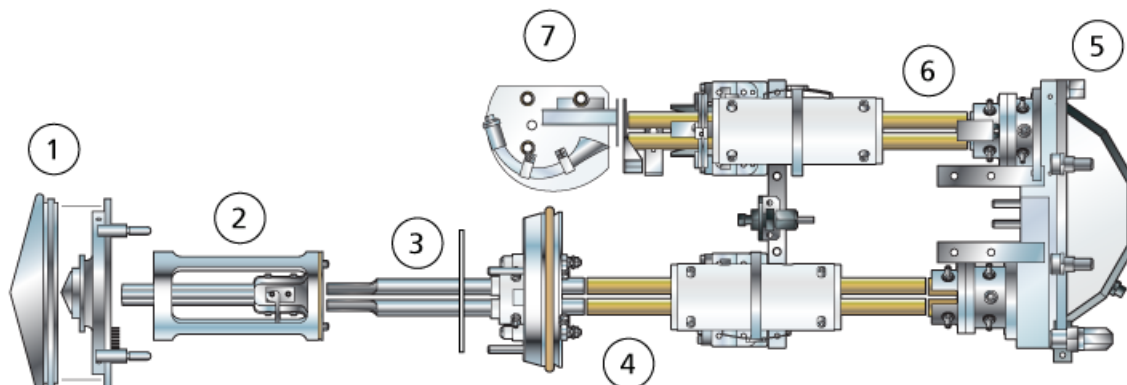Figure 2-3 Ion Path



Item	Description
1	Orifice plate
2	QJet ion guide
3	Q0 region
4	Q1 quadrupole
5	Q2 collision cell
6	Q3 quadrupole
7	Detector

The Q1 quadrupole is a filtering quadrupole that sorts the ions before they enter the Q2 collision cell. The Q2 collision cell is where the internal energy of an ion is increased through collisions with gas molecules to the point that molecular bonds break creating product ions. This technique allows users to design experiments that measure the m/z of product ions to determine the composition of the parent ions.

After passing through the Q2 collision cell, the ions enter the Q3 quadrupole for additional filtering, and then they enter the detector. In the detector, the ions create a current that is converted into a voltage pulse. The voltage pulses leaving the detector are directly proportional to the quantity of ions entering the detector. The system monitors these voltage pulses and then converts the information into a signal. The signal represents the ion intensity for a particular m/z value and the system shows this information as a mass spectrum.

Data Handling

The Analyst[®] software requires a computer running the Windows 7 (32- or 64-bit) operating system or the Windows 10 (64-bit) operating system. The computer and the associated system software work with the system controller and the associated firmware to control the system and data acquisition. During system operation, the acquired

Principles of Operation

data is sent to the Analyst[®] software where it can be shown as either full mass spectra, intensity of single or multiple ions over time, or total ion current over time.

Analyst[®] Software Overview


The Analyst[®] software works with the mass spectrometer and the liquid chromatography (LC) system and the associated firmware to control the system and data acquisition. During system operation, the acquired data is sent to the Analyst[®] software where it can be shown as full mass spectra, intensity of single or multiple ions over time, or total ion count over time.

Different Data Views

The following figures show examples of two types of data views: total ion chromatogram (TIC) and extracted ion chromatogram (XIC).


TIC: The plot of the total ion current as a function of time.

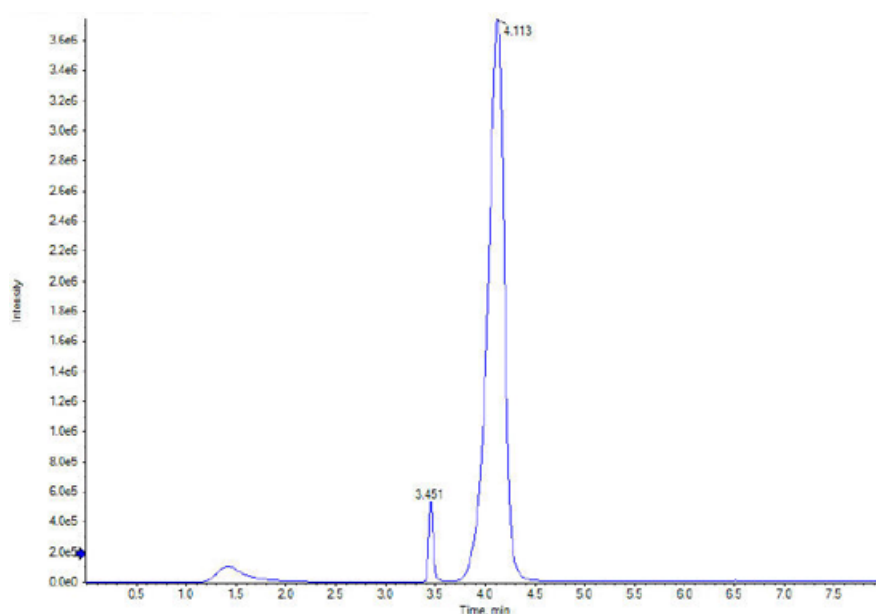Figure 2-4 Example TIC



XIC: An ion chromatogram created by taking intensity values at a single, discrete mass value or a mass range, from a series of mass spectral scans. It indicates the behavior of a given mass or mass range as a function of time.

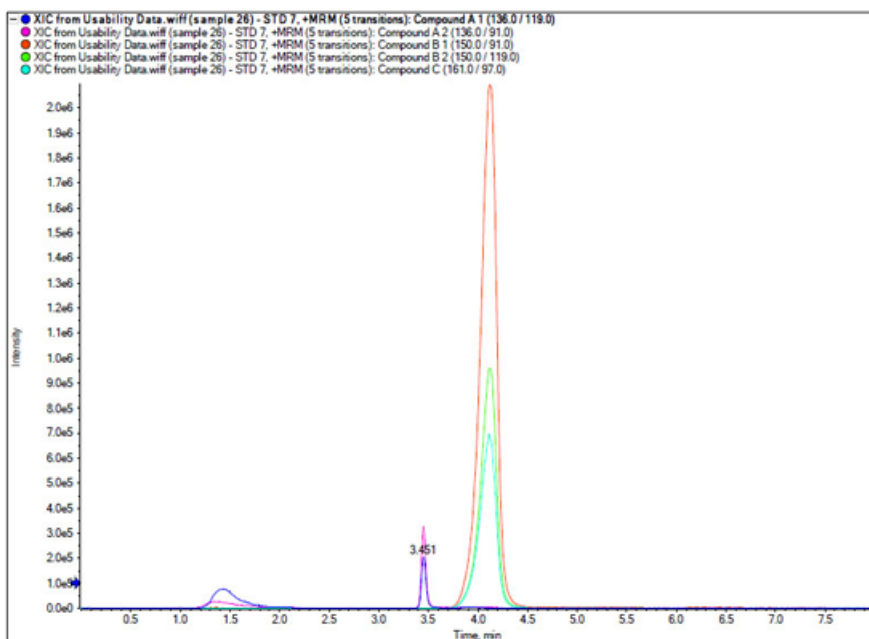Figure 2-5 Example XIC




Parameters

The working parameters are the set of instrument parameters currently being used.

- Source and gas parameters: These parameters can change depending on the ion source used.
- Compound parameters: These parameters consist mostly of voltages in the ion path. Optimal values for compound-dependent parameters vary depending on the compound being analyzed.
- Resolution parameters: These parameters affect the resolution and calibration.
- Detector parameters: These parameters affect the detector.

The following figure shows the location of the parameters on the ion optics path.

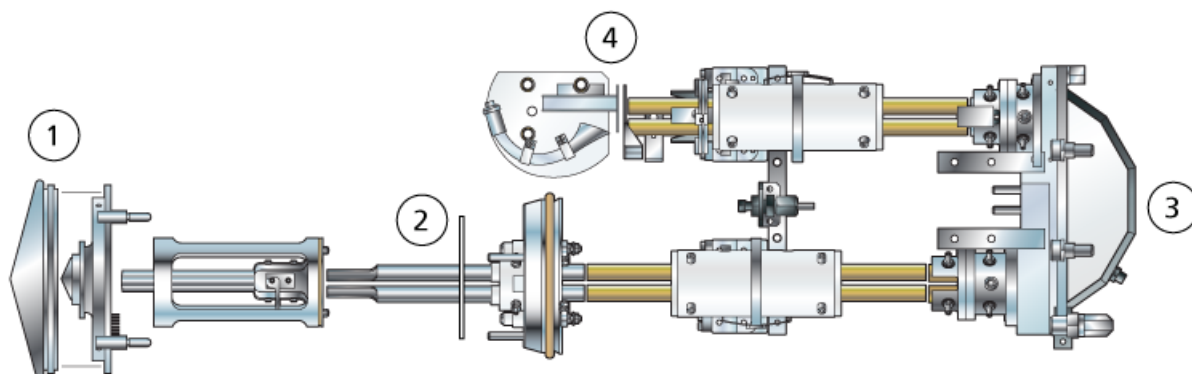Figure 2-6 Ion Optics Path and Parameters



Location	Parameter	Parameter Type	Use	Scan Type
1	IonSpray Voltage (IS)	Source and gas	The IS parameter controls the voltage applied to the electrode in the TurbolonSpray [®] probe that ionizes the sample in the ion source. The parameter depends on the polarity and it affects the spray stability and the sensitivity. This parameter is also used with the NanoSpray [®] ion source. Refer to the applicable ion source <i>Operator Guide</i> . The parameter can be compound-dependent and should be optimized for each compound.	All
1	Nebulizer Current (NC)	Source and gas	The NC parameter controls the current applied to the corona discharge needle in the APCI probe, used in the Turbo V [™] or DuoSpray [™] ion source. The discharge ionizes solvent molecules, which in turn ionize the sample molecules.	All

Location	Parameter	Parameter Type	Use	Scan Type
1	Interface Heater Temperature (IHT)	Source and gas	The IHT parameter controls the temperature of the NanoSpray [®] interface heater and is only available if the NanoSpray [®] ion source and interface are installed. The optimal heater temperature depends on the type of sample being analyzed and the solvent used. If the heater temperature is too high, the signal degrades. Typically, heater temperatures are in the 130 °C to 180 °C range. The maximum heater temperature that can be set is 250 °C, but this is too high for most applications.	All
1	Ion Source Gas 1 (GS1)	Source and gas	The GS1 parameter controls the nebulizer gas for both the TurbolonSpray [®] and APCI probes. The GS1 parameter controls the nebulizer gas for the TurbolonSpray [®] probe.	All
1	Ion Source Gas 2 (GS2)	Source and gas	The GS2 parameter controls the heater gas for the TurbolonSpray [®] probe. The GS2 parameter controls the heater gas for the TurbolonSpray [®] probe and the nebulizer gas for the APCI probe.	All
1	Temperature (TEM)	Source and gas	The TEM parameter controls the temperature of the heater gas for the TurbolonSpray [®] probe or the temperature of the APCI probe.	All
1	Curtain Gas (CUR)	Source and gas	The CUR parameter controls the gas flow of the Curtain Gas [™] interface. The Curtain Gas [™] interface is located between the curtain plate and the orifice. It prevents the contamination of the ion optics.	All
1	sdp	Source and gas	The sdp parameter controls the selection of the DuoSpray [™] ion source probes: TurbolonSpray [®] or APCI probe.	n/a

Principles of Operation

Location	Parameter	Parameter Type	Use	Scan Type
1	Declustering Potential (DP)	Compound	<p>The DP parameter controls the voltage on the orifice, which controls the ability to decluster ions between the orifice and the QJet[®] ion guide. It is used to minimize the solvent clusters that might remain on the sample ions after they enter the vacuum chamber, and, if required, to fragment ions. The higher the voltage, the higher the energy imparted to the ions. If the DP parameter is too high, then unwanted fragmentation might occur.</p> <p>Use the preset value and optimize for the compound.</p>	All
2	Entrance Potential (EP)	Compound	<p>The EP parameter controls the potential difference between the voltage on Q0 and ground. The entrance potential guides and focuses the ions through the high-pressure on the Q0 region.</p> <p>Use the preset value.</p>	All
3	CAD Gas	Source and gas	<p>The CAD parameter controls the pressure of the CAD gas in the collision cell during Q3, MS/MS. For Q3 scans, the collision gas helps to focus the ions as they pass through the Q2 collision cell. The preset value for the CAD parameter is in fixed mode. For MS/MS scan types, the CAD gas helps to fragment the precursor ions. When the precursor ions collide with the collision gas, they dissociate to form product ions.</p> <p>Use the preset value and optimize for the compound.</p>	Q3 MI, Q3 MS, MRM, Prec, NL

Location	Parameter	Parameter Type	Use	Scan Type
3	Collision Energy (CE)	Compound	The CE parameter controls the potential difference between the Q0 region and the Q2 collision cell. It is used only in MS/MS scan types. This parameter is the amount of energy that the precursor ions receive as they are accelerated into the Q2 collision cell, where they collide with gas molecules and fragment. Use the preset value and optimize for the compound.	EPI, MS/MS/MS, MRM, MS2, Prec, NL
3	Collision Cell Exit Potential (CXP)	Compound	The CXP parameter is only used in Q3 and MS/MS scan types. This parameter transmits the ions into the Q3 quadrupole. Use the preset value and optimize for the compound.	Q3, MRM, MS2, Prec, NL
4	CEM	Detector	The CEM parameter controls the voltage applied to the detector. The voltage controls the detector response.	All

Theory of Operation—Analyst[®] Software

Quantitation

The goal of quantitation in LC-MS/MS is to accurately determine the concentration of a compound in an unknown sample. The MRM scan type is primarily used for quantitative analysis in a triple quadrupole or QTRAP[®] mass spectrometer. In an MRM scan, the ability to define the precursor ion and characteristic product ion enables the creation a pair that is highly specific of the analyte. The MRM transition (pair) coupled with the retention time associated with the analyte during liquid chromatography provides the specificity required for quantitation.

Quantitation is accomplished through the use of validated MRM LC-MS/MS acquisition methods, acquisition of calibration standard curves, and the subsequent integration of the peaks associated with the compounds of interest. The calibration curve relationship between signal response and concentration is used to determine the quantity of a particular analyte in an unknown sample.

Integration

In LC-MS/MS data, integration refers to obtaining the area under a curve for the peak associated with a specific compound. Through the development of a method which specifies the transitions, expected retention times, internal standards, integration, and regression parameters, the software is able to automatically integrate peaks for a given set of samples.

About Results Tables

Results tables summarize the calculated concentration of an analyte in each unknown sample based on the calibration curve. Results tables also include the calibration curves as well as statistics for the results. The user can customize the results tables and view the results tables in layouts.

The data from a results tables can be exported to a txt file for use in other applications, such as Microsoft Excel. The user can also export data in the table or just the data in the visible columns.

Calibration Curves

A calibration curve, also known as a standard concentration curve, is a method for determining the concentration of a substance in an unknown sample by comparing the unknown sample to a set of standard samples of known concentration. The calibration curve is a plot of how the instrument responds (the analytical signal) to changes in the concentration of the analyte (the substance to be measured). The operator prepares a series of standards across a range of concentrations near the expected concentration of the analyte in the unknown sample.

Calibration standards are used to build calibration curves. Incorrect readings or missing readings on some of the calibration samples might indicate issues with the analytical run. Follow acceptable methods found in literature and regulatory agency guidances to create a calibration curve. Examples of good practice in the preparation of calibration curves include:

- Preparing calibration standards in blank matrix in which the analyte is to be measured.
- Generating a calibration curve for each analyte to be measured.
- Making sure of the coverage of the expected concentration range of the analyte, including typical and atypical specimens.
- Using six to eight standards to generate the curve.

This is not a comprehensive list and other guidances should be used in determining the best practice in developing a calibration curve for the laboratory.

Note: In some analytical runs, single-point calibration standards are used. Single-point calibrations are performed using a matrix blank sample and a single standard concentration. The relationship between instrument response and analyte concentration is determined by the line created by these two points. Both the acquisition and quantitation methods should be validated before being accepted for their intended use.

Regression

The area of the analyte peaks in the calibration curve standards is plotted against the known concentrations. Subsequently, a line is fitted to the points. This regression line is used to calculate the concentration of the unknown samples.

- Linear ($y = mx + b$)
- Linear through Zero ($y = mx$)
- Quadratic ($y = a_2 + bx + c$)

As well, it is possible to add different types of weighting for the regression, including:

- $1/x$
- $1/x^2$
- $1/y$
- $1/y^2$

Operating Instructions — Sample Workflows

3

Table 3-1 Instrument Setup Workflow

Step	To do this...	Find the information in....	What does it do?
1	Create a hardware profile.	Create a Hardware Profile on page 54	Each hardware profile must include a mass spectrometer and other devices, such as an LC system. Only devices included in the active hardware profile can be used when creating acquisition methods.
2	Create projects to store data.	Create Projects and Subprojects on page 61	Using projects and subprojects facilitates data management and makes comparison of results easier.
3	Optimize the mass spectrometer.	Verify Instrument Performance on page 65	This is the process of optimizing the resolution and mass spectrometer parameters, and calibrating the mass spectrometer to obtain the best sensitivity and performance from the system.

Table 3-2 Sample Acquisition Workflow

Step	To do this...	Find the information in....	What does it do?
1	Create projects to store data.	Create Projects and Subprojects on page 61	Using projects and subprojects facilitates data management and makes comparison of results easier.
2	Create an acquisition method.	Operating Instructions — Acquisition Methods on page 84	To analyze samples, create an acquisition method for the mass spectrometer and any LC devices. An acquisition method indicates which peripheral devices to use, when to use them to acquire data, and the associated parameters.
3	Create and submit a batch.	Add Sets and Samples to a Batch on page 135 and Submit a Sample or Set of Samples on page 139	After creating an acquisition method, run samples by creating an acquisition batch and submitting the batch to the acquisition queue.
4	Run samples to acquire data.	Acquire Data on page 139	Running samples involves managing the acquisition queue and monitoring instrument and device status. To submit samples and acquire data, use the Queue Manager. The Queue Manager shows queue, batch, and sample status, and facilitates management of samples and batches in the queue.

Table 3-2 Sample Acquisition Workflow (continued)

Step	To do this...	Find the information in....	What does it do?
5	Analyze data in Explore mode. —OR— Step 6—Analyze quantitative data.	Operating Instructions — Analyze and Explore Data on page 110	In Explore mode, many tools are available for viewing and processing the acquired data. Graphs can be customized with peak labels and captions, contour plots can be shown, and spectra can be saved in the library.
6	Analyze quantitative data.	Operating Instructions — Analyze and Process Quantitative Data on page 90	Use the various quantitative method creation tools in Quantitate mode to analyze the acquired data and build a quantitative method to generate a Results Table. Use the Results Table to manually review all of the peaks for each analyte and internal standard within a batch and to view calibration curves, sample statistics, and metric plots.

Table 3-3 Example Routine Analysis Workflow

Step	To do this...	Find the information in....	What does it do?
1	Activate a hardware profile applicable to the method.	Create a Hardware Profile on page 54	Each hardware profile must include a mass spectrometer and other devices, such as an LC system. Only devices included in the active hardware profile can be used when creating acquisition methods.
2	Create projects to store data.	Create Projects and Subprojects on page 61	Using projects and subprojects facilitates data management and makes comparison of results easier.
3	Create and submit a batch.	Add Sets and Samples to a Batch on page 135 and Submit a Sample or Set of Samples on page 139	After creating an acquisition method, run samples by creating an acquisition batch and submitting the batch to the acquisition queue.
4	Equilibrate the system.	Equilibrate the System on page 138	Equilibrate the system prior to starting data acquisition. A system that is not equilibrated can result in poor data.
5	Run samples to acquire data.	Acquire Data on page 139	Running samples involves managing the acquisition queue and monitoring instrument and device status. To submit samples and acquire data, use the Queue Manager. The Queue Manager shows queue, batch, and sample status, and facilitates management of samples and batches in the queue.

Table 3-3 Example Routine Analysis Workflow (continued)

Step	To do this...	Find the information in....	What does it do?
6	Analyze data in Explore mode (optional).	Operating Instructions — Analyze and Explore Data on page 110	In Explore mode, many tools are available for viewing and processing the acquired data. Graphs can be customized with peak labels and captions, contour plots can be shown, and spectra can be saved in the library.
<ul style="list-style-type: none"> For quantitative analysis using the MultiQuant™ software follow steps 7 to 8. The MultiQuant™ software is recommended for larger data sets. For quantitative analysis using the Analyst® software follow steps 9 to 10a For qualitative analysis using the Analyst® software follow steps 9 to 10b. 			
7	Analyze quantitative data in the MultiQuant™ software.	<i>MultiQuant™ Software Reference Guide</i> : Chapter 7, 8, 10, 11, 12, 13, 14	Generate and use a Results Table to manually review all of the peaks for each analyte and internal standard within a batch and to view calibration curves, sample statistics, and metric plots.
8	Create a report in the MultiQuant™ software.	<i>MultiQuant™ Software Reference Guide</i> : Appendix C	Generate a report using the provided report templates for the generated and reviewed results.
9	Analyze qualitative (or quantitative) data in the Analyst® software.	Operating Instructions — Analyze and Process Quantitative Data on page 90	Generate a Results Table to manually review all of the peaks for each analyte and internal standard within a batch. For quantitative analysis, also review calibration curves, sample statistics, and metric plots.

Table 3-3 Example Routine Analysis Workflow (continued)

Step	To do this...	Find the information in....	What does it do?
10a	Create a report in the Analyst Reporter.	Generate Reports	Generate a report using the provided report templates for the generated and reviewed results. For reports specific to qualitative analysis, use the set of report templates labeled with library searching.
10b	Select a library and then create a report using the Analyst Reporter.	Generate Reports	Select the appropriate MS/MS spectra library for the results and generate a report using the provided report templates labeled as library search for the generated and reviewed results. For reports specific to qualitative analysis, use the set of report templates labeled with library search.

Table 3-4 Example Method Developer Workflow

Step	To do this...	Find the information in....	What does it do?
1	Create a hardware profile.	Create a Hardware Profile on page 54	Each hardware profile must include a mass spectrometer and other devices, such as an LC system. Only devices included in the active hardware profile can be used when creating acquisition methods.
2	Create projects to store data.	Create Projects and Subprojects on page 61	Using projects and subprojects facilitates data management and makes comparison of results easier.

Table 3-4 Example Method Developer Workflow (continued)

Step	To do this...	Find the information in....	What does it do?
3	Optimize the compound automatically - Or - Step 4 Optimize the compound manually.	Automatic Optimization on page 68	The software automatically optimizes the compound and mass spectrometer parameters for the compounds of interest.
4	Optimize the compound manually.	Operating Instructions — Manual Compound Optimization on page 163	The user manually optimizes the compound and mass spectrometer parameters for the compounds of interest. Manual optimization provides more experienced users with greater control during the optimization process.
5	Create an acquisition method.	Operating Instructions — Acquisition Methods on page 84	To analyze samples, create an acquisition method for the mass spectrometer and any LC devices. An acquisition method indicates which peripheral devices to use, when to use them to acquire data, and the associated parameters.
6	Create and submit a batch.	Add Sets and Samples to a Batch on page 135 and Submit a Sample or Set of Samples on page 139	After creating an acquisition method, run samples by creating an acquisition batch and submitting the batch to the acquisition queue.
7	Equilibrate the system.	Equilibrate the System on page 138	Equilibrate the system prior to starting data acquisition. A system that is not equilibrated can result in poor data.

Table 3-4 Example Method Developer Workflow (continued)

Step	To do this...	Find the information in....	What does it do?
8	Run samples to acquire data.	Acquire Data on page 139	Running samples involves managing the acquisition queue and monitoring instrument and device status. To submit samples and acquire data, use the Queue Manager. The Queue Manager shows queue, batch, and sample status, and facilitates management of samples and batches in the queue.
9	Analyze data in Explore mode (optional).	Operating Instructions — Analyze and Explore Data on page 110	In Explore mode, many tools are available for viewing and processing the acquired data. Graphs can be customized with peak labels and captions, contour plots can be shown, and spectra can be saved in the library.
<ul style="list-style-type: none"> For quantitative analysis using the MultiQuant™ software follow steps 10 to 12. The MultiQuant™ software is recommended for larger data sets. For quantitative analysis using the Analyst® software follow steps 13 to 15. For qualitative analysis, contact Support. 			
10	Create a quantitative method in the MultiQuant™ software.	<i>MultiQuant™ Software Reference Guide: Quantitation Method Editor</i>	Use the various quantitative method creation tools in the software to analyze the acquired data and build a quantitative method to generate a Results Table.

Table 3-4 Example Method Developer Workflow (continued)

Step	To do this...	Find the information in....	What does it do?
11	Analyze quantitative data in MultiQuant™ software.	<i>MultiQuant™ Software Reference Guide</i> : Chapter 7, 8, 10, 11, 12, 13, 14	Generate and use a Results Table to manually review all of the peaks for each analyte and internal standard within a batch and to view calibration curves, sample statistics, and metric plots.
12	Create a report in MultiQuant™ software.	<i>MultiQuant™ Software Reference Guide</i> : Appendix C	Generate a report using the provided report templates for the generated and reviewed results.
13	Create a quantitative method in the Analyst® software.	<i>MultiQuant™ Software Reference Guide</i> : Quantitation Method Editor	Use the various quantitative method creation tools in the software to analyze the acquired data and build a quantitative method to generate a Results Table.
14	Analyze qualitative (or quantitative) data in Analyst® software.	Operating Instructions — Analyze and Process Quantitative Data on page 90	Generate a Results Table to manually review all of the peaks for each analyte and internal standard within a batch. For quantitative analysis, also review calibration curves, sample statistics, and metric plots.
14	Create a report in the Analyst Reporter.	Generate Reports	Generate a report using the provided report templates for the generated and reviewed results. For reports specific to qualitative analysis, use the set of report templates labeled with library searching.

Operating Instructions—Hardware

4



WARNING! Personal Injury Hazard. Follow the instructions in the documentation when using the system. The protection provided by the equipment might be impaired if the equipment is used in a manner not specified by SCIEX.

Start Up the System



WARNING! Electrical Shock Hazard. Make sure that the system can be disconnected from the mains supply outlet in an emergency. Do not block the mains supply outlet.



WARNING! Lifting Hazard. Do not move the system. Risk of personal injury or system damage. If the system must be moved, then contact a Field Service Employee (FSE).

Note: Before operating the instrument, read the safety information in [Operational Precautions and Limitations on page 8](#).


Prerequisites

- The site requirements specified in the *Site Planning Guide* are met. The *Site Planning Guide* includes information on the mains supply and connections, compressed air, nitrogen, roughing pump, ventilation, exhaust, and site clearance requirements. Contact us for a copy of the *Site Planning Guide*, if required. For contact information, go to sciex.com/contact-us.
- The source exhaust gas, compressed air, and nitrogen gases are connected to the mass spectrometer.
- The 4 L source exhaust drain bottle is connected to the exhaust waste connection on the back of the mass spectrometer and to the laboratory ventilation system.
- The source exhaust hoses are securely clamped at the mass spectrometer, drain bottle, and ventilation connections.
- The mass spectrometer convenience switch is turned off and the mains supply cable is plugged into the mass spectrometer.
- The mass spectrometer and roughing pump mains supply cables are plugged into the 200 VAC to 240 VAC mains supply.
- The Ethernet cable is connected to both the mass spectrometer and the computer.

1. Turn on the roughing pump.

The On/Off switch is beside the mains supply input connection on the roughing pump.

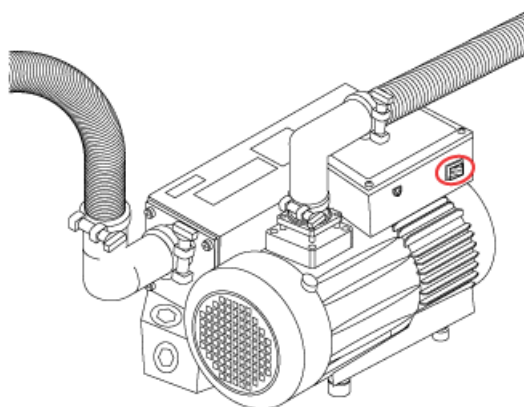Figure 4-1 Roughing Pump—On/Off Switch



2. Wait five minutes, and then turn on the mass spectrometer convenience switch.
3. Turn on the computer.
4. Open the Analyst[®] software.

Reset the Mass Spectrometer

1. Stop any ongoing scans and then turn off sample flow to the mass spectrometer.
2. In the Analyst[®] software, deactivate the hardware profile, if it is active.
3. Press and hold the **Reset** button for five seconds.

A click is heard when the relay activates. After approximately three minutes, the mass spectrometer reaches operating pressure.

Shut Down and Vent the System

Some procedures require that the system be shut down. Others require that it also be vented. Follow these steps to shut down and, if required, vent the system.

Note: If the input gas supply must be disconnected, then relieve the pressure in the gas lines before disconnecting it.

Tip! If the mass spectrometer will not be used for a length of time, then leave it in Standby mode with the ion source in place. If the mass spectrometer must be shut down, then follow these instructions. Do not turn off the roughing pump until after the turbo pumps have spun down.

1. Complete or stop any ongoing scans.

CAUTION: Potential System Damage. Turn off the sample flow before shutting down the system.

2. Turn off the sample flow to the system.
3. In the Analyst[®] software, deactivate the hardware profile, if it is active.
4. Close the software.
5. (If required) Follow these steps to vent the system:

Note: Vent the system before performing a full cleaning of the vacuum interface, before cleaning the Q0 region, and before replacing the roughing pump oil. For more information contact the Qualified Maintenance Person (QMP) or FSE.

- a. Press and hold the **Vent** button for three seconds.

The Vacuum LED begins flashing quickly (more quickly than during pump down). The turbo pump spins down gradually.
- b. Allow the system to vent for 15 minutes and then turn off the roughing pump.
6. Turn off the mass spectrometer convenience switch.
7. Disconnect the mass spectrometer mains supply cable from the mains supply outlet.
8. (If venting the system) Disconnect the roughing pump mains supply cable from the mains supply outlet.

Adjust the Integrated Syringe Pump Position




WARNING! Puncture Hazard. Take care when handling the syringe. The tip of the syringe is extremely sharp.



WARNING! Puncture Hazard. Make sure that the syringe is seated properly in the syringe pump and that the automatic syringe pump stop is adjusted properly to avoid damaging or breaking the glass syringe. If the syringe breaks, follow established safety procedures for sharps disposal.


1. Press the **Release** button on the right side of the syringe pump to lower the base and then insert the syringe. Refer to [Figure 4-2](#).

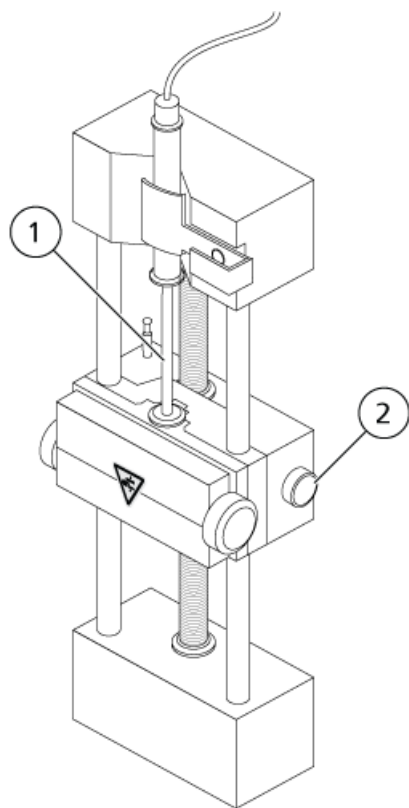Figure 4-2 Lowering the Syringe

Item	Description
1	Syringe plunger
2	Release button. Press to raise or lower the base.

2. Make sure that the end of the syringe is flush against the base and that the shaft of the syringe rests in the cutout.
3. Adjust the post so that it triggers the automatic syringe stop before the syringe plunger hits the bottom of the glass syringe. Refer to [Figure 4-3](#).

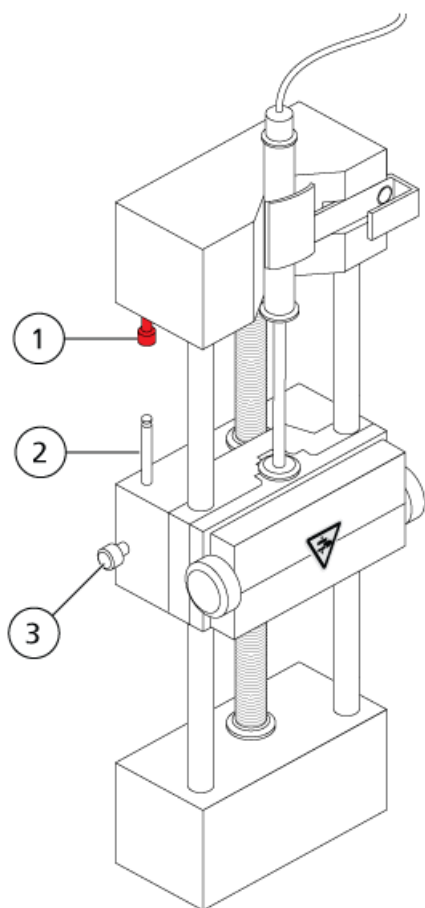
Figure 4-3 Automatic Syringe Stop



Item	Description
1	Automatic syringe stop. After the post hits the automatic syringe stop, the syringe pump stops.
2	Post. Adjust the height to prevent the syringe plunger from hitting the syringe during sample infusion.
3	Post lock screw. Tighten the screw after the height of the post is adjusted.

4. Turn the side screws as shown in [Figure 4-4](#) to secure the syringe.

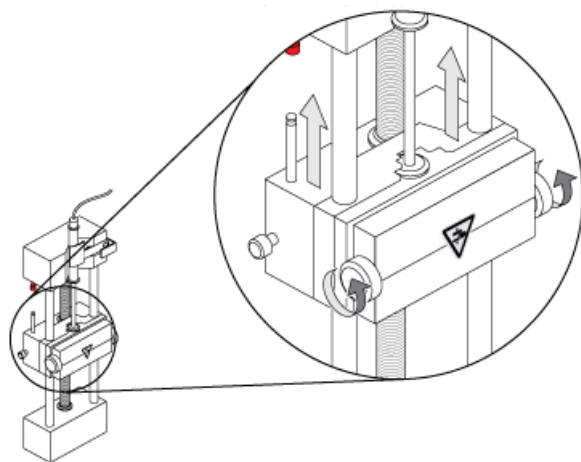
Figure 4-4 Syringe Pump Screws



5. Make sure that the mass spectrometer and integrated syringe pump are activated in the software.

Note: For subsequent manual use, on the mass spectrometer, press the button on the right side of the syringe pump to start the flow. Refer to [Figure 4-5](#). The LED beside the button flashes when the syringe pump is in use. The syringe pump flow can also be controlled automatically by the Analyst[®] software.

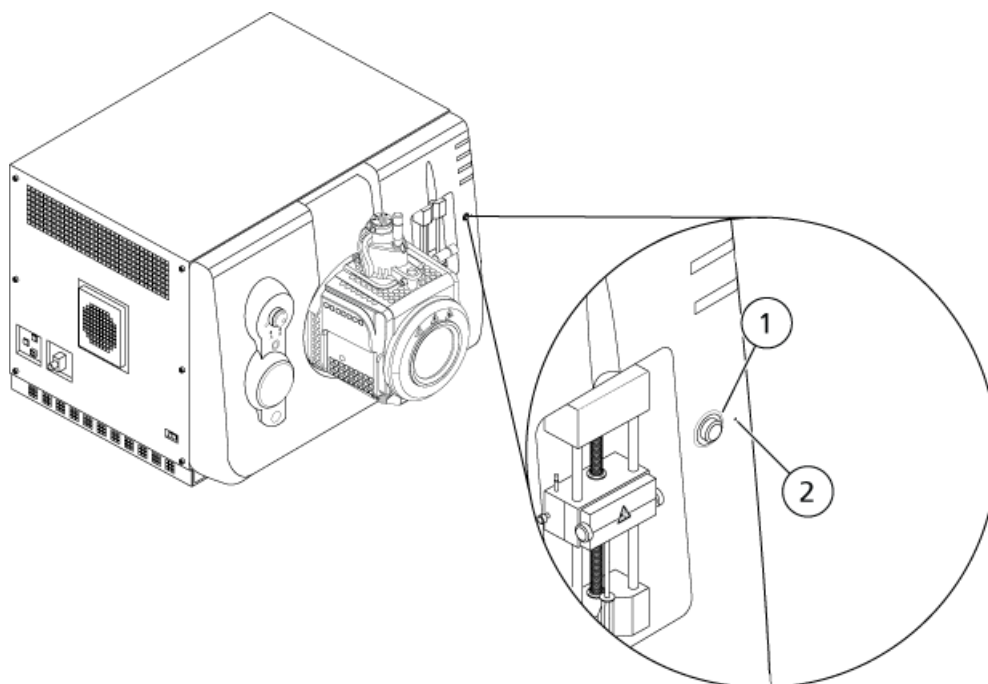Figure 4-5 Syringe Pump LED



Item	Description
1	Syringe pump on and off button
2	Syringe pump status LED

6. In the Analyst[®] software, on the Navigation bar, double-click **Manual Tuning**.
7. Click **Start Syringe**.
8. To stop the syringe pump, click **Stop Syringe**.

Plumb the Diverter Valve

The integrated diverter valve, which is located next to the ion source, can be plumbed in injector mode or diverter mode. To configure the valve, access the **Configuration** tab and then make sure that the **Use integrated injector/diverter valve** check box is selected. Refer to [Add Devices to a Hardware Profile on page 58](#).

CAUTION: Potential Wrong Result. Do not press the diverter valve button during a run. Doing so might result in incorrect data.

Plumb the Diverter Valve in Injector Mode

When the valve is in Position A, the sample flows through the external loop. When the valve switches to Position B, the sample is injected.

- Plumb the valve for injector mode.

Figure 4-6 Diverter Valve—Injector Mode Position A


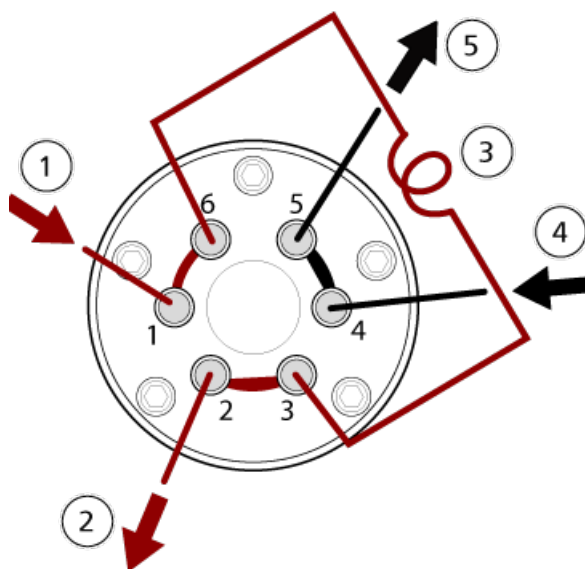


Figure 4-7 Diverter Valve—Injector Mode Position B

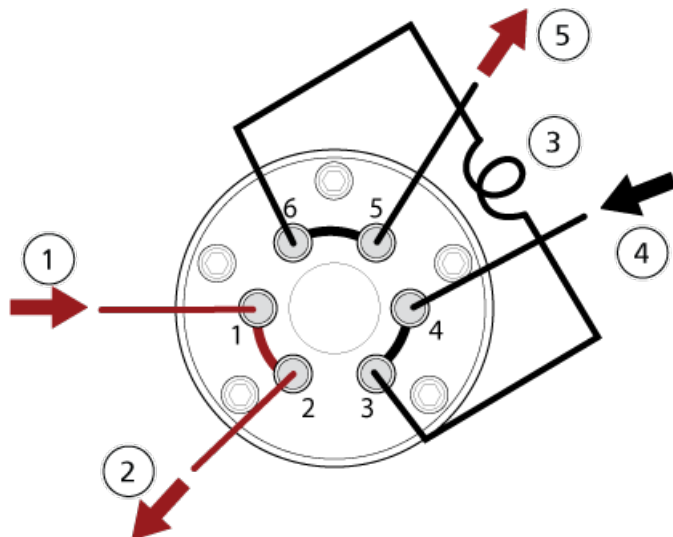

Item	Description
1	Sample in
2	Waste out
3	Sample loop (ports 3 and 6)
4	Mobile phase in
5	To column (or to the mass spectrometer, if a column is not installed)

Plumb the Diverter Valve in Diverter Mode

When the valve is in Position A, the flow goes to waste. When the valve switches to Position B, the flow goes to the mass spectrometer.

- Plumb the valve for diverter mode.

Figure 4-8 Diverter Valve—Diverter Mode Position A


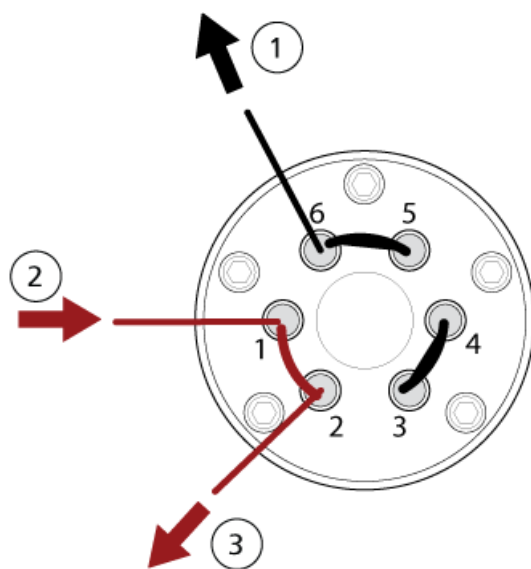


Figure 4-9 Diverter Valve—Diverter Mode Position B

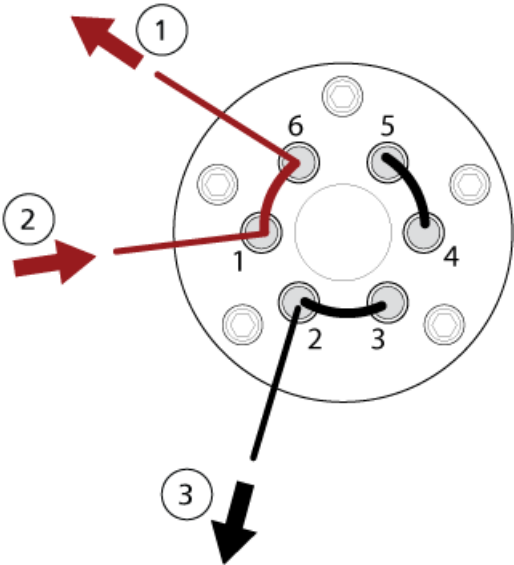

Item	Description
1	To mass spectrometer
2	From column
3	Waste out

Operating Instructions — Hardware Profiles and Projects

5

Hardware Profiles

A hardware profile tells the software how the mass spectrometer and the devices are configured and connected to the computer. Multiple hardware profiles can be set up, but only one profile can be active at any time.

When a hardware profile is created in the Hardware Configuration Editor, the peripheral devices must be configured so that the software can communicate with them. Configuring the peripheral devices requires two procedures: setting up the physical connections and configuring the software to communicate with the peripheral devices. When the software is installed, the driver required for each peripheral device is also installed. After the peripheral devices are physically connected to the computer, set up the appropriate configuration information.

Each hardware profile must include a mass spectrometer. Before creating an acquisition method, make sure that all devices to be used in the method are included in the hardware profile, including the syringe pump. The devices configured in the active hardware profile and selected in the Add/Remove Device Method dialog are shown as icons in the Acquisition method pane. Only peripheral devices included in the active hardware profile can be used when creating acquisition methods.

Each hardware profile must include a mass spectrometer. Before creating an acquisition method, make sure that all devices to be used in the method are included in the hardware profile. The devices configured in the active hardware profile and selected in the Add/Remove Device Method dialog are shown as icons in the Acquisition method pane. Only peripheral devices included in the active hardware profile can be used when creating acquisition methods.


For information about setting up the physical connections to the devices, refer to the *Peripheral Devices Setup Guide*. For a list of the supported devices, refer to the *Software Installation Guide* for the Analyst® software.

Create a Hardware Profile

The user can create multiple hardware profiles, but only one profile can be active at any time.


1. On the Navigation bar, under **Configure**, double-click **Hardware Configuration**.

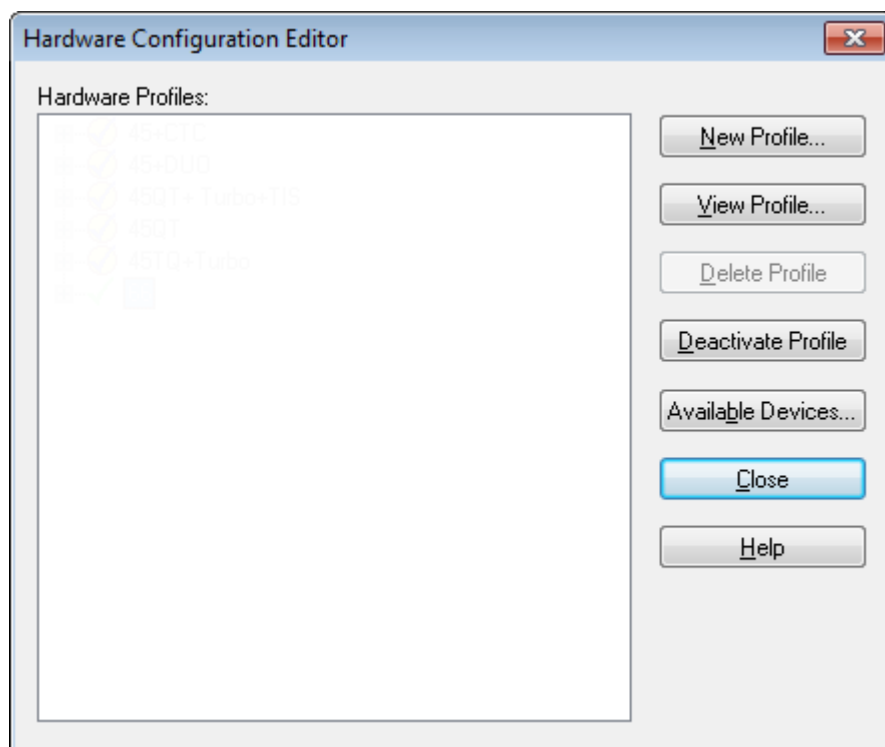Figure 5-1 Hardware Configuration Editor Dialog




2. In the Hardware Configuration Editor dialog, click **New Profile**.

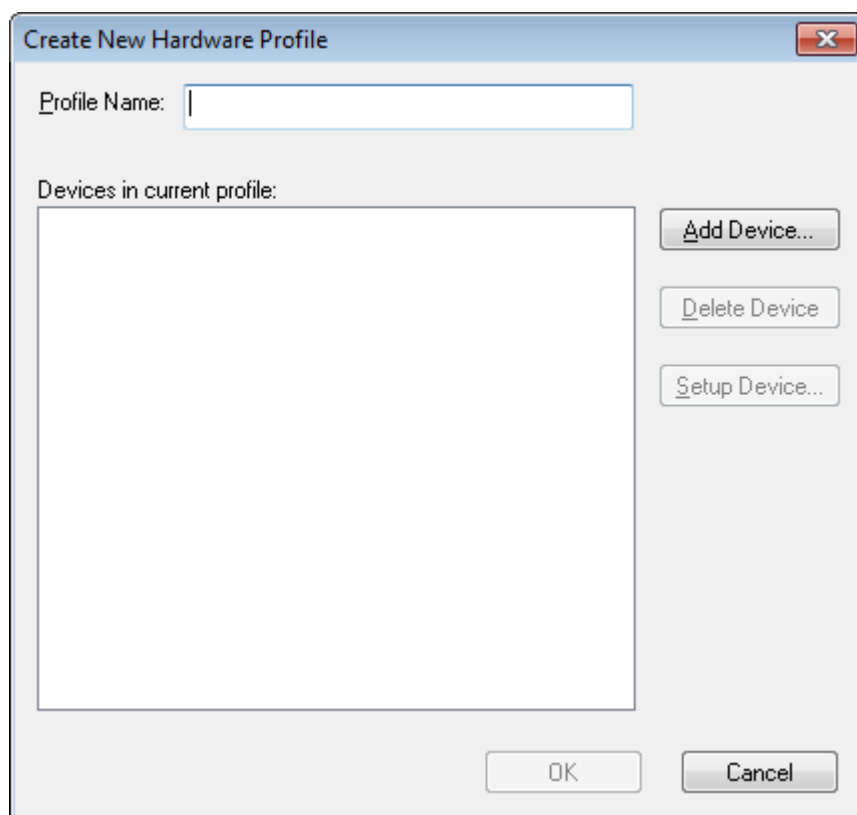Figure 5-2 Create New Hardware Profile Dialog



3. Type a name in the **Profile Name** field.
4. Click **Add Device**.

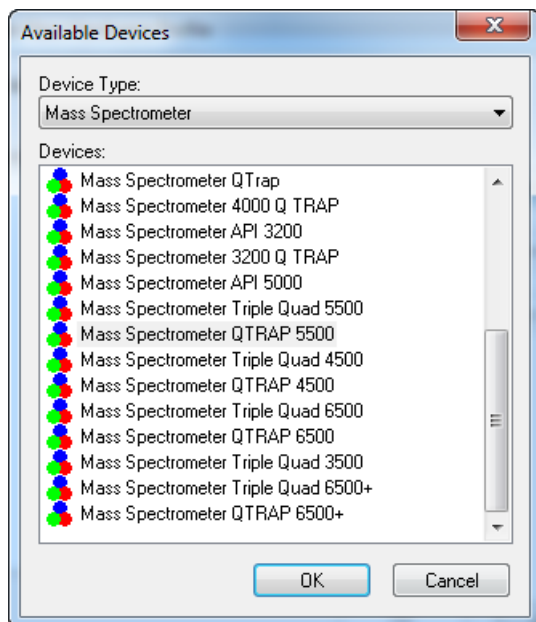Figure 5-3 Available Devices Dialog



In the Available Devices dialog, in the **Device Type** field, **Mass Spectrometer** is the preset value.

5. In the **Devices** list, select the appropriate mass spectrometer and then click **OK** to return to the Create New Hardware Profile dialog.
6. Click **Setup Device**.
7. (Optional) To configure mass spectrometers that use the integrated syringe pump, on the **Configuration** tab, select the **Use integrated syringe pump** check box.
8. (Optional) To configure the mass spectrometer for the diverter valve, on the **Configuration** tab, select **Use integrated injector/diverter valve**.
9. (Optional) Select additional features on the **Configuration** and **Communication** tabs as required.
10. Click **OK** to return to the Create New Hardware Profile dialog.
11. Click **Add Device** and then add and set up each device that is used with the mass spectrometer. Refer to [Add Devices to a Hardware Profile on page 58](#).
12. Click **OK** in the **Create New Hardware Profile** dialog.
13. Click the hardware profile to be activated in the **Hardware Configuration Editor**.
14. Click **Activate Profile**.

The check mark turns green. If a red x is shown, then there is an issue with the hardware profile activation.

Tip! A hardware profile need not be deactivated before another is activated. Click a hardware profile and then click **Activate Profile**. The other profile is deactivated automatically.

15. Click **Close**.

Add Devices to a Hardware Profile

Devices must be configured to enable the software to communicate with them. When the software is installed, the driver required for each device is also installed. After the devices are physically connected to the computer, configure them.


Only the devices configured in the active hardware profile and selected in the Add/Remove Device Method dialog are shown as icons in the Acquisition Method Browser pane.

1. Open the Hardware Configuration Editor.
2. In the **Hardware Profiles** list, deactivate the hardware profile.
3. Click **Edit Profile**.
4. Click **Add Device**.

Note: Remember to add a mass spectrometer. For more information, refer to [Create a Hardware Profile on page 54](#).

5. In the Available Devices dialog, in the **Device Type** list, select the device and then click **OK**.

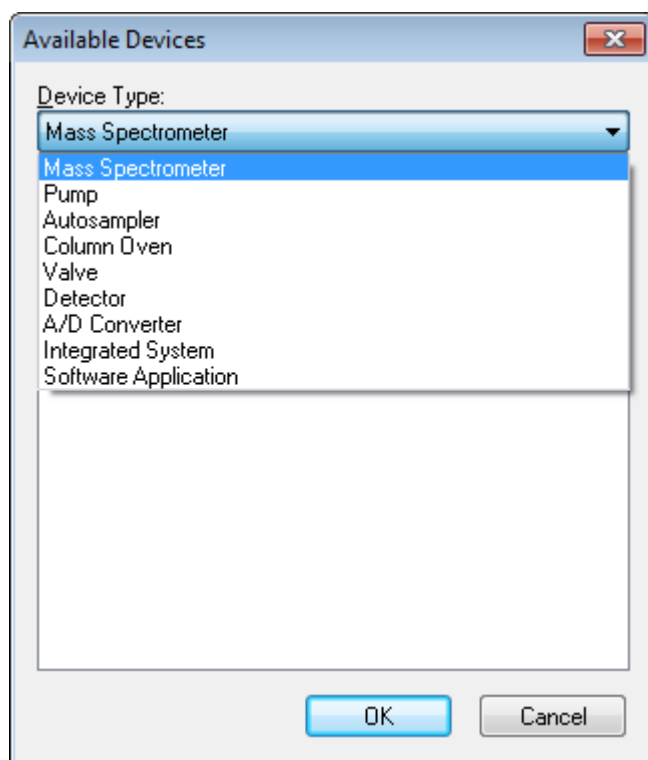Figure 5-4 Available Devices Dialog



6. Click **OK**.
7. Select the device from the **Devices** list and then click **OK**.
8. Click **Setup Device**.

A dialog containing configuration values for the device opens.

9. (Optional) On the Communication tab, in the **Alias** field, type a name or other identifier for the device.

Note: For devices using serial communication, make sure that the serial port selected matches the serial port to which the device is physically connected. When the serial expansion cable is used, the number selected in the profile is the number on the cable plus two.

Note: The **Alias** field might also be referred to as the **Name** box and might be found on another tab under **Alias**.

- If the device uses a **Serial Port** as a communication interface, then in the **COM Port Number** list, select the COM port to which the device is connected.

- If the device uses **Ethernet** as a communication interface, then type the **IP Address** assigned to the device by the administrator or use the corresponding **Host Name** for the address.
- If the device uses **GPIB Board** as a communication interface, then do not change the settings for the GPIB board.

The rest of the preset values for the device are probably appropriate. Do not change them. For information about the Configuration and Communication tabs, refer to the Help.

10. To restore the device preset values, on the Communication tab, click **Set Defaults**.
11. To save the configuration, click **OK**.
12. Repeat step 4 to step 11 for each device.
13. Click **OK** in the Create New Hardware Profile dialog.
14. To activate the hardware profile, in the Hardware Configuration Editor, click the hardware profile.
15. Click **Activate Profile**.

The check mark turns green. If a red x is shown, then there is an issue with the hardware profile activation. For more information, refer to [Troubleshoot Hardware Profile Activation on page 60](#).

Tip! An active hardware profile does not have to be deactivated before another one is activated. Click an inactive hardware profile and then click **Activate Profile**. The other profile is deactivated automatically.

16. Click **Close**.

Troubleshoot Hardware Profile Activation

If a hardware profile fails to become active, then a dialog opens indicating which device in the profile failed. A device might fail to activate because of communications errors.

1. Read the error message generated. Depending on the message, there might be an issue with a device or how the communication is set up.
2. Verify that the device has power and is turned on.
3. Verify that the COM port or IP address assigned to the device is correct.

Tip! On computers with two built-in serial ports, the first port on the serial port expansion card is usually COM3, even though the cable indicates P1.

4. Verify that the communication settings for the device (for example, dip switch settings) are set correctly and match the settings on the Communication tab.
5. Turn off the device.
6. Wait 10 seconds.

7. Turn on the device.

Wait until all device power-up activities are complete before trying to activate the hardware profile again. Some devices might require 30 seconds or more to complete the power-up activities.

8. Activate the hardware profile.
9. If the issue persists, then delete the failing profile and create a new one.
10. If the issue still persists, then contact technical support.

Projects and Subprojects


Before beginning an experiment, decide where to store the files related to the experiment. Use projects and subprojects for each experiment to better manage data and compare results. For example, use subprojects to store the results for specific dates.

Create Projects and Subprojects

To use a subproject structure within a project, create the subproject structure when the project is created.

1. Click **Tools > Project > Create Project**.

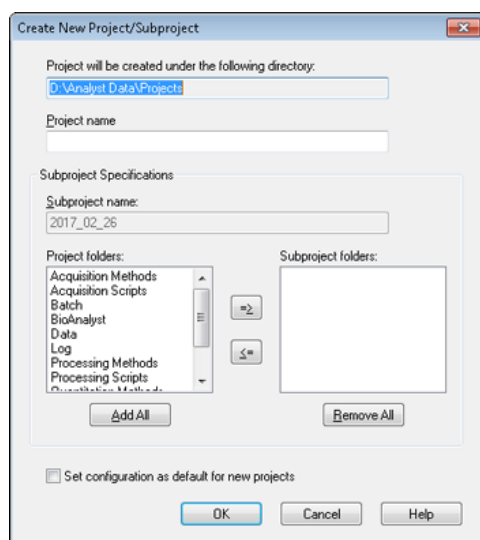Figure 5-5 Create New Project/Subproject Dialog



Note: A new subproject cannot be created for a project that was not originally created with a subproject.

2. Type a project name in the **Project name** field.

3. (Optional) To use subprojects, select the required folders and then use the arrow buttons to move them to the **Subproject folders** list.
4. (If subprojects are used) In the **Subproject name** field, type a name for the first subproject or use the existing date.
5. (Optional) To use this project and subproject folder organization for all new projects, select the **Set configuration as default for new projects** check box.

All new projects are created with this folder configuration.

6. Click **OK**.

Create Subprojects

Subprojects can only be created in a project that has an existing subproject structure.

1. On the **Project** tool bar, from the **Project** list, select the project.
2. Click **Tools > Project > Create Subproject**.
3. In the **Subproject name** box, type a name for the subproject or use the existing date.
4. Click **OK**.

Copy Subprojects

A subproject can be copied from another project that has existing subprojects. If the copied subprojects contain folders that also exist in the project folder, then the software uses the project level folders.

1. Click **Tools > Project > Copy Subproject**.

The Copy Subproject dialog is shown.

2. Click **Browse** to navigate to the subproject source.
3. Click **OK**.
4. Select the subproject from the **Source Subproject** list.
5. Click **Browse** to navigate to the subproject destination.
6. Type the name in the **Target Subproject** field.
7. Click **OK**.
8. Do one of the following:
 - To copy all folders and files from the **Subproject Source** into the **Subproject Destination**, select the **Copy Contents** check box.


- To copy only the folders in the same structure into the **Subproject Destination**, make sure that the **Copy Contents** check box is cleared.

9. Click **Copy**.

Switch Between Projects and Subprojects

- On the software tool bar, from the project list, click the required project or subproject.

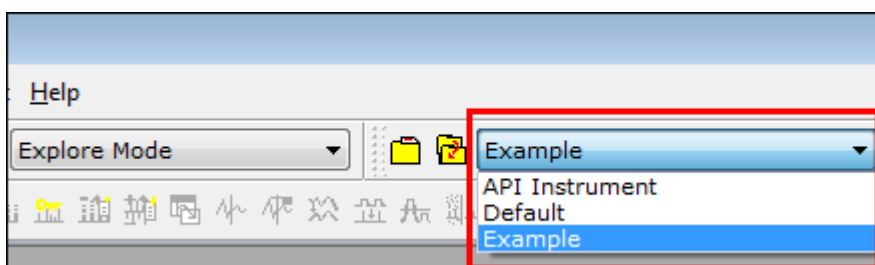Figure 5-6 Project List



The project list in this figure shows the **API Instrument**, **Default**, and **Example** folders.

Installed Project Folders

Three project folders are installed with the software: **API Instrument**, **Default**, and **Example**.

API Instrument Folder

The API Instrument folder is unique and very important to the correct functioning of the mass spectrometer. The API Instrument folder contains the information required for tuning and calibrating the mass spectrometer. This information includes parameter settings files, reference files, instrument data files that contain calibration and resolution information, and the acquisition methods used during automatic tuning. The API Instrument folder also contains data files for manual tuning runs that were performed using the Start button rather than the Acquire button. These data files are saved automatically in the API Instrument folder in the Tuning Cache folder and named with the date and time they were created. The Tuning Cache folder is automatically cleared periodically.

Default Folder

The Default folder contains folders that are present in new projects and serves as a template for new projects.

Example Folder

The Example folder contains sample methods and data files. Users can practice working with the Explore or Quantitate modes using the example data files. The example files are sorted into subfolders by mass spectrometer type and application area.

Back Up the API Instrument Folder

Back up the API Instrument folder regularly and after routine maintenance has been performed.

- Copy the API Instrument folder, paste it to a different location, preferably to another computer, and then rename the folder. Use the date and a mass spectrometer reference if there is more than one mass spectrometer when the folder is named. For example, API Instrument_*instrument model*3_010107

Recover the API Instrument Folder

Back up the **API Instrument** folder regularly and after routine maintenance has been performed. To recover the API Instrument folder, do the following:

1. Rename the current **API Instrument** folder.
2. Copy the backup folder into the **Projects** folder.
3. Change the name of the backup folder to **API Instrument**.

Tune and Calibrate

6

Run the **Verify instrument performance** option weekly or after the mass spectrometer is cleaned to confirm that the system is working properly. In general, for triple quadrupole systems, calibration and resolution is maintained for three to six months unless the system loses vacuum. For QTRAP® systems, the resolution should also be maintained for three to six months but the system should be calibrated approximately monthly. If the system loses vacuum, then check the calibration and resolution before using the system. For more information about tuning and calibration, refer to the *Advanced User Guide* and the *Manual Tuning Tutorial*.

Tip! Perform maintenance tasks regularly to make sure that the mass spectrometer is performing optimally.

Prerequisites

- The spray is stable and the correct tuning solution is being used.
- A printer is configured.

Required Materials

- Tuning solutions that are supplied in the Standards Chemical Kit shipped with the system. If required, a new kit can be ordered from SCIEX.
- 5 mL, 1 mL, and 250 µL serial gas-tight syringes.
- Red PEEK sample tubing.

Verify Instrument Performance

The following procedure describes how to verify or adjust the performance of the mass spectrometer. For information about using the other instrument performance options, refer to the Help. Refer to [Calibration Ions and Solutions](#).

1. On the Navigation bar, under **Tune and Calibrate**, double-click **Manual Tuning**.
2. Start the syringe pump and then run a calibration method. Confirm that the Total Ion Chromatogram (TIC) is stable and that the peaks of interest are present in the spectrum.
3. On the Navigation bar, under **Tune and Calibrate**, double-click **Instrument Optimization**.

The **Instrument Optimization** dialog opens.

4. Click **Verify instrument performance**.
5. Click **Next**.

6. Click **Approved Tuning**.
7. Click **Next**.
8. Select a **Tuning Solution** from the list.

Depending on the solution selected, different modes are available:

- a. Click a polarity.
 - b. (If available) Click **Q1** and **Q3** in the **Quad** section.
 - c. (If available) Click the required scan speeds.
 - d. (If available) Click the scan speeds in the **LIT** section.
 - e. (If available) Click **Excitation** in the **MS/MS/MS** section
9. Click **Next**.
 10. If the **Select a mode** page opens, then select **Automatic**.
 11. Click **Next**.
 12. Click **GO**.

The Verifying or Adjusting Performance dialog opens. After the process has completed, the Results Summary opens. For more information, refer to the Help.

13. If applicable depending on the options selected. Change solutions for different scan types and polarities when prompted.

About the Verifying or Adjusting Performance Dialog

The top left corner shows the part of the instrument that is being tuned.


The Current Spectrum graph shows the spectrum of the current scan, the optimal scan selected by the software, or the scan at the current parameter value when the software results are viewed in interactive mode.

The Instrument Optimization Decision Plots, in the top right graph, dynamically show the intensity versus voltage curves of the parameters that are currently being optimized.

Results Summary

The Results Summary is a record of any instrument settings changes that were made by the Instrument Optimization wizard.

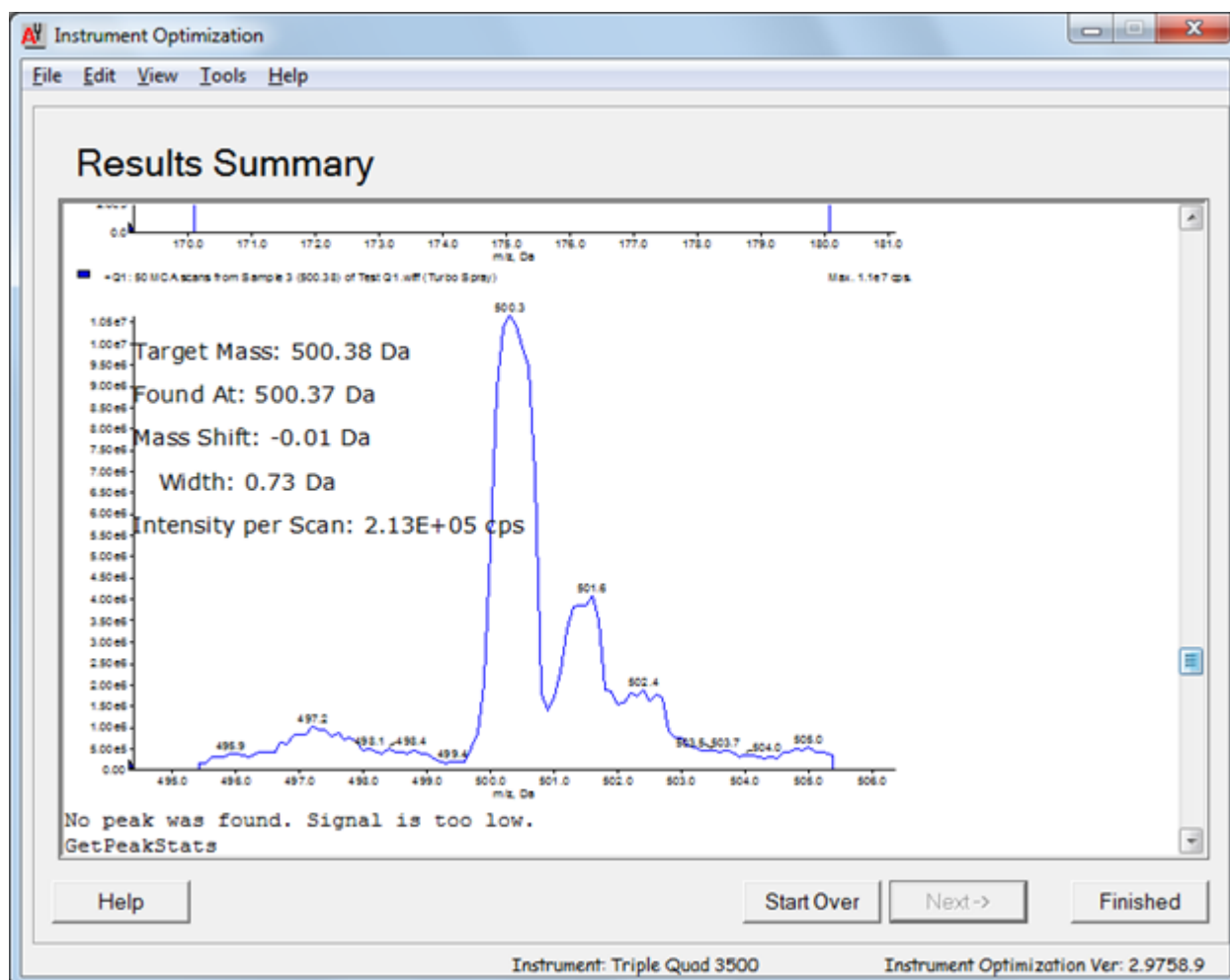Figure 6-1 Results Summary



The Results Summary includes the location of data files and instrument settings backups, as well as step-by-step changes and results during optimization.

The Results Summary also shows a verification report. This report contains a snapshot of the mass spectrum for each relevant mass for the scan modes being verified. The spectrum is labeled with the target mass, where the mass was found, mass shift, peak width, and peak intensity. The spectrum can be used as a visual record of peak shape or scan mode performance. A summary table of results follows the spectra.

The Results Summary is automatically saved in the following path: <drive>\Analyst Data\Projects\API Instrument\Data\Instrument Optimization\yyyy-mm-dd\results.doc, where *yyyy-mm-dd* is the date on which the report was created. Users can print the Results Summary or open a previously saved Results Summary.

Automatic Optimization

7

Users will learn how to:

- Automatically optimize for the analyte using the Compound Optimization wizard.
- Choose between infusion and flow injection analysis (FIA).
- Optimize parameters:
 - If infusion analysis is chosen, then use an infusion method to optimize compound-dependent parameters.
 - If FIA analysis is chosen, then use FIA to optimize compound- and ion source-dependent parameters.

This tutorial uses minoxidil, tolbutamide, reserpine, and rescinnamine as example compounds. Other available compounds can be used, but the methods must be adjusted accordingly.

The user can also optimize compounds manually.

Prerequisites

- The mass spectrometer is tuned and calibrated.
- (For FIA analysis) An acquisition method template is available.
- (For FIA analysis) An LC pump and an autosampler are connected and configured in the hardware profile.
- If the system has an integrated syringe pump, a syringe pump is configured in the hardware profile.
- All of the required peripheral devices, including the LC components, if needed, are configured in the hardware profile.

Required Materials

- Four-compound mix (50 ng/mL), consisting of reserpine, minoxidil, tolbutamide, and rescinnamine. The solution can be used for infusion and FIA. The concentration is system-dependent. Use a solution that is 49.9% acetonitrile, with 50% deionized water and 0.1% formic acid as diluent.
- A syringe, preferably a 1.0 mL syringe.
- (For FIA analysis) Mobile phase: 1:1 acetonitrile:water + 2 mM ammonium acetate + 0.1% formic acid.

Note: Users can choose a different mobile phase based on the experimental properties of the compound.

- LC pump and autosampler.
- (For FIA Analysis) Autosampler vials.

Table 7-1 Compounds and Molecular Weights

Compound	<i>m/z</i>
Minoxidil	210.2
Tolbutamide	271.3
Reserpine	609.4
Rescinnamine	635.4

Table 7-2 Estimates for Starting Concentrations

System	Concentration
SCIEX Triple Quad™ 3500 system	10 ng/mL

About Automatic Optimization

Automatic optimization first checks for the presence of the compounds. The voltages of the various ion path parameters are gradually increased or decreased to determine the maximum signal intensity (Q1 scan) for each ion. A text file is generated and shown during the optimization process. This file records the various experiments performed and the optimal values for each parameter. A file folder containing all the experiments performed is also generated and can be found by opening the data file folder in Explore mode. For each experiment performed, an acquisition method is also generated and saved in the Acquisition Method folder.

During the optimization process, select how the precursor ion and the corresponding product ions are to be chosen.

Types of Sample Introduction

Infusion

Infusion is the continuous flow of the sample at low flow rates into the ion source using a syringe pump. During the infusion optimization process, the software can select precursor and product ions and optimize for declustering potential, collision energy, and collision cell exit potential. The voltages of the ion path parameters are gradually increased or decreased to determine the maximum signal intensity for the precursor and product ions.

Use infusion optimization to optimize compound-dependent parameters only at much lower flow rates than those used during LC-MS/MS analysis.

FIA

FIA is the injection of a sample by the autosampler into the mass spectrometer using LC. During the FIA optimization process, multiple sample injections are performed for various ion source- or compound-dependent parameter types

Automatic Optimization

that are changed between injections. FIA compound optimization optimizes parameters by performing looped experiments in succession. One compound-dependent parameter is optimized first, followed by the next compound-dependent parameter. FIA optimizes for ion source-dependent parameters by making one injection for each value.

Compound parameters must be narrowed down using at least two more FIA cycles. Use FIA optimization to optimize both compound-dependent and source-dependent parameters using LC at higher flow rates.

Table 7-3 Differences Between Sample Introduction Methods

Method	Required Devices	Parameters	Typical Flow Rate Range
Infusion	Syringe pump	Compound-dependent	5 µL/min to 25 µL/min
FIA	LC pump and autosampler	Source- and compound-dependent	25 µL/min to 1000 µL/min

During optimization, a text file is generated and then shown. This file records the various experiments performed and the optimal values for each parameter. A file folder containing all the data files is also generated. For each experiment performed, an acquisition method is also generated and saved in the Acquisition Methods folder.

Automatically Optimize for an Analyte Using Infusion

In this section, users will perform automatic MS/MS optimization using infusion with a known precursor ion and an unknown product ion.

Confirm the Presence of Compounds

Confirm the presence of compounds of interest before continuing with automatic optimization.

1. In the Analyst[®] software, create a project.
2. Activate the hardware profile.
3. Infuse the compound in solution at a rate of 5 µL/min. to 10 µL/min.
4. On the Navigation bar, under **Tune and Calibrate**, double-click **Manual Tuning**.
5. In the method list field, click **Syringe Pump Method**.
6. On the Syringe Pump Method Properties tab, type the appropriate parameter values. Refer to [Table 7-4](#).

Table 7-4 Syringe Pump Method Properties Tab

Parameter	Typical Value
Syringe Diameter	Syringe dependent; 1.0 mL syringe is 4.610 mm
Flow Rate	10
Unit	μL/min.

Figure 7-1 Syringe Pump Method Properties Tab

Integrated Harvard Syringe Pump Method Properties

Syringe Diameter (mm): 4.610


Flow Rate: 10.000 Unit: uL/min

7. Click **Start Syringe Pump**.
8. Click **MS Method** from the method list.
9. Click **Start**.
10. Wait until an even TIC is shown on the left and peaks are shown on the right and then click **Stop**.
11. Select the **MCA** check box.
12. Type **10** in the **Cycles** field.
13. Click **Start**.

When the ten scans are complete, the graph should show the masses of the four compounds as ions.

Note: The intensities of the compounds should be much higher than the smallest noise peaks but not so high that any noise peaks are not seen. In the first case, the peak might not be a real compound. In the second case, the concentration might be too high for the software to optimize properly.

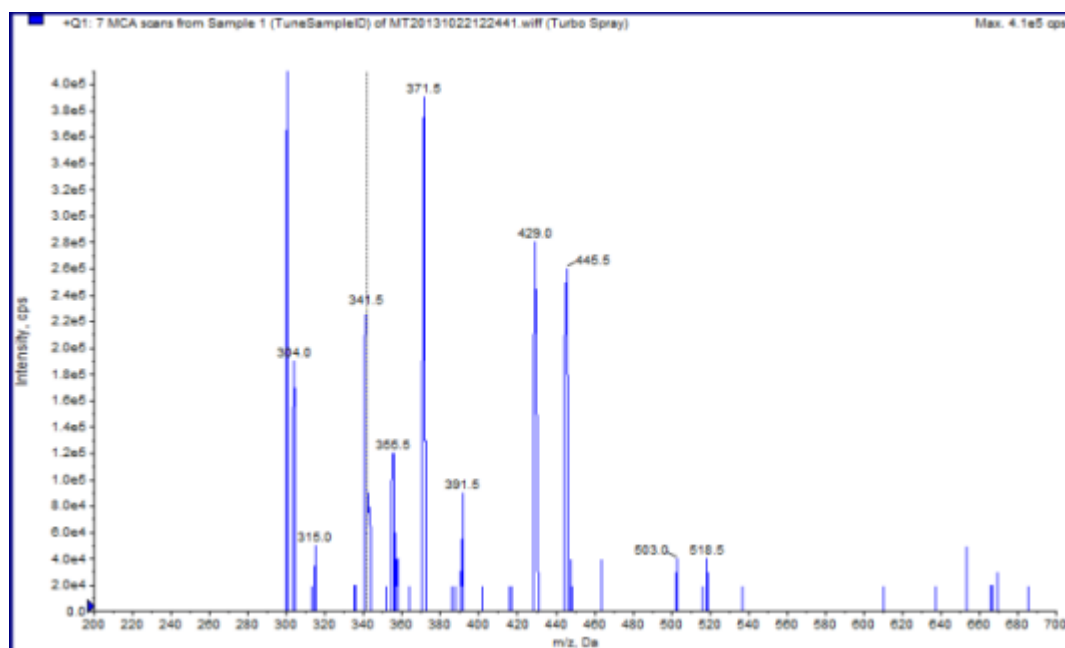Figure 7-2 Compound Ions



Perform Automatic MS and MS/MS Optimization Using Infusion with a Known Precursor Ion and an Unknown Product Ion

Automatic optimization for MS/MS analysis optimizes certain compound-dependent parameters for one or more MRM transitions. The software finds the ion of interest and then optimizes the compound-dependent parameters to get the maximum sensitivity for the compound. The software ramps CE and then selects the most intense fragments that meet all of the product ion selection criteria.

If the initial Q1 scan signal is too high, then the software attempts to reduce the CEM to keep ions within the detector range. If the signal is still too high after the CEM is reduced, then the process stops and an error message is shown. Dilute the solution and then restart the optimization. Be sure to purge the infusion line. The parameters from the last quantitative optimization are stored.

1. Make sure that the proper concentration of solution has been installed in the syringe pump, and that the syringe pump has been started. If an integrated syringe pump is available and has been started, then the Syringe Pump Status LED is flashing.

The syringe pump must be started in the Manual Tune window before compound optimization begins.

2. On the Navigation bar, under **Tune and Calibrate**, double-click **Compound Optimization**.
3. On the Instrument Settings page, in the **Inlet** section, click **Infusion**.
4. Click **MS/MS Analysis** in the **Mass Spectrometer** section.

5. Click **Next**.
6. On the Ions to use in MS/MS Analysis page, select the appropriate parameter values. Refer to [Table 7-5](#).

Table 7-5 Example Parameters to use in MS/MS Analysis Page

Parameter	Value
MW Ion: Search Window	2.500
Resolution	Unit
Polarity	Positive
Product Ion	Auto Select
Resolution	Unit

Note: The optimization algorithm looks for the most intense peak in the specified search window. If the most intense peak in that window is not the mass of interest, then the software optimizes on the wrong ion.

7. Click **Criteria** next to the **Auto Select** option.
8. In the Product Ion Auto Selection Criteria dialog, type the appropriate parameters. Refer to [Table 7-6](#).

Table 7-6 Example Product Ion Auto Selection Criteria Dialog Parameters

Parameter	Parameter	Value
From the Most Intense (peaks)	3	The number of fragment peaks to be optimized. The algorithm generates a product ion scan spectrum while ramping the CE in MCA mode. In this example, the algorithm takes the three most intense fragment ions from the spectrum and then continues the MS/MS optimization on those fragments only. For unknown compounds, select more peaks.
Build final method using (most intense peaks)	2	The number of fragment ions per precursor ion (target compound) to be automatically included in the acquisition method. The specified number defines the number of MRM transitions to be included for each target compound in the method. The order of preference is based on the intensity of the fragment ion. Two is a better starting value than one because typically two product ions are required for quantitation. Start with three in case there is a issue with one of the two best. Go back and the third is already identified. For unknown compounds, select more peaks for use in case of interference.
Exclude Product Ions within \pm (Da of Precursor Ion m/z)	20.000	The Da value that defines the exclusion window around the precursor ion so that fragment ions that fall within this window are not selected for MRM optimization. For example, if the user types ± 5 Da for a precursor ion that has an m/z of 500, then any fragment ions within an m/z of 495 to 505 are excluded. This prevents the precursor ion from being optimized as a product ion.
Min. Mass for Product Ion (Da)	60.000	The lowest fragment mass to be considered for optimization. Use this option to narrow or widen the window of fragment ions to be considered from the precursor mass.
Threshold for Product Ion (cps)	100.000	Minimum number of counts for a product ion to be considered.

9. Click **OK** to save the changes to the selection criteria.
10. Click **Next**.
11. In the Target Components dialog, type the appropriate parameter values. Refer to [Table 7-7](#).

Note: The compound name must be unique for each compound or transition.

Table 7-7 Example Target Compounds Dialog Parameters

Target Compound	Field	Value
Reserpine	Compound Name	Reserpine
	MW (Da)*	609.3
	No. Charges	1
Minoxidil	Compound Name	Minoxidil
	MW (Da)*	210.2
	No. Charges	1
Tolbutamide	Compound Name	Tolbutamide
	MW (Da)*	271.3
	No. Charges	1
Rescinnamine (IS)	Compound Name	Rescinnamine
	MW (Da)*	635.4
	No. Charges	1
*Type the exact ion mass.		

- Click **Finish** to begin the optimization process.

The screen shows two windows, a text file window and an acquisition window. The user might need to minimize one of them to see the other. The experiment being run is shown on the top of the acquisition window. The X-axis shows the parameter that is being optimized for each experiment. The text file window is updated as results are generated.

After optimization is complete, an MRM acquisition file is created and named [compound]_QOpt_FinalMRM_Pos.dam, where [compound] is the first compound in the Target Components page.

Review the Optimization Results

At the end of the optimization, the optimized parameters are saved in an acquisition method. All dam and wiff files generated in the optimization process are saved in the Acquisition Methods folder and in a subfolder in the Data folder, respectively, in the project. The name of the subfolder is generated using the name of the compound and the date.

- After completing the optimization, print the text file containing the optimized parameters for each compound.
- Click **File > Open** and then select the Reserpine_QOpt_FinalMRM.POS.dam file.

Automatic Optimization

3. Compare the values in the text file to those in the dam file.
4. Inspect the contents of the following folders:
 - **Data:** Look through all of the runs executed during optimization. Compare a wiff file with the optimized value in the method or printed parameters.
 - **Acquisition Method:** Reserpine_QOpt_FinalMRM.POS.dam file and other dam files created during optimization.
 - **Log:** Report file (rtf) shown during the optimization process.

Automatically Optimize for an Analyte Using FIA

Prerequisites
<ul style="list-style-type: none">• Identify the ions for the compounds and save the basic acquisition method.• Add an autosampler and LC pump to the basic acquisition method. Using FIA for optimization requires that these devices be active in the hardware profile.• Create an LC-MS/MS acquisition method based on the Reserpine_QOpt_FinalMRM.POS.dam file and then name the new method FIA_AutoOpt_Tutorial.dam.• Make sure that the selected project contains the acquisition method that is being used.

Note: Although FIA can be used to optimize compound-dependent parameters, this is typically not done because of the number of cycles required to obtain the optimal parameter values.

1. Put a dilution of the four-compound mix in an autosampler.

Enough sample is required to review each variable of each parameter and have sample left over. For example, for the temperature to run at 300 °C, 400 °C, and 500 °C, if the injection volume is 10 µL, then more than a 30 µL (3 × 10 µL injection) is required.

2. Confirm that **LC Sync** is selected in the method.

<p>Note: In LC Sync mode, the mass spectrometer coordinates with the action of the LC system to make sure that data is acquired properly.</p>
--

Figure 7-3 Acquisition Method with LC Sync Selected

Acquisition Method Properties

Comment:

Duration (min):

Synchronization Mode: **LC Sync** ▼

Auto-Equilibration

☐ Auto-Equilibration

Auto-Equilibration Duration (min):

Original Configuration

Instrument signature: Triple Quad 3500

Ion Source: Turbo Spray

Device methods:

Agilent 1290 G4226A Autosampler Me

Agilent 1290 LC Binary SL Pump Meth

3. Make sure that the ion source and gas parameters are set to reasonable levels to prevent the contamination of the mass spectrometer while optimizing. For more information, refer to the ion source *Operator Guide*.
4. Set the horizontal micrometer to 5.
5. Set the vertical micrometer on the ion source for the flow rate. As a starting point, use the parameters in the following table: For more information, refer to the ion source *Operator Guide*.

Table 7-8 Ion Source Vertical Parameters

Flow rate	Initial Vertical Parameters
1 µL/min to 20 µL/min	10 mm
20 µL/min to 250 µL/min	5 mm
250 µL/min to 500 µL/min	2 mm
500 + µL/min	0 mm

6. Set the values for the HPLC system and use an autosampler injection volume of 10 µL. Use the same concentration or lower as for the infusion experiment.

The LC pumps must be configured for an isocratic run with no column. The MS and LC times must be the same to collect the proper data.

Automatic Optimization

The flow rate and percent of mobile phases used should be based on the LC column used, the general chromatography, and the approximate mobile phase concentration at which the compounds of interest elute.

7. On the Navigation bar, under **Tune and Calibrate**, double-click **Compound Optimization**.
8. On the Instrument Settings page, depending on the LC system being used, type the appropriate parameter values. Refer to [Table 7-9](#).

Table 7-9 Example Instrument Settings Parameters

Parameter	Value
Inlet	FIA
Default Acq. Method	FIA_AutoOpt_Tutorial.dam
Rack Code	Autosampler specific
Rack Position	Autosampler specific
Injection Volume	10 µL
Mass Spectrometer	MS/MS Analysis

9. Click **Next**.
10. Make sure that the **Int. Std.** check box is cleared.

Selecting the check box indicates which MRM corresponds to the internal standards. Internal standards are not optimized during the optimization process.

11. In the **Resolution** section, select **Unit** in both the **Q1 Resolution** and **Q3 Resolution** fields.

Figure 7-4 Q1 and Q3 Resolution Fields

FIA Target Compounds

These target compounds will be optimized. You may change the compound name. Please specify which one of them are used as internal standards.

	Compound Name	Q1 Mass (Da)	Q3 Mass (Da)	Int. Std.	Vial Pos.
1	Compound 609.300-195.00	609.300	195.000	<input type="checkbox"/>	1
2	Compound 210.200-164.20	210.200	164.200	<input type="checkbox"/>	1
3	Compound 271.100-91.100	271.100	91.100	<input type="checkbox"/>	1
4	Compound 635.300-221.20	635.300	221.200	<input type="checkbox"/>	1

Note: All compounds identified as I.S. (internal standard) will not be used to determine optimum Source / Gas Parameter conditions.

Resolution

Q1 Resolution:

Q3 Resolution:

< Back Next > Cancel Help

12. Click **Next**.
13. On the FIA Source Parameters page, type numbers that are lower or higher than the original value as long as they are still within specifications.

Be sure not to go too low with any of the settings to keep the system clean. Refer to [Table 7-10](#) for parameters that can be used as a starting point.

Tip! Type the values before selecting the check box.

Automatic Optimization

Table 7-10 Example Parameters for FIA Source Parameters Page

Parameter	Select the Optimize Check Box?	Values for Optimization
Curtain Gas	Yes	20;40;55
Collision Gas	No	—
IonSpray Voltage	Yes	1500;2000;3000;4000;5000
Temperature	Yes	300;400;500;600;700
Ion Source Gas 1	Yes	40;50;60;70;80;90
Ion Source Gas 2	Yes	40;50;60;70;80;90
Interface Heater	No	—

14. Select **1** or **2** beside **Replicate Injection for each Parameter**.

The total number of injections and the total sample volume are calculated based on the specifications here. Note the total sample volume needed. Sample volume might be high depending on how many variables for each parameter are being optimized, as each variable is a separate method.

Figure 7-5 Example Replicate Injection for each Parameter Field

Please select the Source Parameters to optimize in FIA:

	Parameter Name	Optimize	Current Val.	Values for Optimization
1	Curtain Gas	<input checked="" type="checkbox"/>	10.0	20.0;40.0;55.0;
2	Collision Gas	<input type="checkbox"/>	9.00	
3	IonSpray Voltage	<input checked="" type="checkbox"/>	5500.0	1500.0;2000.0;3000.0;4000
4	Temperature	<input checked="" type="checkbox"/>	0.0	300.0;400.0;500.0;600.0;70
5	Ion Source Gas 1	<input checked="" type="checkbox"/>	15.0	40.0;50.0;60.0;70.0;80.0;90
6	Ion Source Gas 2	<input checked="" type="checkbox"/>	0.0	40.0;50.0;60.0;70.0;80.0;90

Replicate Injection for each Parameter: ☐ 1 ☒ 2 ☐ 3 ☐ 4

Total # of injections: 50

Total Sample Volume: 50 (μl)

< Back Next > Cancel Help

15. Click **Next**.
16. On the FIA Compound Parameters page for each analyte, use the example parameter values provided as a starting point. Refer to [Table 7-11](#).

Note: The values in [Table 7-11](#) are suggested values. For more information, refer to the Help.

Table 7-11 Example FIA Compound Parameters Page

Parameter	Select the Optimize Check Box?	Values for Optimization
Decustering Potential	Yes	60;80;100;120;200
Entrance Potential	No	—

Table 7-11 Example FIA Compound Parameters Page (continued)

Parameter	Select the Optimize Check Box?	Values for Optimization
Collision Energy	Yes	20;30;40;50;70;80;100
Collision Cell Exit Potential	Yes	2;4;6;8;10;12

The total number of injections and dependent sample volume update automatically. In contrast to ion source parameters, which require one injection per value per replicate, compound-dependent parameters only require one injection per parameter. A looped experiment is performed for each parameter. The values are alternated scan-by-scan within one injection.

Note: Do not enter too many values. Doing so will prevent proper evaluation of the parameter.

17. Using the Compound list, move to another compound and type the parameters to be optimized.
18. Repeat step 19 until all of the parameters for all of the compounds are provided.
19. Type **1.5** in the **Mass Spec. Duration** field.

Note: This value should be at least the required length of time for each injection.

Figure 7-6 Example Mass Spec. Duration Field

Please select the Compound Parameters to optimize by FIA:

Compound: MRM: 609.300 - 195.000

	Parameter Name	Optimize	Current Val.	Values for Optimization
1	Declustering Potential	<input checked="" type="checkbox"/>	120.0	60.0;80.0;100.0;120.0;200.
2	Entrance Potential	<input type="checkbox"/>	10.0	
3	Collision Energy	<input checked="" type="checkbox"/>	30.0	20.0;30.0;40.0;50.0;70.0;80
4	Collision Cell Exit Potential	<input checked="" type="checkbox"/>	15.0	2.0;4.0;6.0;8.0;10.0;12.0;

Total # of Injections: 53
 Total Sample Volume: 53 (μl)
 Mass Spec. Duration: (min)

< Back Finish Cancel Help

20. Click **Finish** to begin the optimization process.

The software optimizes the specified ion source- and compound-dependent parameters to get the maximum sensitivity for the MRM transition of the compound. As the software proceeds through the optimization, it creates a **Compound Optimization** report.

21. To obtain optimized parameters, repeat this routine.

Note: Typically, the ion source and gas parameters must be narrowed using one more FIA cycle.

22. Open the final optimized FIA method called *_FIA_sample_1.

Note: The software generates several acquisition methods.

23. Save this method using a simpler name.

Operating Instructions — Acquisition Methods

8

An acquisition method consists of experiments and periods. Use the Acquisition Method Editor to create a sequence of periods and experiments for the instrument and devices.

We recommend that only users who are proficient in method development create or modify acquisition and quantitation methods. Refer to the *Laboratory Director's Guide* for more information about roles and security.

Create an Acquisition Method Using the Acquisition Method Editor

Tip! If users are creating a new acquisition method file from an existing file, then some or all of the peripheral device methods in the acquisition method might be used.

Only devices configured in the active hardware profile appear in the Acquisition method pane. Any devices added to the hardware profile must also be added to existing acquisition methods. For more information about devices, refer to the *Peripheral Devices Setup Guide*.


1. Make sure that a hardware profile containing the mass spectrometer and peripheral devices is active.
2. On the Navigation bar, under **Acquire**, double-click **Build Acquisition Method**.
3. Select a **Synchronization Mode** on the Acquisition Method Properties tab.
4. (Optional) Select the **Auto-Equilibration** check box and then type the required equilibration time, in minutes.
5. Click the **Mass Spec** icon in the Acquisition method pane.
6. Select a **Scan type** on the MS tab.
7. Type values in the fields as required. Refer to [Parameters on page 27](#).
8. On the Advanced MS tab, type values in the fields as required.
9. On the MS tab, click **Edit Parameters**.
10. On the Source/Gas tab, specify values in the fields as required.
11. On the Compound tab, specify values in the fields as required and then click **OK**.
12. Click a device icon and then select the parameters for the device.
13. Add any additional periods and experiments. Refer to [Add an Experiment on page 85](#) and [Add a Period on page 86](#).
14. Click **File > Save**.

Configure the Integrated Syringe Pump

1. Confirm that the integrated syringe pump is selected in the hardware device profile.
2. Click the **Syringe Pump** icon in the Acquisition method pane.

The Syringe Pump Method Properties tab opens in the Acquisition Method editor.

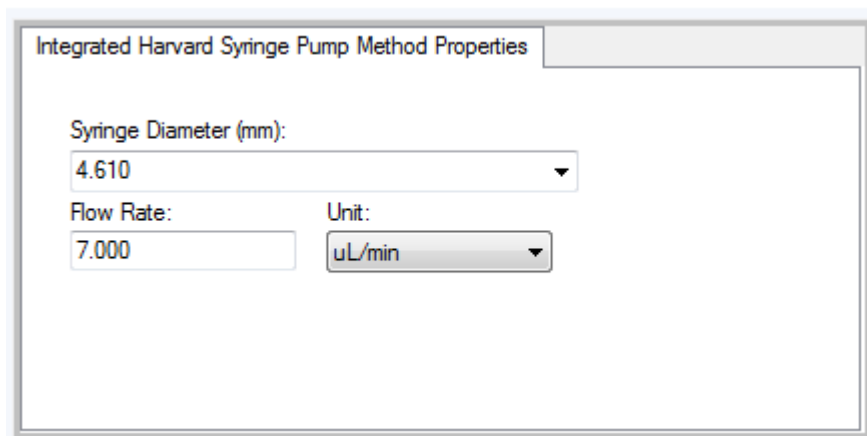Figure 8-1 Syringe Pump Properties tab




The screenshot shows a dialog box titled "Integrated Harvard Syringe Pump Method Properties". It contains three input fields: "Syringe Diameter (mm):" with a dropdown menu showing "4.610", "Flow Rate:" with a text box containing "7.000", and "Unit:" with a dropdown menu showing "uL/min".

3. Type the syringe diameter in the **Syringe Diameter (mm)** field.
4. Type the flow rate in the **Flow Rate** field.
5. Select the units of flow from the **Unit** list.

Add an Experiment

1. In the Acquisition method pane, on the period where the experiment is to be added, right-click and then click **Add experiment**.

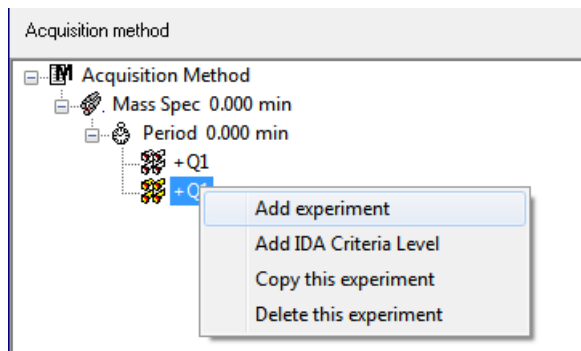Figure 8-2 Add Experiment



An experiment is added below the last experiment in the period.

Note: An experiment cannot be inserted between experiments, IDA criteria, or periods. Users can only add an experiment at the end of the period.

2. In the MS tab, select the appropriate parameters.

Add a Period

- In the Acquisition method pane, right-click the **Mass Spec** icon, and then click **Add period**.

A period is added below the last period created.

Note: Users cannot use multiple periods in an IDA experiment.

Copy an Experiment into a Period

1. Open a multi-period method.
2. In the Acquisition method pane, press **Ctrl**, and then drag the experiment to the period.

The experiment is copied below the last experiment in the period.

Copy an Experiment within a Period

Use this procedure to add the same or similar experiments to a period if most or all of the parameters are the same.

- Right-click the experiment and then click **Copy this experiment**.

A copy of the experiment is added below the last experiment created. This is useful when the same or similar experiments are added to an acquisition method.

Scan Techniques

MS: In MS scans, also referred to as single MS scans, ions are separated according to their mass-to-charge ratio (m/z). A single MS scan might be used to determine the molecular weight of a compound. Single MS scans can also be referred to as survey scans. MS scans do not provide any information about the chemical make-up of the ions other than the mass. Perform MS/MS scans to obtain more information about the ions.

MS/MS: MS/MS scans are used to determine a molecular species.

- For MS/MS scans in triple quadrupole systems, precursor ion fragmentation occurs in the collision cell.

If enough energy is used, then the precursor ion fragments to produce characteristic product ions.

Quadrupole-Mode Scan Types

Triple quadrupole instruments have high-sensitivity Multiple Reaction Monitoring (MRM) capabilities required for quantitation experiments. In addition, they have highly specific scan types, such as precursor ion and neutral loss scans, which allow a more advanced search to be performed on the components of the samples.

Q1 MS (Q1): A full scan type using the first quadrupole (Q1). The ion intensity is returned for every mass in the scan range.

Q1 Multiple Ions (Q1 MI): A zero-width scan type using the Q1 quadrupole. The ion intensity is returned for the specified masses only.

Q3 MS (Q3): A full scan type using the third quadrupole (Q3). The ion intensity is returned for every mass in the scan range.

Q3 Multiple Ions (Q3 MI): A zero-width scan type using the Q3 quadrupole. The ion intensity is returned for the specified masses only.

MRM (MRM): An MS/MS scan in which a user-selected ion is passed through the Q1 quadrupole and then fragmented in the Q2 collision cell. The Q3 quadrupole then selects the fragment ion that enters the detector. This scan mode is used primarily for quantitation.

Product Ion (MS2): An MS/MS full scan where the Q1 quadrupole is fixed to transmit a specific precursor ion and the Q3 quadrupole scans a defined mass range. Used to identify all of the products of a particular precursor ion.

Precursor Ion (Prec): An MS/MS scan where the Q3 quadrupole is fixed at a specified mass-to-charge ratio to transmit a specific product ion and the Q1 quadrupole scans a mass range. Used to confirm the presence of a precursor ion or, more commonly, to identify compounds sharing a common product ion.

Neutral Loss (NL): An MS/MS scan where both the Q1 quadrupole and the Q3 quadrupole scan a mass range, a fixed mass apart. A response is observed if the ion chosen by the Q1 quadrupole fragments by losing the neutral

loss (the fixed mass) specified. Used to confirm the presence of a precursor ion or, more commonly, to identify compounds sharing a common neutral loss.

LIT-Mode Scan Types

The LIT-mode scans use the Q3 quadrupole as a linear ion trap. Ions are trapped and stored in the Q3 quadrupole before being scanned out, giving increased sensitivity. In addition, MS/MS/MS analysis can be performed in the linear ion trap, providing more information about the sample. LIT-mode scan types are typically used for qualitative measurements.

Enhanced MS (EMS): Ions are scanned in the Q1 quadrupole to the linear ion trap where they are collected. These ions are scanned out of the Q3 quadrupole to produce single MS type spectra.

Enhanced Multi-Charge (EMC): This scan type is similar to the EMS scan except that before ions are scanned out of the linear ion trap, there is a delay period during which low-charge state ions (primarily singly-charged ions) are allowed to preferentially escape from the linear ion trap. When the retained ions are scanned out, the multiply-charged ion population dominates the resulting spectrum.

Enhanced Product Ion (EPI): This scan type is used to obtain a high quality MS/MS spectrum on a specific ion. The fragmentation is done in the Q2 collision cell and thus provides the information-rich MS/MS spectrum typical of collisionally activated dissociation (CAD) fragmentation. In this scan mode, the precursor ion to be fragmented is first selected in the Q1 quadrupole with a mass window that is 1 Da to 4 Da wide, filtering out all other ions. The precursor ion is fragmented by CAD gas in the Q2 collision cell. The fragment ions generated are captured in the linear ion trap and then scanned out at one of three scan speeds, depending on the required fragment ion resolution.

For IDA experiments, the **Product Of** field is set to 30 Da by default, and this value should not be changed.

Enhanced Resolution (ER): This scan type is similar to the EMS scan except that a small 30 Da mass around the precursor mass is scanned out of the linear ion trap at the slowest scan rate to produce a narrow window of the best-resolved spectra.

MS/MS/MS (MS3): A precursor ion is selected by the Q1 quadrupole and fragmented with collisionally activated dissociation in the Q2 collision cell. The resulting product ions are all transmitted to the linear ion trap, where a single product ion is then isolated. The isolated ion is further fragmented in the linear ion trap, and the resulting product ions are scanned out of the trap at one of three scan speeds. As with any in-trap Collision Induced Dissociation (CID) technique, there is a low mass cut-off for the second MS/MS step due to the condition that the lowest mass fragment and precursor must be simultaneously stable in the trap. For QTRAP® systems, this results in the loss of ions lower than 28 percent of the mass of the precursor ion during MS3 experiments. This phenomenon is often referred to as the one-third cut-off rule.

About Spectral Data Acquisition

Spectral data can be acquired in one of the modes described in [Table 8-1](#).

Table 8-1 Spectral Data

Mode	Description
Profile	The preset value is 0.1 Da. Profile data is the data generated by the mass spectrometer and corresponds to the intensity recorded at a series of evenly spaced discrete mass values. For example, for a mass range from 100 Da to 200 Da and a step size of 0.1 Da, the mass spectrometer scans 99.95 to 100.05 (recorded as value 100), 100.05 to 101.15 (recorded as value 101)...199.95 to 200.05 (recorded as value 200).
Peak Hopping	The preset value is 1.0 Da. Peak Hopping is a mode of operating a mass spectrometer in which large steps (approximately 1 Da) are made. It has the advantage of speed (fewer data steps are made) but with the loss of peak shape information.
Centroid	The mass spectrometer scans as in Profile mode, but creates a centroid of the data, replacing found peaks with the intensity-weighted center of gravity for each peak. Centroid data has the advantage of significantly reducing file size. The disadvantage is that peak shape information is lost and if data has been collected as a centroid then it cannot be altered. We recommend the use of profile mode and centroiding of the data post-acquisition.

Operating Instructions — Analyze and Process Quantitative Data

9

This section describes how to use the Analyst[®] software to analyze and process quantitative data. Data can also be processed using the MultiQuant[™] software. We suggest that the MultiQuant[™] software is used to quantitate data. Refer to the documentation that comes with the MultiQuant[™] software.

Use the sample files found in the Example folder to learn how to select samples for quantitation, how to select preset queries and create table-specific queries, and how to analyze the acquired data. For more information about the following topics, refer to the *Advanced User Guide*.

- Metric Plots
- Layout of a Results Table

Quantitative Analysis

Quantitative analysis is used to find the concentration of a specific substance in a sample. By analyzing an unknown sample and comparing it to other samples containing the same substance with known concentrations (standards), the software can calculate the concentration of the unknown sample. The process involves creating a calibration curve using the standards and then calculating the concentration for the unknown sample. The calculated concentrations of each sample are then available in a Results Table.

Quantitation Methods

A quantitation method is a set of parameters used to generate peaks in a sample. The quantitation method can include parameters used to locate and integrate peaks, generate standard curves, and calculate unknown concentrations. A previously saved quantitation method can be selected from the Quantitation menu in the batch. The user can create a quantitation method before data acquisition and then apply the method to the quantitative data automatically after the batch is complete. Alternatively, a quantitation method can be created and applied post-acquisition.

Three tools can be used to create a quantitation method: the Quantitation Wizard, the Build Quantitation Method, and Quick Quant.

Quantitation Wizard

With the Quantitation Wizard, a Results Table is generated at the same time as the quantitation method. Alternatively, an existing quantitation method can be used to quantitate different sets of data.

Build Quantitation Method

The Build Quantitation Method does not generate a quantitation Results Table although the method can subsequently be used in the Quantitation Wizard to create a Results Table. The Build Quantitation Method can also be used to change existing quantitation methods. This is the most flexible way of creating a quantitation method. Refer to [Create a Method Using the Quantitation Method Editor on page 91](#).

Quick Quant

Quick Quant is not recommended for quantitation of results.

Quick Quant is part of the Batch Editor. Use Quick Quant to add compound concentrations prior to data acquisition. Because a sample has not been acquired, a representative sample cannot be selected nor can peaks be reviewed. With this process, only the method components are defined.

To use a previously saved quantitation method, select it from the Quantitation menu in the batch. For instructions on creating a batch, refer to [Create and Submit a Batch on page 135](#).

About Results Tables

Results tables summarize the calculated concentration of an analyte in each unknown sample based on the calibration curve. Results tables also include the calibration curves as well as statistics for the results. The user can customize the results tables and view the results tables in layouts.

The data from a results tables can be exported to a txt file for use in other applications, such as Microsoft Excel. The user can also export data in the table or just the data in the visible columns.

Quantitation Methods and Results Tables

For the following procedures, use the sample data that is installed in Example folder. The **Triple Quad folder** contains the data files, Mix_Batch_1 and Mix_Batch_2. These sample files are used to demonstrate the usefulness of metric plots to isolate problematic samples. The ions scanned were reserpine (609.3/195.0), minoxidil (210.2/164.2), tolbutamide (271.1/91.1) and rescinnamine (635.3/221.2), which is the internal standard. Mix_Batch_1 contains no errors in terms of sample preparation, whereas Mix_Batch_2 contains a QC sample where the internal standard was added twice (sample QC2).

Create a Method Using the Quantitation Method Editor

Prerequisites
<ul style="list-style-type: none"> Select the project that contains the data to be quantified. Switch Between Projects and Subprojects on page 63

- Select the project that contains the data to be quantified.
- [Switch Between Projects and Subprojects on page 63](#)

1. Make sure that the Example folder is selected.

2. On the Navigation bar, under **Quantitate**, double-click **Build Quantitation Method**.

The Select Sample dialog opens.

3. Double-click the **Triple Quad** folder in the **Data Files** list.
4. Select **Mix_Batch_2. wiff**.

The samples in the selected data file are shown in the **Samples** list.

Note: If the **Compound ID** field was populated for the samples and internal standards in the acquisition method, then in the **Internal Standards** table, when a value is selected in the **Q1/Q3** field, the **Name** field is automatically populated.

5. Select a sample that provides a detectable signal to select integration parameters that fit the entire batch and then click **OK**.
6. In the **Internal Standards** table, in the **Name** column, select **rescinnamine**. In **Q1/Q3** column, select **635.3/221.2**.
7. In the **Analytes** table, do the following:
 - a. In the **Name** column, select **minoxidol** for the **Q1/Q3** column masses of **210.2/164.188**, **tolbutamide** for **271.3/91.146**, and **reserpine** for **609.4/195.039**.
 - b. In the **Internal Standard** column, from the list, select the **rescinnamine** as internal standard to be associated with each analyte.
 - c. Delete **635.4/221.185** from the **Q1/Q3** column in the **Analytes** table.

Note: If the Compound ID field was populated for the samples and internal standards in the acquisition method, then in the **Analytes** table, the **Name** field and **Q1/Q3** field are populated.

8. Click the **Integration** tab.

The preset integration parameters are suitable for most peaks.

9. If the integration is not suitable, then change the algorithm. Refer to [Manually Integrate Peaks on page 103](#).
10. Click the **Show or Hide Parameters** icon to show the additional integration algorithms.
11. Click the **Calibration** tab.

The preset parameters are suitable for these samples. User can change the fit, weighting, and regression parameter depending on the specific applications.

12. Save the quantitation method.

The new method can be used when a batch is created in the Batch Editor or when the Quantitation Wizard is used to generate a Results Table.

Tip! The quantitation method can only be used in the current project unless it is copied to another project. To do this, click **Tools > Project > Copy Data**. A new project must be created and selected to be available for use.

Create a Results Table Using the Quantitation Wizard

Prerequisites

- Select the project that contains the data to be quantified.
- [Switch Between Projects and Subprojects on page 63](#)

1. On the Navigation bar, under **Quantitate**, double-click **Quantitation Wizard**.

The Create Quantitation Set - Select Samples page opens.

2. Double-click the **Triple Quad** folder in the **Available Data Files** list.
3. Select **Mix_batch_2.wiff**.
4. Click **Add All**.

Note: We recommend that users do not process or report results from any sample for which acquisition was abnormally or unexpectedly terminated.

5. Click **Next**.

The Create Quantitation Set - Select Settings & Query page opens.

6. Click **Select Existing: Query** in the **Default Query** section.
7. Select **Accuracy 15%** from the **Query** list.

Note: To create a query at the same time, refer to [Create a Standard Query \(Optional\) on page 94](#).

Note: It is the responsibility of the user to evaluate and validate the query to be used for specific applications.

8. Click **Next**.

The Create Quantitation Set - Select Method page opens.

9. Click **Choose Existing Method**.
10. Select **PK Data_Mix.qmf** from the **Method** list.

11. Click **Finish**.

The Results Table opens.

Tip! To add or remove samples in the Results Table, click **Tools > Results Table > Add/Remove Samples**.

12. Review the sample type, actual concentration, peak integration, calibration curves, statistics pane, metric plot for the internal standard, and other information related to data quantitation.
13. Save the Results Table.

Note: We recommend that users do not change datafile (wiff) names if a Results Table includes samples from that file.

Tip! Well-formatted reports can be created from a Results Table using the Reporter software. We recommend that the user validate the results if a Reporter template that contains a query is used. Refer to [Reporter Software on page 129](#).

Create a Standard Query (Optional)

Advanced users can create a query and a standard query numerous ways. The following is one example. For more information about creating queries, refer to the Help.

1. On the Navigation bar, under **Quantitate**, double-click **Quantitation Wizard**.
2. Select samples in the Create Quantitation Set - Select Samples page.
3. Click **Next**.
4. In the Select Settings & Query page, in the **Default Query** section, select **Create New Standard Query**.
5. Type a query name.

Figure 9-1 Create Quantitation Set — Select Settings & Query Page

Create Quantitation Set - Select Settings & Query

Please select the settings for the new results table and the default query (if any). Integration Algorithm: IntelliQuan

Settings to Use: Default

Default Query

☐ None

☐ Select Existing:

Query: Accuracy 15% ☐ Execute Query as Standard Query

☒ Create New Standard Query

Name:

< Back Next > Finish Cancel Help

6. Click **Next**.

Figure 9-2 Create Quantitation Set — Create Default Query Page

Please specify the concentrations/sample names and the corresponding allowed accuracy variations (in percent). You can leave any of the "variation" fields empty as desired.

Maximum Allowed Accuracy Variation for QCs (%) Maximum Allowed Accuracy Variation for Standards (%)

Concentration	Max. Variation
4.000000	
40.000000	
400.000000	
4000.000000	
12000.000000	

Concentration	Max. Variation
0.120000	
0.240000	
0.490000	
0.980000	
1.950000	
3.910000	
7.810000	
15.630000	
31.250000	
62.500000	
125.000000	

☐ Query By Name

< Back Next > Finish Cancel Help

7. In the **Maximum Allowed Accuracy Variation for QCs (%)** table in the **Max. Variation** column, type the maximum allowable percent of variation for each QC, for example 5 is $\pm 5\%$, in the same row as the corresponding concentration. If the concentrations were not specified during acquisition, then they are not shown here. In this case, type them in the **Concentration** column.
8. In the **Maximum Allowed Accuracy Variation for Standards (%)** table, in the **Max. Variation** column, type the maximum allowable percent of variation for each standard, for example 10 is $\pm 10\%$, in the same row as the corresponding concentration. If the concentrations were not specified during acquisition, then they are not shown here. Type the concentrations in the **Concentration** column.
9. Click **Next**.

Figure 9-3 Create Quantitation Set — Select Method Page

10. Select or create a method.

11. Click **Finish**.

The query is applied as a standard query. The query results are shown as a Pass or Fail entry in the **Standard Query Status** column of the Results Table.

Tip! To return to the full view, right-click and then click **Full**.

Results Table Right-click Menu

Right-click in the Results Table to access the options shown in [Table 9-1](#).

Table 9-1 Results Table Right-click Menu

Menu	Function
Full	Shows all the columns.
Summary	Shows specific columns.
Analyte	Shows a specific analyte.

Table 9-1 Results Table Right-click Menu (continued)

Menu	Function
Analyte Group	Creates an analyte group.
Sample Type	Shows samples of a specific type or all samples.
Add Formula Column	Adds a formula column. We recommend that the user validate the results if a formula column is used.
Table Settings	Edits or selects a table setting.
Query	Creates or selects a query.
Sort	Creates a sort or sorts by index.
Metric Plot	Creates a metric plot.
Delete Pane	Deletes the active pane.
Fill Down	Copies the same data into the selected cells.
Add Custom Column	Adds a custom column.
Delete Custom Column	Deletes the selected custom column.

Peak Review and Manual Integration of Peaks

Use peak review to survey the peaks that the software has identified and then redefine the peak or the start and end points where required.

After identifying the analytes and internal standards that the software must find, the software searches for the peaks in the samples. When the software identifies a peak, it shows the chromatograms for each analyte and internal standard in the Create Quantitation Method: Define Integration page of the Standard Wizard or on the Integration tab of the Full Method Editor. The user can confirm the peaks that are found or change the quantitation method to better define the peaks. We recommend that users manually review all integration results.

Review Peaks

During peak review, the user might want to view a peak in its entirety or to examine the baseline to find out how well the software found the start and end points of the peak. The automatic zooming feature can be used to do either.

To help the software find a peak, define the exact start and end points of the peak and background manually. These changes will apply only to that individual peak unless the global method is updated.


Note: We recommend that manually integrated results be validated.

Tip! To review an individual peak, right-click on a point on the curve and then click **Show Peak**. The software opens the Peak Review window with the selected peak.

1. Right-click in the **Results Table** and then click **Analyze**.
2. Select an analyte.
3. Click **Tools > Peak Review > Pane**.


The peaks are shown below the **Results Table** with only the peaks listed in the Results Table.

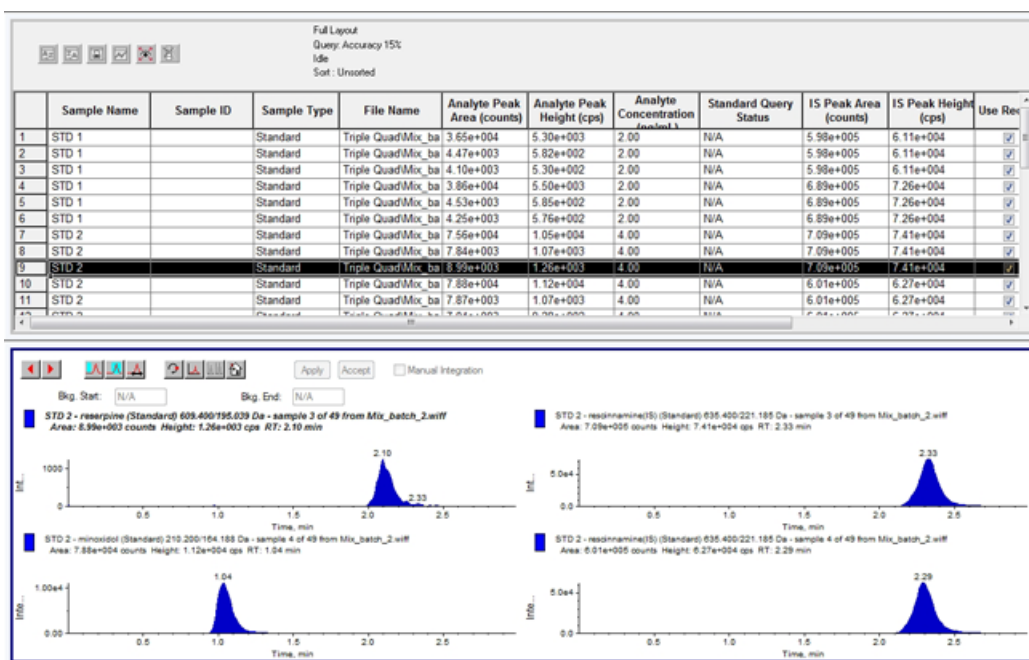Figure 9-4 Peak Review



4. Right-click in the pane and then click **Options**.
5. In the **Peak Review Options** dialog, in the **Appearance** section, change **Num. rows** to **1** and **Num. columns** to **2**.
6. In the **Automatic Zooming** section, click **Zoom Y axis to: 100% of largest peak** to show the entire peak.

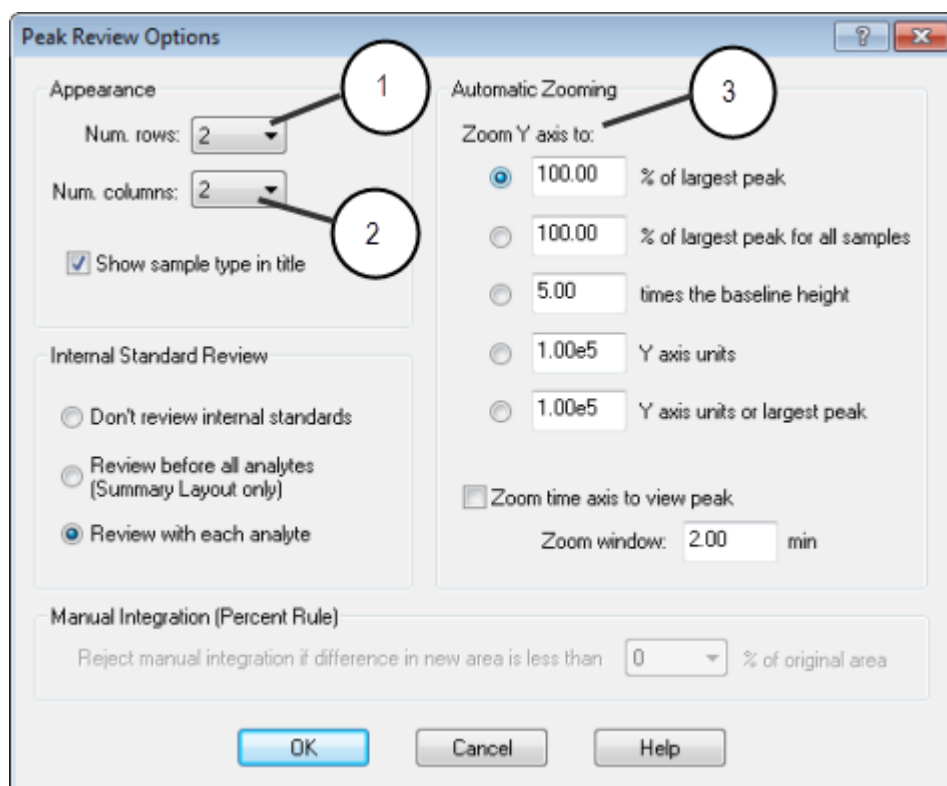Figure 9-5 Peak Review Options Dialog



Item	Definition
1	Number of rows
2	Number of columns
3	Zoom Y-axis to 100% of largest peak to show the entire peak

- Click **OK**.
- To move through the peaks, click the right-pointing arrow. Refer to [Figure 9-6](#) and [Figure 9-7 on page 102](#).
- Go to the second injection of standard 3.

In this example, the peak can be integrated closer to the baseline by selecting the **Specify Parameters** option.

Tip! To move to a specific peak in the Peak Review pane, select the corresponding row in the Results Table.

Figure 9-6 Peak Review Pane


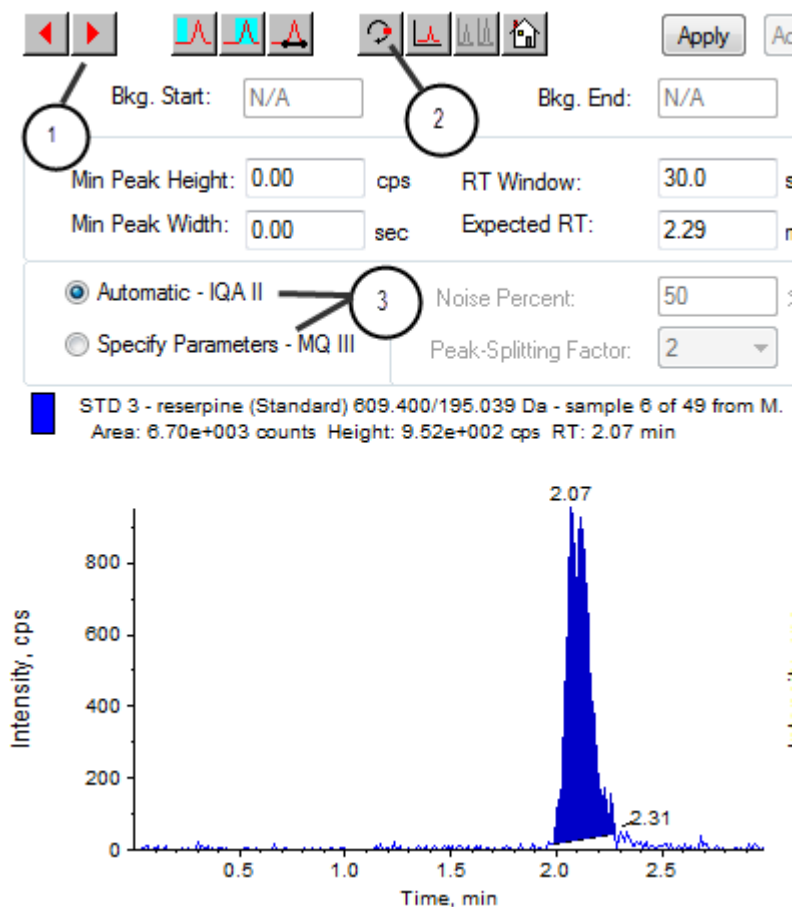


Figure 9-7 Peak Review Pane

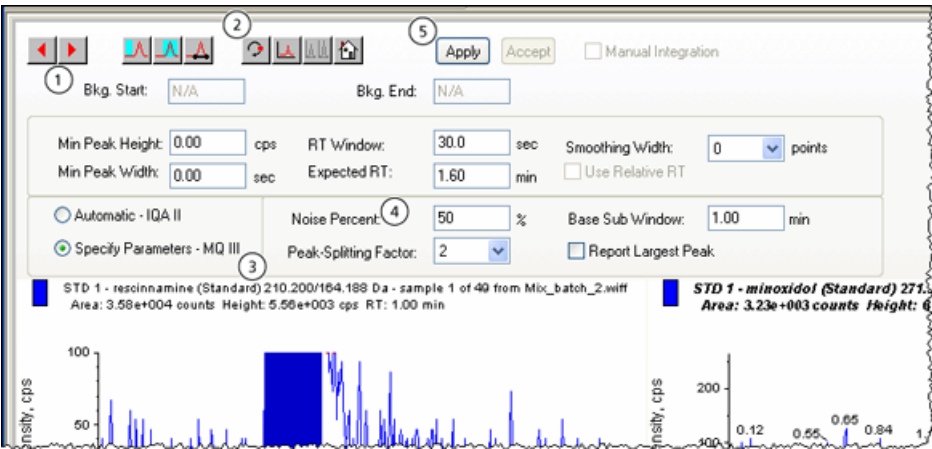

Item	Definition
1	Arrows: Click to move through the peaks.
2	Show or Hide Parameters: Click to show the integration parameters.
3	Integration parameters: Click to change the parameters.
4	Noise Percentage: Type a noise percent.
5	Apply: Click to integrate the parameters.

10. Click **Show or Hide Parameters** twice.

11. Click **Specify Parameters - MQ III**.

12. Change the **Noise Percent** value.

13. Click **Apply**.


The peak is integrated closer to the baseline.

14. If the change does not improve the peak integration, then adjust the **Noise Percent** parameter until the optimal value is found.

Note: The **Update Method** function only updates the algorithm values for that specific analyte (or internal standard) and not all analytes.

15. To update the algorithm for all peaks, right-click in the pane and then click **Update Method**.

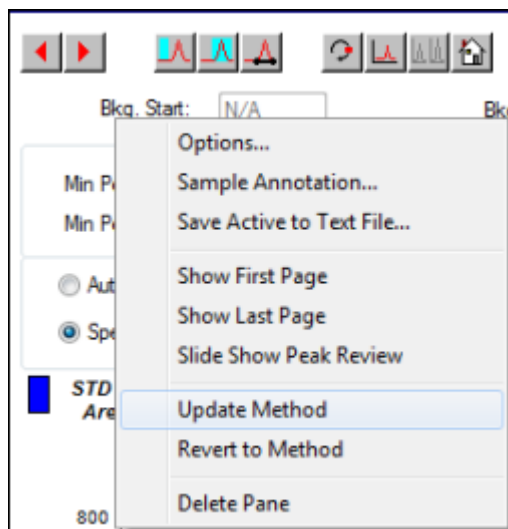Figure 9-8 Update Method



Manually Integrate Peaks

Manually integrating peaks should be done last, to limit person-to-person variability. Manually integrate peaks only if all the peaks have not been found after the algorithm parameters have been adjusted and updated. We recommend that users validate the results to determine whether manual integration is acceptable for specific applications.

Note: Peaks that are manually integrated, or where the algorithm was changed for only that peak, are identified in the Record Modified column of the Results Table, as are peaks that have algorithm parameter changes for a sample that are not applied to the entire analyte group.


1. In the **Peak Review** pane, click **Manual Integration Mode**.

Figure 9-9 Peak Review Pane: Manual Integration



2. Zoom in on the lower 10% of the peak.

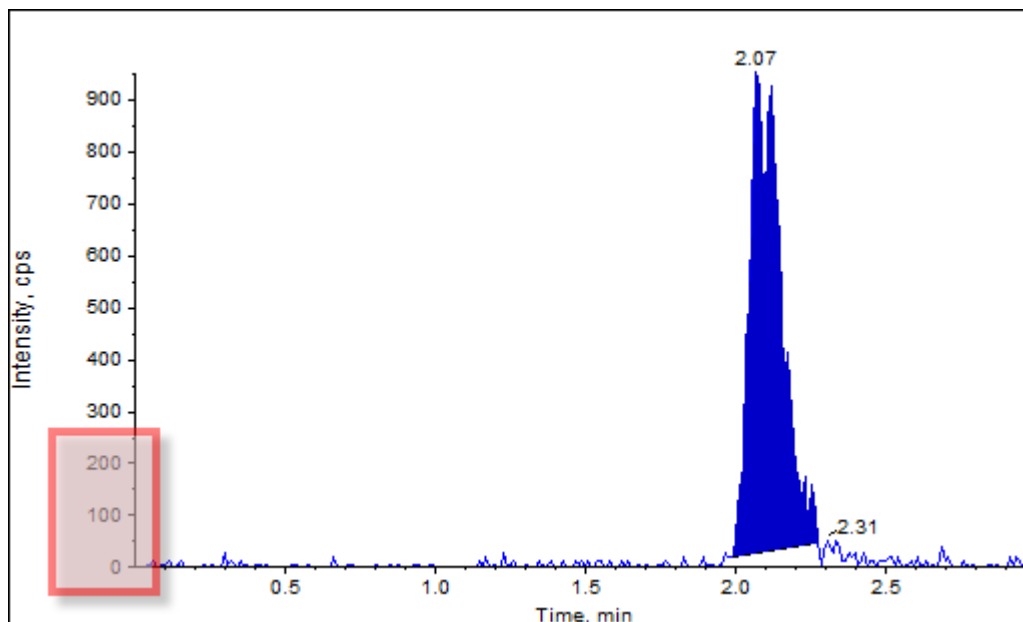Figure 9-10 Peak Review Pane: Zooming in on a Peak



3. Move the cross-hair to where the start of the peak is to be defined and then drag the cross-hair to where the end of the peak is to be defined.

The software shades the area bounded by the base and sides of the peak. Peak parameters are gray as they are no longer applicable because the peak was drawn manually.

4. Do one of the following:
 - To make this change permanent, click **Accept**.
 - To discard the changes, clear the **Manual Integration** check box.

Tip! If a peak was correct as originally selected, right-click the peak and then click **Revert to Method**.

Peak Review Right-Click Menu

Right-click in the **Peak Review** window or pane to access the options shown in [Table 9-2](#).

Table 9-2 Peak Review Right-click Menu

Menu	Function
Options	Opens the Peak Review Options dialog.
Sample Annotation	Opens the Sample Annotation dialog.
Save Active to Text File	Saves the selected peak as a text file.
Show First Page	Goes to the first sample.
Show Last Page	Goes to the last sample.
Slide Show Peak Review	Opens the slide show.
Update Method	Updates the algorithm for all peaks.
Revert to Method	Selects a redefined peak based on the current quantitation method.
Delete Pane	Deletes the active pane.

Calibration Curves

Use calibration curves to find the calculated concentration of samples, including quality control (QC) samples. QC samples are added to a batch to estimate the data quality and accuracy of standards in the batch. QC samples have known analyte concentrations but are treated as unknowns so that the measured concentrations can be compared to the actual value.

The calibration curve is generated by plotting the concentration of the standard versus its area or height. If an internal standard is used, then the ratio of the standard concentration or internal standard versus the ratio of the standard peak height or area to the internal standard peak height or area is plotted. The area or height ratio of a sample is then applied to this curve to find the concentration of the sample, as shown in the Results Table. A regression equation is generated by this calibration curve according to the regression that was specified. The regression equation is used to calculate the concentration of the unknown samples.

View Calibration Curves

The user can view the calibration curve and change the regression options in an open Results Table. If two or more Results Tables are open, then the calibration curves can be overlaid. To overlay curves, make sure that the method used to create the tables is the same.

Plot a calibration curve to see the curve used for regression. The Calculated Concentration field in the Results Table reflects any changes resulting from the fit of the curve to the points of the standard.


Operating Instructions — Analyze and Process Quantitative Data

Note: This option is available only when a Results Table is open in the workspace.

1. Open a Results Table.
2. Click **Tools > Calibration > Pane**.


The Calibration Curve pane containing the calibration curve opens.

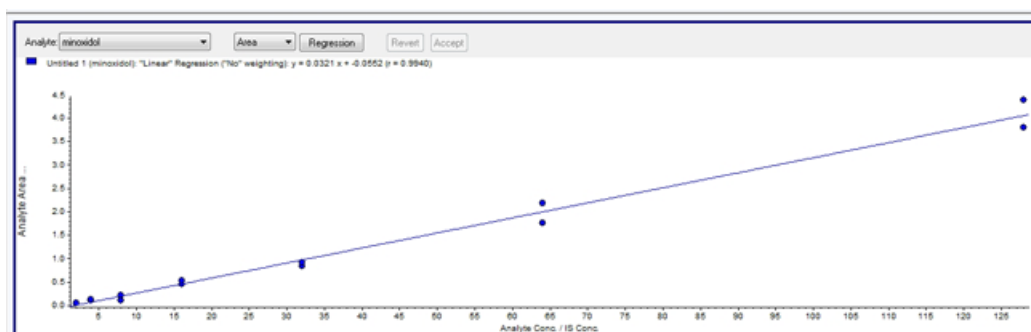Figure 9-11 Calibration Curve



3. If there is more than one analyte, then use the following steps to view the calibration curve for another analyte:
 - a. From the **Analyte** list, select an analyte.
 - b. If required, from the next list, select **Area** or **Height**.
4. To change the regression options for the calibration curve, do the following:
 - a. Click **Regression**.

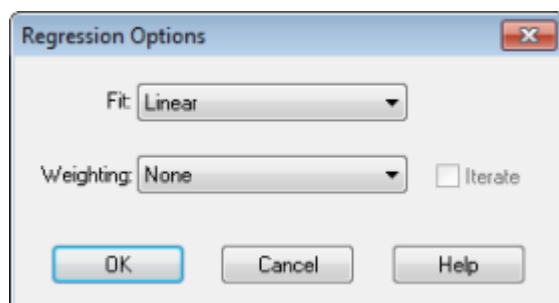Figure 9-12 Regression Options Dialog



- b. Select **Linear** in the **Fit** list.
- c. Select **1 / x** in the **Weighting** list.
- d. Click **OK**.

The calibration curve opens. The user can review individual peaks on the curve or exclude points from the curve to produce a better curve.


5. If required, repeat these steps to create a more appropriate curve.
6. To save the changes, click **Accept**.

Overlay Calibration Curves

Tip! To examine the curve for one table more closely, right-click the curve and then click **Active Plot**. Select the curve to be plotted on top.

1. With two or more open Results Tables, view a calibration curve for one of the tables.
2. Right-click the calibration curve and then click **Overlay**.

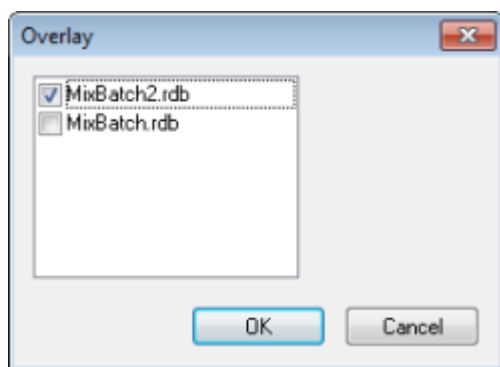Figure 9-13 Overlay Dialog



3. Select the tables to overlay with the current curve.
4. Click **OK**.

The software plots the curves for all selected tables on the same graph.

Calibration Curve Right-Click Menu

Right-click in the Calibration window or pane table to access the options shown in [Table 9-3](#).

Table 9-3 Calibration Curve Right-click Menu

Menu	Function
Exclude (Include)	Right-click the point and then click Exclude to exclude the point from the curve. Right-click the point and then click Include to include the point.
Exclude All Analytes (Include All Analytes)	Right-click a point and then click Exclude All Analytes to exclude all the analytes from the curve. Right-click a point and then click Include All Analytes to include the points.
Show Peak	Reviews an individual peak.
Overlay	Overlays two graphs.
Active Plot	Determines which plot is active.
Legend	Shows the graph legend.
Log Scale X Axis*	Uses a log scale for the X-axis.
Log Scale Y Axis*	Uses a log scale for the Y-axis.
Delete Pane	Deletes the active pane.
Home Graph	Scales the graph to its original size
* A log scale arranges the data points in a more manageable view so that the effect of all points can be monitored simultaneously. For this view, select Log Scale Y Axis versus Log Scale X and not just the log of one axis.	

Sample Statistics

Use the Statistics window to view the statistics samples, typically for standards and QCs (quality controls). The data from each available batch in the Results Table opens in tabular form in the grid and a row of data is shown for each standard or QC concentration.

View the Statistics for Standards and QCs

When more than one **Results Table** is open, statistical information about the standards and QCs for additional batches in the **Statistics** window can be obtained. This facilitates comparison of results between batches and identification of trends in the standards or QCs.

1. Open a Results Table.
2. Click **Tools > Statistics**.
3. Select **Concentration** from the **Statistics Metric** list.
4. Select an analyte in the **Analyte Name** field.

5. Select **Standard** in the **Sample Type** field.

The results are shown.

6. Look at the **%CV** and **Accuracy** columns.

The **%CV** shows the coefficient of variance between the measurements of a single parameter, for example the area. **Accuracy** shows how close the plotted point is to the interpolated value.

7. If required, select the **Display Low/High values** check box and then examine the **Low**, **High** values, and **Mean** for each row in the grid. Each row represents standards that have the same concentration levels.
8. Select another analyte.

The results are shown on a per-analyte basis.

9. To check for **Quality Control** variations at the same concentration levels, select **QC** in the **Sample Type** field.

Compare Results Between Batches

When more than one Results Table is shown, obtain statistical information about the standards and QCs for additional batches in the Statistics window. Normally results are compared between batches to look for trends in the standards or QCs or to verify that the method is valid.

For two or more open Results Tables, compare results in the Statistics window. Both sets of statistics are shown in the Statistics window.

The number of analytes and the analyte names must be the same for the data to be combined in the Statistics pane.

1. Open a Results Table.
2. Click **Tools > Statistics**.
3. Do one of the following:
 - To arrange the results by **Results Table**, select **Group By Batch** in the **Conc. as Rows** list.
 - To arrange the results in order of concentration, select **Group By Concentration** in the **Conc. as Rows** list.
 - To arrange the results in order of concentration without a row showing the statistics for each group or batch, select **Group By Concentration (no All)** in the **Conc. as Rows** list.

The software sorts the results. At the end of each batch or group, one or two additional rows are shown: **All** (statistics for all Results Tables in that group) and **Average** (statistics on the statistics for that batch or group).

Operating Instructions — Analyze and Explore Data

10

Use the sample files installed in the Example folder to learn how to view and analyze data using the most common analysis and processing tools. For more information about the following topics, refer to the *Advanced User Guide*:

- Labeling graphs
- Overlaying and summing spectra or chromatograms
- Performing background subtractions
- Smoothing algorithms
- Working with smoothed data
- Working with centroid data
- Working with contour plots
- Working with the fragment interpretation tool
- Working with library databases and library records

Open Data Files

Tip! To turn off the automatic update on the mass spectrum, right-click the mass spectrum and then click **Show Last Scan**. If there is a check mark beside **Show Last Scan**, then the spectrum will update in real-time.

1. On the Navigation bar, under **Explore**, double-click **Open Data File**.

The Select Sample dialog is shown.

2. In the **Data Files** list, navigate to the data file to open, select a sample, and then click **OK**.

The data acquired from the sample is shown. If data is still being acquired, then the mass spectrum, DAD/UV trace, and TIC continue to update automatically.

Navigate Between Samples in a Data File

Note: If samples were saved in separate data files, then open each file individually.

Table D-5 on page 180 shows the navigation icons used in this procedure.

- Open a data file that contains multiple samples and then do one of the following:
 - Click the icon with the arrow pointing to the right to skip to the next sample in the data file.
 - Click the icon with the arrow curving to the right to skip to a non-sequential sample.
 - In the Select Sample dialog, from the **Sample** list, select the sample to view.
 - Click the icon with the arrow pointing to the left to go to the previous sample in the data file.

View Experimental Conditions

The experimental conditions used to collect data are stored in the data file with the results. The information contains the details of the acquisition method used: the MS acquisition method (that is, the number of periods, experiments, and cycles) including instrument parameters and the HPLC device method (LC pump flow rate). In addition, it also contains the MS resolution and mass calibration tables used for the sample acquisition. Table 10-1 shows the software functionality available when the user views the file information.

Note: If data is acquired from more than one sample into the same wiff file, then the file information pane does not refresh automatically while scrolling through the samples. Close the file information pane and then reopen it to view the details for the next sample in the wiff file.

- Click **Explore > Show > Show File Information**.

The File Information pane opens below the graph.

Tip! To create an acquisition method from the **File Information** pane, right-click the **File Information** pane and then click **Save Acquisition Method**.

Table 10-1 Right-click Menu for Show File Information Pane

Menu	Function
Copy	Copies the selected data.
Paste	Pastes data.
Select All	Selects all the data in the pane.
Save To File	Saves data as an rtf file.
Font	Changes the font.
Save Acquisition Method	Saves the acquisition method as a dam file.

Table 10-1 Right-click Menu for Show File Information Pane (continued)

Menu	Function
Save Acquisition Method to CompoundDB	Opens the Specify Compound Information dialog. Select the IDs and molecular weights to be saved in the compound database.
Delete Pane	Deletes the selected pane.

Show Data in Tables

1. Open a data file.
2. Click **Explore > Show > Show List Data**.

The data is shown in a pane below the graph.

Figure 10-1 Peak List Tab

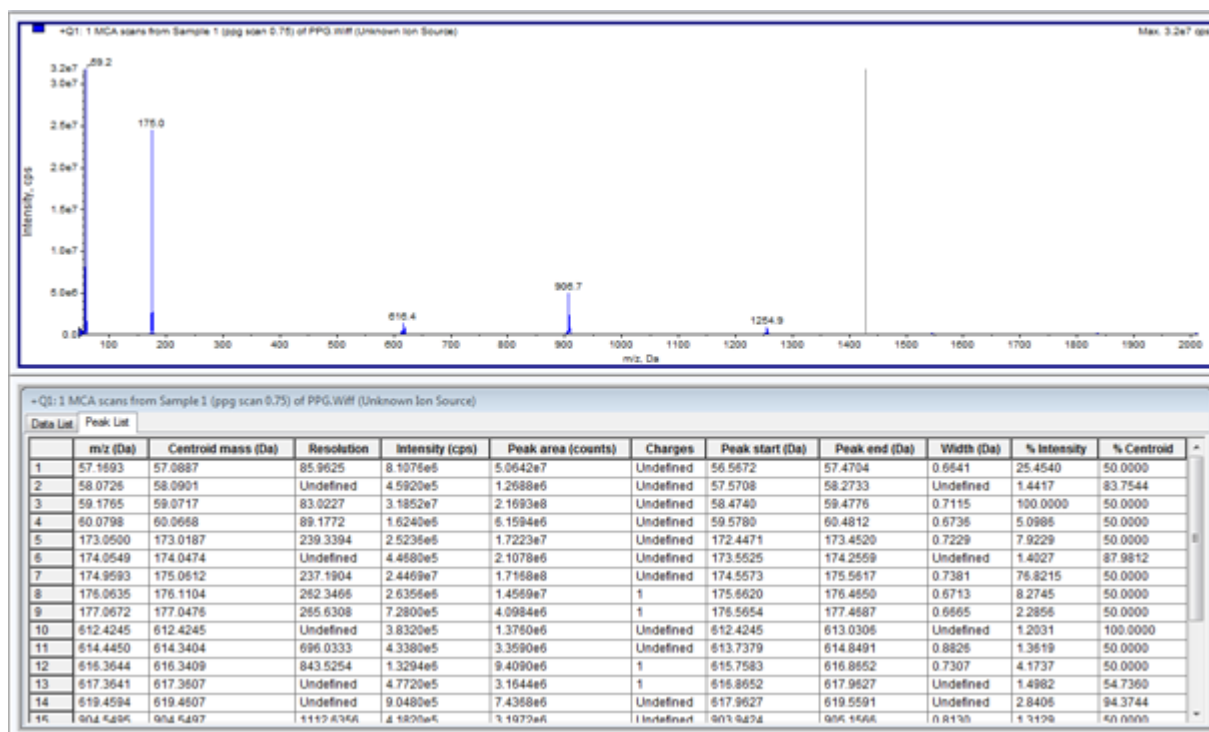


Table 10-2 Right-click Menu for Spectral Peak List Tab

Menu	Function
Column Options	Opens the Select Columns for Peak List dialog.
Save As Text	Saves the data as a .txt file.
Delete Pane	Deletes the selected pane.

Table 10-3 Right-click Menu for Chromatographic Peak List Tab

Menu	Function
Show Peaks in Graph	Show the peaks in two colors in the graph.
IntelliQuan Parameters	Opens the IntelliQuan dialog.
Save As Text	Saves the data as a txt file.
Delete Pane	Deletes the selected pane.

Show ADC Data


ADC (analog-to-digital converter) data is acquired from a secondary detector (for example from a UV detector through an ADC card), and is useful for comparison with mass spectrometer data. To have ADC data available, acquire the data and the mass spectrometer data simultaneously and then save it in the same file.

1. Make sure that the **Example** folder is selected.
2. On the Navigation bar, under **Explore**, double-click **Open Data File**.

The Select Sample dialog opens.

3. In the **Data Files** field, double-click **Devices** and then click **Adc16chan.wiff**.
4. In the **Samples** list, select a sample, and then click **OK**.
5. Click **Explore > Show > Show ADC Data**.

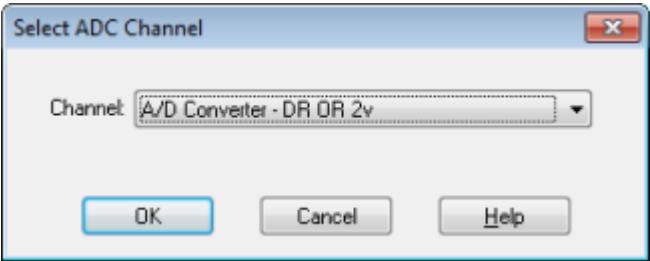Figure 10-2 Select ADC Channel Dialog




6. In the **Channel** list, select a channel, and then click **OK**.

The ADC data is shown in a new pane beneath the active pane.

Show Basic Quantitative Data

1. Open a data file.
2. Click **Explore > Show > Show List Data**.

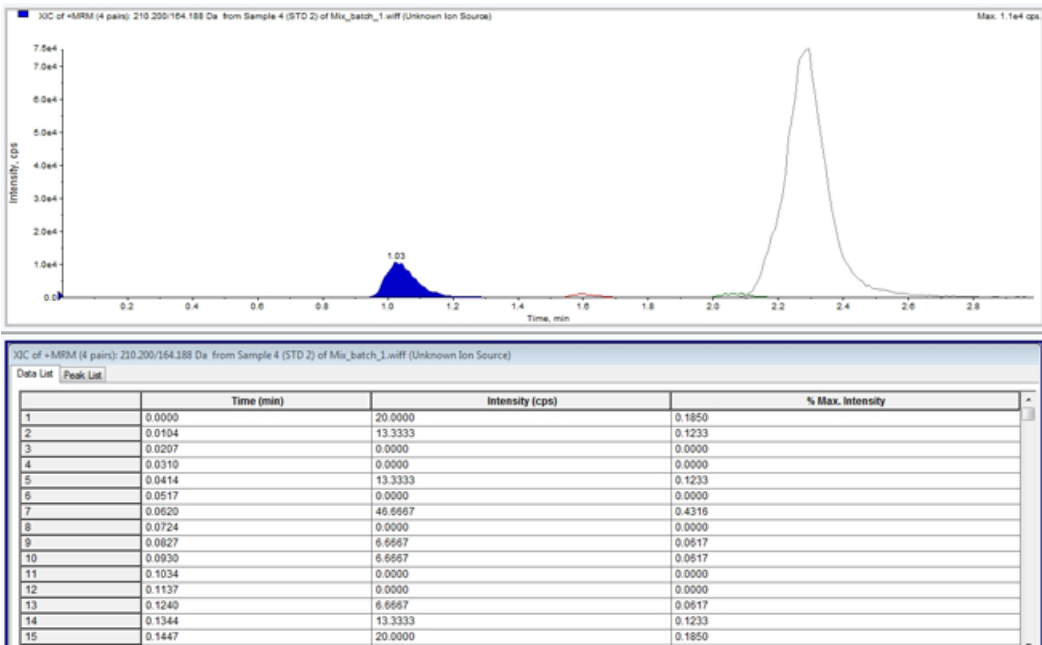Figure 10-3 List Data



3. In the Peak List tab, right-click and then select **Show Peaks in Graph**.

Peaks are shown in two colors.

4. To change the peak finding algorithm settings, right-click and then select either **Analyst Classic Parameters** or **IntelliQuan Parameters**, whichever is active.
5. (Optional) To remove the colored peaks, right-click in the Peak List tab and then clear **Show Peaks in Graph**.

Chromatograms

A chromatogram is a graphical view of the data obtained from the analysis of a sample. It plots the signal intensity along an axis that shows either time or scan number. For more information about software functionality available for chromatograms, refer to [Table 10-6 on page 124](#).

The software plots intensity, in counts per second (cps), on the Y-axis against time on the X-axis. Peaks above a set threshold are labeled automatically. In the case of LC-MS, the chromatogram is often shown as a function of time. [Table 10-4](#) contains the a description of the types of chromatograms.

Refer to [Table 10-8 on page 127](#) for more information about using the available icons.

Table 10-4 Types of Chromatograms

Types of Chromatograms	Purpose
TIC (Total Ion Chromatogram)	<p>A chromatographic view generated by plotting the intensity of all ions in a scan against time or scan number.</p> <p>When a data file is opened, it is preset to open as a TIC. If the experiment contains only one scan, then it is shown as a spectrum.</p> <p>If the MCA check box is selected during acquisition of the data file, then the data file opens to the mass spectrum. If the MCA check box is not selected, then the data file opens as the TIC.</p>
XIC (Extracted Ion Chromatogram)	<p>An ion chromatogram created by taking intensity values at a single, discrete mass value, or a mass range, from a series of mass spectral scans. It indicates the behavior of a given mass, or mass range, as a function of time.</p>
BPC (Base Peak Chromatogram)	<p>A chromatographic plot that shows the intensity of the most intense ion within a scan versus time or scan number.</p>
TWC (Total Wavelength Chromatogram)	<p>A chromatographic view created by summing all of the absorbance values in the acquired wavelength range and then plotting the values against time. It consists of the summed absorbances of all ions in a scan plotted against time in a chromatographic pane.</p>

Table 10-4 Types of Chromatograms (continued)

Types of Chromatograms	Purpose
XWC (Extracted Wavelength Chromatogram)	A subset of TWC. An XWC shows the absorbance for a single wavelength or the sum of the absorbance for a range of wavelengths.
DAD (Diode Array Detector)	A UV detector that monitors the absorption spectrum of eluting compounds at one or more wavelengths.

Show TICs from a Spectrum

To see an example data file, make sure that the Example project is selected. Open the Triple Quad folder, and then open the Mix_batch_1.wiff file.

- Click **Explore > Show > Show TIC**.

The TIC opens in a new pane.

Tip! Right-click inside a pane containing a spectrum and then click **Show TIC**.

Show a Spectrum from a TIC

A TIC is created by summing the intensity contributions of all ions from a series of mass scans. Use the TIC to view an entire data set in a single pane. It consists of the summed intensities of all ions in a scan plotted against time in a chromatographic pane. If the data contains results from multiple experiments, then the user can create individual TICs for each experiment and another TIC that represents the sum of all experiments.


When a data file is opened, it is preset to appear as a TIC. However, if the experiment contains only one scan, it is shown as a spectrum. If the user selects the MCA check box before acquiring the data file, then the data file opens to the mass spectrum. If the MCA check box is not selected, then the data file opens with the TIC.

1. In a pane containing a TIC, select a range.
2. Click **Explore > Show > Show Spectrum**.

The spectrum opens in a new pane.

Tip! Double-click in the TIC pane at a particular time to show the spectrum.

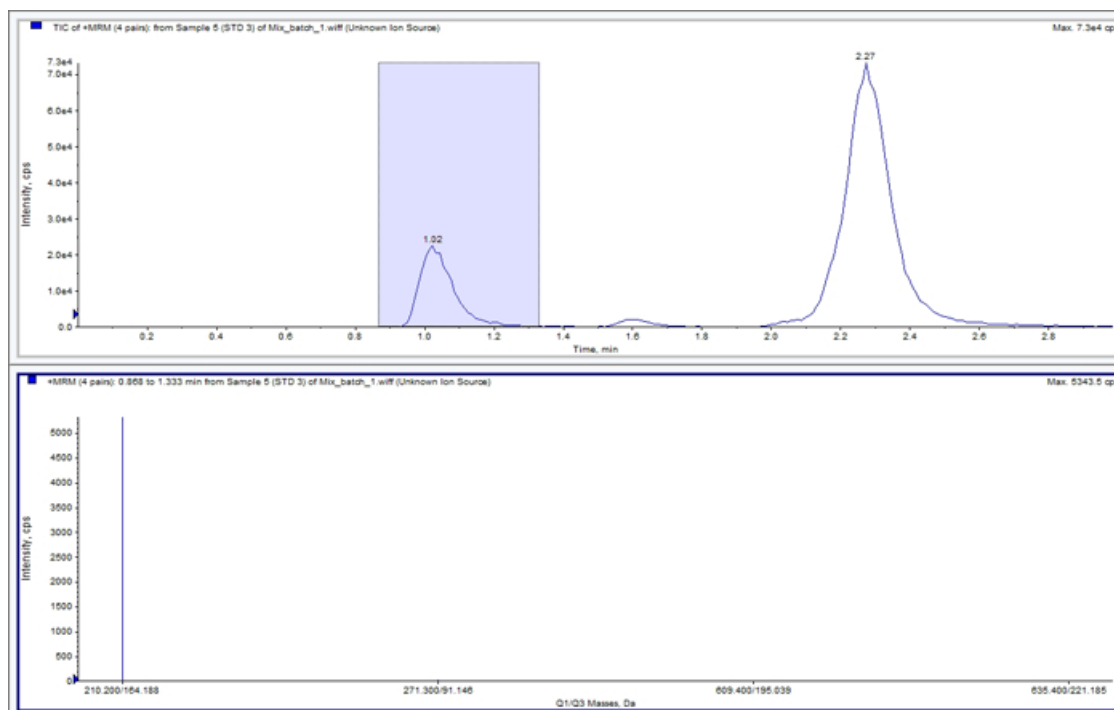Figure 10-4 Example of a TIC



About Generating XICs

XICs can be generated only from single-period, single-experiment chromatograms or spectra. To obtain an XIC from multi-period or multi-experiment data, split the data into separate panes by clicking the triangle under the X-axis. Refer to [Table 10-8 on page 127](#) for more information about using the available icons.

Several methods are available for extracting ions to generate an XIC, depending on whether chromatographic or spectral data is used. [Table 10-5](#) contains a summary of methods that can be used with chromatograms and spectra.

Table 10-5 Summary of XIC Generation Methods

Method	Use with Chromatogram	Use with Spectrum	Extraction
Selected range	No	Yes	Extracts ions from a selected area in a spectrum.
Maximum	No	Yes	Extracts ions from a selected area in a spectrum using the most intense peak in the selected area. This option creates an XIC using the maximum mass from the selected spectral range.
Base peak masses	Yes	Yes	Can be used only with Base Peak Chromatograms (BPCs). Use the Use Base Peak Masses command to extract ions results in an XIC with a different colored trace for each mass. If the selection includes multiple peaks, then the resulting XIC will have an equal number of colored traces, one for each mass.
Specified masses	Yes	Yes	Extracts ions from any type of spectrum or chromatogram. Select up to ten start and stop masses for which to generate XICs.

Generate an XIC Using a Selected Range

1. Open a data file containing spectra.
2. Select a range by pressing the left mouse button at the start of the range, dragging the cursor to the stop point, and then releasing the left mouse button.

The selection is indicated in blue.

3. Click **Explore > Extract Ions > Use Range**.

An XIC of the selection opens in a pane below the spectrum pane. The experiment information at the top of the pane contains the mass range and the maximum intensity in counts per second.

Generate an XIC Using the Maximum Peak

1. Open a data file containing spectra.
2. Select a range in a spectrum.

The selection is indicated in blue.

3. Click **Explore > Extract Ions > Use Maximum.**

An XIC of the maximum peak specified selection opens below the spectrum pane. The experiment information at the top of the pane contains the mass range and the maximum intensity in counts per second.

Generate an XIC Using Base Peak Masses

1. Open a data file containing spectra.
2. In a BPC, select the peak from which to extract ions.

The selection is indicated in blue.


3. Click **Explore > Extract Ions > Use Base Peak Masses.**

An XIC of the specified selection opens below the spectrum pane. The experiment information at the top of the pane shows the mass range and the maximum intensity in counts per second.

Extract Ion by Selecting Masses

1. Open a spectrum or chromatogram.
2. Click **Explore > Extract Ions > Use Dialog.**

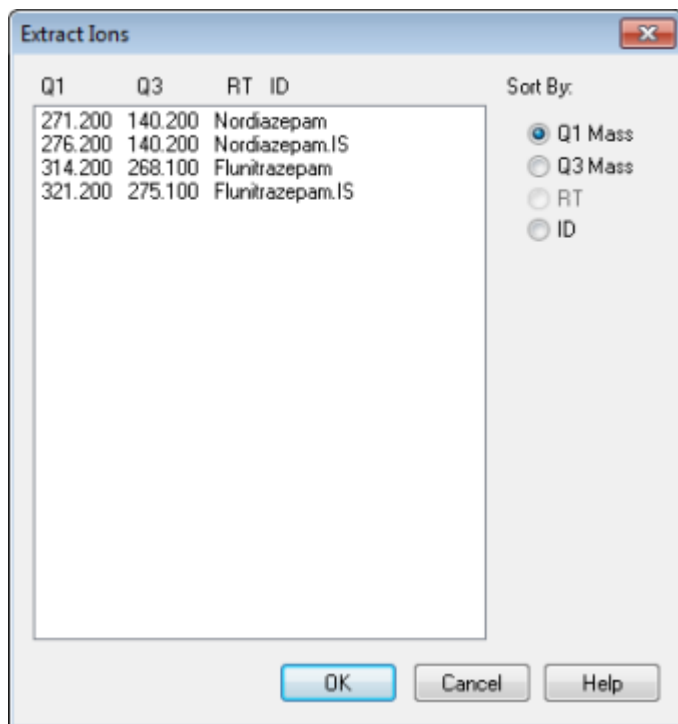Figure 10-5 Extract Ions Dialog



3. Type the values for each XIC to be created. If a stop value is not typed, then the range is defined by the start value.
 - In the **Start** field, type the start value (lower value) for the mass range.
 - In the **Stop** field, type the stop value (higher value) for the mass range.
4. Click **OK**.

An XIC of the selection opens below the chromatogram pane. The experiment information at the top of the pane includes the masses and the maximum intensity in counts per second.

Generate BPCs

BPCs can be generated only from single-period, single-experiment data.


1. Open a data file.
2. Select an area within a TIC.

The selection is indicated in blue.

3. Click **Explore > Show > Show Base Peak Chromatogram**.

The selections are shown in the **Start Time** and **End Time** fields.

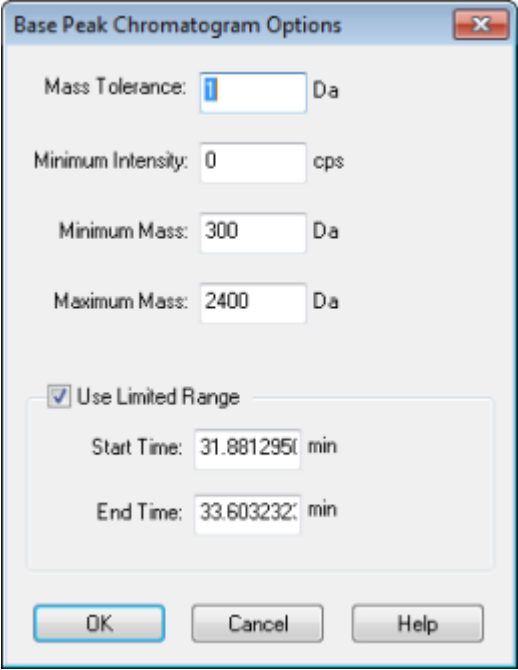Figure 10-6 Base Peak Chromatogram Options

A screenshot of the 'Base Peak Chromatogram Options' dialog box. The dialog has a title bar with a close button. It contains several input fields: 'Mass Tolerance' with a value of 1 and unit 'Da'; 'Minimum Intensity' with a value of 0 and unit 'cps'; 'Minimum Mass' with a value of 300 and unit 'Da'; and 'Maximum Mass' with a value of 2400 and unit 'Da'. Below these is a checked checkbox labeled 'Use Limited Range'. Under this checkbox are two more input fields: 'Start Time' with a value of 31.881295 and unit 'min', and 'End Time' with a value of 33.603232 and unit 'min'. At the bottom are three buttons: 'OK', 'Cancel', and 'Help'.

4. In the **Mass Tolerance** field, type the value to indicate the mass range used to find a peak. The software finds the peak using a value twice the typed range (\pm the mass value).
5. Type the intensity below which peaks are ignored by the algorithm in the **Minimum Intensity** field.
6. Type the mass that determines the beginning of the scan range in the **Minimum Mass** field.
7. Type the mass that determines the end of the scan range in the **Maximum Mass** field.
8. To set the start and end times, select the **Use Limited Range** check box and do the following:
 - In the **Start Time** field, type the time that determines the start of the experiment.
 - In the **End Time** field, type the time that determines the end of the experiment.
9. Click **OK**.


The BPC is generated in a new pane.

Generate XWCs

An XWC is a wavelength chromatogram created by taking intensity values at a single wavelength, or by the sum of the absorbance for a range of several wavelengths. Up to three ranges can be extracted from a DAD spectrum to generate the XWC. Refer to [Table 10-8 on page 127](#) for more information about using the available icons.

1. Open a data file that contains a DAD spectrum.
2. Right-click anywhere in the pane and then click **Extract Wavelengths**.

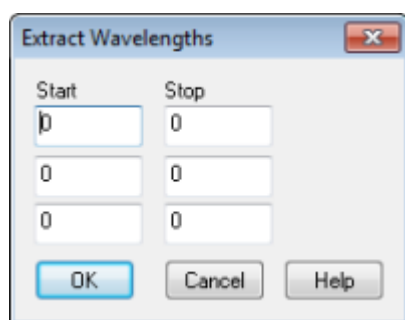Figure 10-7 Extract Wavelengths Dialog



3. Type **Start** and **Stop** values.
4. Click **OK**.

The XWC is shown in a pane below the DAD spectrum.

Generate DAD Data

Like mass spectrometer data, DAD data can be viewed in chromatogram or spectrum form. Users can view the DAD spectrum for a single point in time, or for a range of time as a Total Wavelength Chromatogram (TWC).

1. Open a data file containing data acquired with a DAD.

The TWC, which is analogous to a TIC, opens in a pane below the TIC.

2. In the TWC pane, click a point to select a single point in time, or highlight an area of the spectrum to select a range of time.
3. Click **Explore > Show > Show DAD Spectrum**.

The DAD spectrum opens in a pane below the TWC. The Y-axis shows the absorbance and the X-axis shows the wavelength.

Tip! If the pane with the TWC is closed, then click a point anywhere in the TWC to open it again. Click **Explore > Show > Show DAD TWC**.

Generate TWCs

A TWC is a less commonly used chromatogram. It shows the total absorbance (mAU) as a function of time. The TWC provides a way of viewing an entire data set in a single pane. It consists of the summed absorbances of all ions in a scan plotted against time in a chromatographic pane. If the data contains results from multiple experiments, then create individual TWCs for each experiment and another TWC that represents the sum of all experiments.

A TWC shows total absorbance (mAU) on the Y-axis plotted against time on the X-axis. Refer to [Table 10-8 on page 127](#) for more information about using the available icons.

1. Open a data file that contains a DAD spectrum.
2. Click **Explore > Show > Show DAD TWC**.

The TWC is shown in a pane below the DAD spectrum.

Tip! Right-click inside the pane containing the DAD spectrum and then click **Show DAD TWC**.

Adjust the Threshold

The threshold is an invisible line drawn parallel to the X-axis of a graph that sets a limit below which the software will not include peaks in a spectrum. The line has a handle, represented by a blue triangle to the left of the Y-axis. Click the blue triangle to view a dotted line that represents the threshold. The threshold can be raised or lowered, but changing the threshold value does not change the data. The software does not label any peaks in the region that lies below the threshold.

1. Open a data file.
2. Do one of the following:
 - To raise the threshold, drag the blue triangle up the Y-axis. To lower the threshold, drag the blue triangle down.
 - Click **Explore > Set Threshold**. In the Threshold Options dialog that opens, type the threshold value and then click **OK**.
 - Click **Explore > Threshold**.

The graph updates to show the new threshold. Peak labeling and the peak list are also updated.

Tip! To view the current threshold value, move the pointer over the threshold handle.

Chromatogram Panes

Table 10-6 Right-click Menu for Chromatogram Panes

Menu	Function
List Data	Lists the data points and integrates the peaks found in chromatograms.
Show Spectrum	Generates a new pane containing the spectrum.
Show Contour Plot	Shows a color-coded plot of a data set, where the color represents the intensity of the data at that point. Only certain MS modes are supported.
Extract Ions	Extracts a specific ion or set of ions from a selected pane and then generates a new pane containing a chromatogram for the specific ions.
Show Base Peak Chromatogram	Generates a new pane containing a base peak chromatogram.
Show ADC Data	Generates a new pane containing the ADC data trace, if acquired.
Show UV Detector Data	Generates a new pane containing the UV data trace, if acquired.
Spectral Arithmetic Wizard	Opens the Spectral Arithmetic Wizard.
Save to Text File	Generates a text file containing the data in a pane, which can be opened in Microsoft Excel or other programs.
Save Explore History	Saves information about changes to processing parameters, also called processing options, that were made when a wiff file was processed in Explore mode. The processing history is stored in a file with an eph (Explore Processing History) extension.
Add Caption	Adds a caption at the cursor point in the pane.
Add User Text	Adds a text box at cursor point in the pane.
Set Subtract Range	Sets the subtract range in the pane.
Clear Subtract Range	Clears the subtract range in the pane.
Subtract Range Locked	Locks or unlocks the subtract ranges. If the subtract ranges are not locked, then each subtract range can be moved independently. The subtract ranges are preset to locked.
Delete Pane	Deletes the selected pane.

Spectra Panes

Table 10-7 Right-click Menu for Spectra Panes

Menu	Function
List Data	Lists the data points and integrates chromatograms.
Show TIC	Generates a new pane containing the TIC.
Extract Ions (Use Range)	Extracts a specific ion or set of ions from a selected pane and then generates a new pane containing a chromatogram for the specific ions.
Extract Ions (Use Maximum)	Extracts ions using the most intense peak in a selected area.
Save to Text File	Generates a text file of the pane, which can be opened in Microsoft Excel or other programs.
Save Explore History	Saves information about changes to processing parameters, also called Processing Options, that were made when a wiff file was processed in Explore mode. The processing history is stored in a file with an eph (Explore Processing History) extension.
Add Caption	Adds a caption at the cursor position in the pane.
Add User Text	Adds a text box at the cursor position in the pane.
Show Last Scan	Shows the scan prior to the selection.
Select Peaks For Label	In this dialog, select the parameters to reduce peak labeling.
Delete Pane	Deletes the selected pane.
Add a Record	Adds records and compound-related data, including spectra, to the library. An active spectrum is required to perform this task.
Search Library	Searches the library without constraints or with previously saved constraints.
Set Search Constraints	Searches the library using the criteria typed in Search Constraints dialog.

Graphical Data Processing

Graphical data can be processed many ways. This section provides information and procedures for using some of the most commonly used tools.

The user can zoom in on part of a graph to view a particular peak or an area in greater detail in both spectra and chromatograms. The user can also zoom in repeatedly to view smaller peaks.

Graphs

The same data can be examined in different ways. Data can also be kept for comparison purposes before performing processing operations such as smoothing or subtraction.

A window contains one or more panes arranged in such a way that all the panes are fully visible and that they do not overlap.

Panes might be of a variable or fixed size. Panes are automatically tiled within the window and are arranged into column and row format. If the size of a window is changed, then the panes within the window change in size to accommodate the new size. A window cannot be sized to the point where any of the panes become smaller than their minimum size.

Two or more windows or panes containing similar data can be linked, for example, spectra with similar mass ranges. As one pane or window is zoomed in, the other pane zooms in simultaneously. For example, the user can link an XIC to the BPC from which the XIC was extracted. Zooming in the BPC also zooms the XIC, so that both chromatograms show the same magnification.

Manage Data

Data can be compared or examined in different ways. Users might want to keep the data for comparison purposes before performing processing operations such as smoothing or subtraction.

A window contains one or more panes, arranged in such a way that all the panes are fully visible and they do not overlap.

Panes can be of variable or fixed size. Panes are automatically tiled within the window and are arranged into column and row format. If window size is changed, then the panes within the window change in size to accommodate the resizing. A window cannot be resized to the point where any of the panes would become smaller than its minimum size.

Two or more windows or panes containing similar data can be linked, for example, spectra with similar mass ranges. When the user zooms in one pane or window, the other pane zooms simultaneously. For example, the user can link an XIC to the BPC from which it was extracted. Zooming in the BPC also zooms the XIC, so that both chromatograms are shown with the same magnification.

- Use the following menu options or icons to manage data in graphs.

Table 10-8 Graph Options











To do this...	use this menu option...	...or click this icon
Copy a graph to a new window	Select the graph to copy. Click Explore > Duplicate Data > In New Window .	
Rescale a graph to its original size	Select the graph. Click Explore > Home Graph .	
Move a pane	<ul style="list-style-type: none"> Select the graph. Click Window > Move Pane. Select the pane or window and then drag it to the new position. This position can be inside the same window or within another window. <p>A four-headed arrow is shown when the cursor is on the boundary of the active window or pane.</p> <ul style="list-style-type: none"> If the pane is at the top or bottom of the target pane, then the pane moves above or below that pane, respectively. If the pane is at the left or right of the target pane, then the pane moves to the left or right of that pane, respectively. If the pane is at any other position, then the pane moves to the target row. The drop shadow of the pane as the pane is moved indicates its new position. 	
Link panes	<ol style="list-style-type: none"> With the two graphs open, click one to make that pane active. Click Explore > Link and then click the other pane. 	
Remove linking	Close one of the panes. Click Explore > Remove Link .	
Delete a pane	Select the graph. Click Window > Delete Pane .	
Lock a pane	Select the graph. Click Window > Lock Panes .	
Hide a pane	Select the graph. Click Window > Hide Pane .	

Table 10-8 Graph Options (continued)

To do this...	use this menu option...	...or click this icon
Maximize a pane	Select the graph. Click Window > Maximize Pane .	
Tile panes	Select the graph. Click Window > Tile all Panes .	

Zoom In on the Y-axis

1. Move the pointer to the left of the Y-axis to either side of the area to be expanded and then drag away from the starting point in a vertical direction while holding the left mouse button.

A box is drawn along the y-axis representing the new scale.

Note: Take care when zooming in on the baseline. Zoom in too far and the zoom-in box closes.

2. Release the mouse button to draw the graph to the new scale.

Zoom In on the X-axis

Tip! To return the graph to the original scale, double-click either axis. To restore the entire graph to the original scale, click **Explore > Home Graph**.

1. Move the pointer under the X-axis to either side of the area to be expanded and then drag away from the starting point in a horizontal direction while holding the left mouse button.
2. Release the mouse button to draw the graph to the new scale.

The Reporter software extends the reporting functionality available in the Analyst[®] software.

We recommend that users validate the results if a modified Reporter template or one that contains a query is used.

The Reporter software can be used to create custom reports with Microsoft Word and Excel (2010, 2013, or 2016).

The Reporter software has the following features:


- Provides a variety of reports that use the data available in a Results Table, in file information, and in quantitative peak review windows.
- Provides a variety of reports that present MS/MS library search results. The user can configure the Reporter software to search against any MS/MS spectral library that uses the Analyst[®] (mdb) format.
- Uses Microsoft Word templates to provide the format information needed when generating reports. These templates can be created or modified to provide customized report formats. Refer to the Help for information about creating or editing the Report Template Editor.
- Contains a blank starting template that can be used in the Analyst[®] software Reporter editing environment to design report templates to meet most reporting requirements.
- Automates report generation through the use of the Autoquan Reporter batch script.
- Automatically prints, exports to Adobe Portable Document Format (pdf), and delivers results by e-mail.
- Attaches processing scripts to report templates to expand both the content and automation level for various workflow requirements.
- Generates reports from custom software applications that use the available Analyst[®] software programming libraries.

Reporter software can be used as follows:

- Within the Analyst[®] software to manually generate a report or set of reports.
- By a batch script to automate report generation within a batch. Users can generate reports on a sample-by-sample basis, either during or after batch acquisition.
- By applications that do not use the Analyst[®] software.

Analyst Reporter User Interface

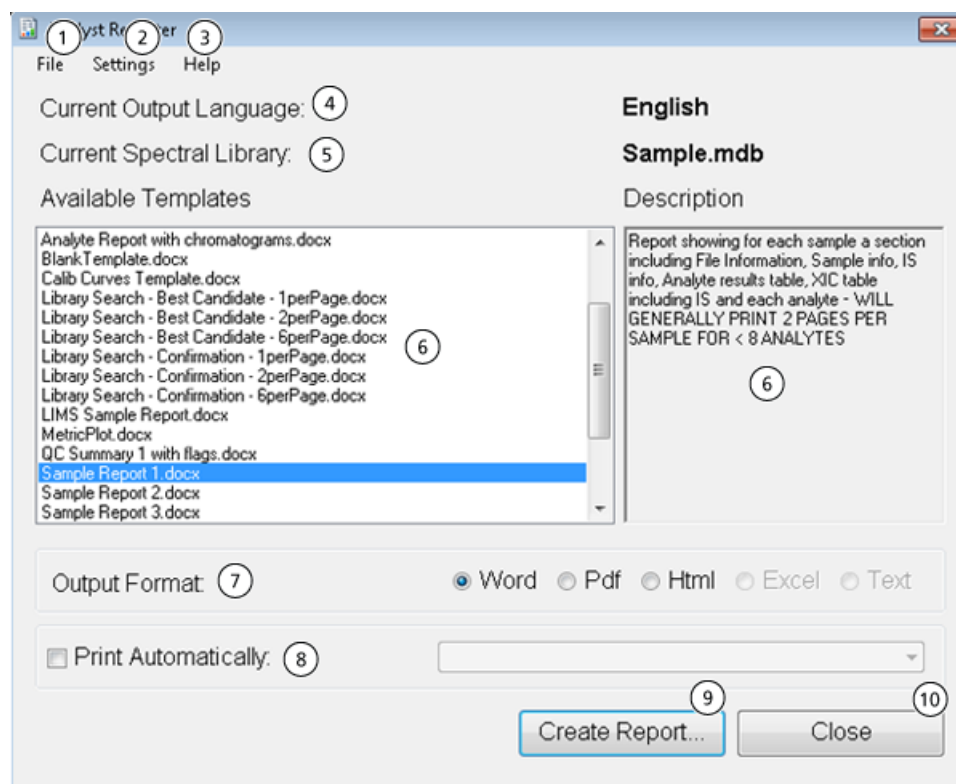Figure 11-1 Analyst Reporter



Item	Option	Description
1	File > Exit	Exits the program and releases all resources.
2	Settings > Select Output Language	Sets the language dictionary that will be used to replace language tags within a report template. Templates that contain language tags can be used to generate reports in any language. The language tags are replaced with text from a matching tag in the dictionary file for the selected language. These dictionary files are contained in the folder: C:\Program Files\AB SCIEX\AnalystReporter\Resources\Languages on Windows 7, 32-bit operating systems or C:\Program Files (x86)\AB SCIEX\AnalystReporter\Resources\Languages on Windows 7, 64-bit or Windows 10, 64-bit operating systems.
2	Settings > Select Library	Browse to a spectral library. This library will be used for matching and scoring MS/MS data from Results Tables that contain data from information dependent acquisition (IDA) triggered MS/MS scan types.
2	Settings > Select Template Folder	Sets the folder from which the available templates will be read. To return to the default template folder, select the Default option.
3	Help > About	Shows information about the version of Reporter software currently installed.
4	Current Output Language	Shows the currently selected language dictionary used for replacing language tags within a report template. The language dictionary can be selected using Settings > Select Output Language .
5	Current Spectral Library	Shows the currently selected spectral library. The spectral library can be selected using Settings > Select Library .
6	Available Templates and Description	Shows a list of available report templates. Selecting a template will show a description of the template. To change the folder from which available templates are read, select Settings > Select Template Folder > Browse .

Item	Option	Description
7	Output Format	<p>Shows the output formats that are supported by the Reporter software. Only formats that are compatible with the selected report template are enabled.</p> <ul style="list-style-type: none">• Word: Microsoft Word document (docx) is produced. This document can be viewed by Microsoft Word 2010 and above.• PDF: A report is created directly in PDF format.• HTML: Microsoft Word is used to generate an HTML file. Associated image files are stored in a folder with the same name as the HTML file.• Excel: A plain text file (csv) is produced. Report templates that contain values separated by commas can be opened in Microsoft Excel, where each value will be shown in a separate cell. Only templates that are specifically marked as text-compatible can be used for this output format.• Text: A plain text document (txt) is produced. Only templates that are specifically marked as text compatible can be used for this output format.
8	Print Automatically	After the report has been created it is printed to the selected printer. Select any available printer.
9	Create Report	Creates the report in the selected output format using the selected report template.
10	Close	Exits the program and releases all resources.

Generate Reports

The Reporter software extracts numerical data from the Results Table and sample and graphical information from the wiff file.

Select a template in the Available Template field.

Tip! For reports that can be generated on a sample-by-sample basis, it may be more efficient to generate the reports automatically using a batch script during acquisition to avoid long processing times at the end of the acquisition. For more information about batch scripts, refer to the *Scripts User Guide*.

1. Open a **Results Table**.

2. Under **Companion Software**, double-click **Reporter**.
3. In the **Analyst Reporter** dialog, in the **Available Templates** field, select the applicable report template.
4. Click the **PDF** output format.

The Word option is pre-selected and the report is automatically saved in the current project Results folder. If this option is not selected, then the report is created and opened in Word or printed as selected, but the report is not saved. This lets the user edit the report in Word prior to saving the original report.

5. Select either one document containing all samples or multiple documents with one sample in each.
6. (Optional) Select the **Print Automatically** check box to print the reports automatically on a pre-selected printer.

The Default Printer set in Windows is used unless a different printer is selected. The Reporter tool retains the selected printer between operations. If the printer is set to a PDF printer driver, then the Reporter generates PDF file versions of the created reports automatically.

7. Click **Create Report**.

The screen shows various progress indicators as the tool opens the template and populates it with data from the Results Table. Some reports can take seconds to generate, others can take longer. A large data set with many MRM transitions or a large number of graphics might result in reports of several hundred pages that take hours to generate.

CAUTION: Potential System Damage. If the HPLC system connected to the mass spectrometer are not controlled by the software, then do not leave the mass spectrometer unattended while in operation. The liquid stream from the HPLC system can flood the ion source when the mass spectrometer goes into Standby mode.


Set Queue Options

The queue goes one-by-one through the list, acquiring each sample with the selected acquisition method. After all of the samples have been acquired, the queue stops and the mass spectrometer goes into Standby mode. In Standby mode, the LC pumps and some instrument voltages are turned off.

The user can change the length of time the queue runs after the last acquisition has finished, before the Analyst[®] software puts the mass spectrometer into Standby mode. For information about the other fields in the Queue Options dialog, refer to the Help.

1. On the Navigation bar, click **Configure**.
2. Click **Tools > Settings > Queue Options**.

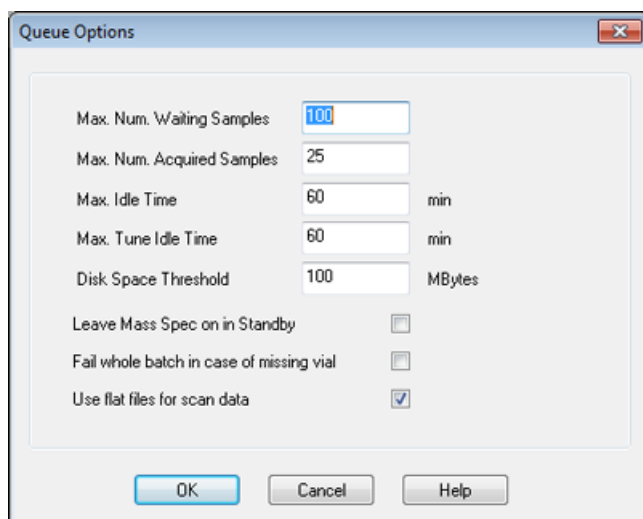Figure 12-1 Queue Options Dialog



3. In the **Max. Num. Waiting Samples** field, set the maximum number of samples to a value that is greater than the number of samples that will be submitted to the queue
4. In the **Max. Idle Time** field, type the length of time the queue will wait after acquisition is completed before going to Standby mode. The preset value is 60 minutes.

If gas cylinders are used, then adjust this time to make sure that the gas in the cylinders is not depleted.

If an LC method is used, then before the run is started, make sure that there is enough solvent in the reservoirs for all of the sample runs at the primary flow rate and the maximum idle time.

5. Select the **Leave Mass Spec on in Standby** check box to keep the mass spectrometer running after analysis has been completed. This feature allows the heaters and gases to continue running, even after devices have entered Idle state, so that the ion source and entrance to the mass spectrometer are kept free of contaminants.
6. Select the **Fail Whole Batch in Case of Missing Vial** check box to fail the entire batch when a missing vial is encountered. If this option is not selected, then only the current sample will fail and the queue will continue to the next sample.

Create and Submit a Batch

Use this workflow to create a batch. In this example, use the MRM scan type that was created previously. Go through the workflow twice more for practice, once using the Q1MS method and the second time using the Q1MI method.

Add Sets and Samples to a Batch

A set can consist of a single sample or multiple samples.

Note: For more information about adding quantitation information to a batch, refer to the *Advanced User Guide*.

1. On the Navigation bar, under **Acquire**, double-click **Build Acquisition Batch**.

Figure 12-2 Batch Editor Dialog

The Batch Editor Dialog box has a tabbed interface with 'Sample', 'Locations', 'Quantitation', and 'Submit' tabs. The 'Sample' tab is active. It contains a 'Select Method for Sample Set' section with a 'Set' dropdown menu (showing 'SET1'), 'Add Set', 'Remove Set', 'Add Samples', and 'Del Samples' buttons. To the right are 'Quantitation' and 'Acquisition' sections, each with a dropdown menu (both showing 'none') and buttons for 'Quick Quant' and 'Method Editor'. Below these is a 'Batch Script' text field and a 'Select Script' button. At the bottom is a table with columns: Sample Name, Rack Code, Rack Position, Plate Code, Plate Position, Vial Position, Data File, and Inj. Volume (µl).

2. In the Sample tab, in the **Set** list, type a name.
3. Click **Add Set**.
4. Click **Add Samples** to add samples to the new set.

Figure 12-3 Add Sample Dialog

The Add Sample Dialog box has a title bar with a close button. It contains three main sections: 'Sample name' with a 'Prefix' text field (containing 'Sample') and a 'Sample number' checkbox (checked) with a 'Number of digits' text field (containing '3'); 'Data file' with a 'Prefix' text field (containing 'Data'), a 'Set name' checkbox (checked), an 'Auto Increment' checkbox (checked), a 'Sub Folder' text field, and a 'Browse' button; and 'New samples' with a 'Number' text field (containing '1'). At the bottom are 'OK', 'Cancel', and 'Help' buttons.

5. In the **Sample name** section, in the **Prefix** field, type a name for the samples in this set.
6. To add incremental numbering to the end of the sample name, select the **Sample number** check box.
7. If the **Sample number** check box is selected, then in the **Number of digits** field, type the number of digits to include in the sample name.

For example, if 3 is typed, then the sample names would be samplename001, samplename002, and samplename003.

8. In the **Data file** section, in the **Prefix** field, type a name for the data file that will store the sample information.
9. Select the **Set name** check box to use the set name as part of the data file name.
10. Select the **Auto Increment** check box to increment the data file names automatically.

Note: The data for each sample can be stored in the same or a separate data file. The names of the data file will have numerical suffixes starting from 1.

11. Type a name in the **Sub Folder** field.

The folder is stored in the Data folder for the current project. If the **Sub Folder** field is left blank, then the data file is stored in the **Data** folder and a subfolder is not created.

12. In the **New samples** section, in the **Number** field, type the number of new samples to add.
13. Click **OK**.

The sample table fills with the sample names and data file names.

Tip! Fill Down and **Auto Increment** options are available in the right-click menu after a single column heading or several rows in a column are selected.

14. On the Sample tab, in the **Acquisition** section, select a method from the list.

Depending on how the system is set up, specific information for the autosampler must be entered. Even if the injection volume is set in the method, the user can change the injection volume for one or more samples by changing the value in the injection volume column.

Note: To use different methods for some of the samples in this set, select the **Use Multiple Methods** check box. The **Acquisition Method** column is shown in the **Sample** table. Select the acquisition method for each sample in this column.

15. To change the injection volumes from the volumes listed in the method, in the **Inj. Volume (μL)** column, type the injection volume for each sample.
16. To set sample locations, do one of the following:

- [Set Sample Locations in the Batch Editor on page 140](#)
 - [Select Vial Positions Using the Locations Tab \(Optional\) on page 141](#)
17. (Optional) To define quantitation details prior to submitting the batch, refer to [Set Quantitation Details in the Batch Editor \(Optional\) on page 142](#).
 18. Click the **Submit** tab.

Note: The order of samples can be edited before the samples are submitted to the queue. To change the order of samples, on the **Submit** tab, double-click any of the numbers at the far left of the table (a very faint square box is shown), and then drag them to the new location.

19. If the **Submit Status** section contains a message about the status of the batch, then do one of the following:
 - If the message indicates that the batch is ready for submission, then proceed to step [20](#).
 - If the message indicates that the batch is not ready for submission, then make the changes as indicated by the message.
20. After confirming that all of the batch information is correct, click **Submit**.

The batch is submitted to the queue and can be viewed in the Queue Manager.

21. Save the file.

Equilibrate the System

Equilibrate the system before submitting a batch. Equilibration warms up and prepares the mass spectrometer for the next sample or batch.

1. Click .

The Equilibrate dialog opens.

2. Select the acquisition method used for the submitted batch.
3. Type the equilibration time in the **Time (min)** field, in minutes.
4. Select **OK**. The system starts the equilibration.

After the equilibration is completed, the system will change to Ready status.

Tip! If the equilibration does not finish as expected, or if the system status does not change to Ready after the equilibration finishes, then make sure that

- The activated hardware profile is appropriate for the acquisition method.
 - The HPLC system is turned on.
 - The HPLC system has correctly established communication with the software.
-

Submit a Sample or Set of Samples

Note: Run the sample again in the event of an abnormal termination during sample acquisition. If the abnormal termination is caused by a power failure, then the temperature of the autosampler tray is not maintained and sample integrity might be compromised.

1. Select one sample or a set of samples.
2. Click the **Submit** tab in the **Batch Editor**.
3. If the **Submit Status** section contains a message about the status of the batch, then do one of the following:
 - If the message indicates that the batch is ready for submission, then proceed to the next step.
 - If the message indicates that the batch is not ready for submission, then make the changes as indicated by the message.
4. Click **Submit**.

Change Sample Order


The order of the samples can be edited before the samples are submitted to the **Queue**.

- On the **Submit** tab, double-click any of the numbers at the far left of the table (a very faint square box is visible), and then drag them to the new location.

Acquire Data


The system must not be in Tune and Calibrate mode when sample acquisition is started. Also, if the system has been previously run that day and has not yet been set to Standby mode, then sample acquisition will start automatically.

1. Make sure that the column oven temperature is reached.

2. Make sure that the  icon is not pressed in.
3. On the Navigation bar, click **Acquire**.
4. Click **View > Sample Queue**.

The Queue Manager opens with all submitted samples.

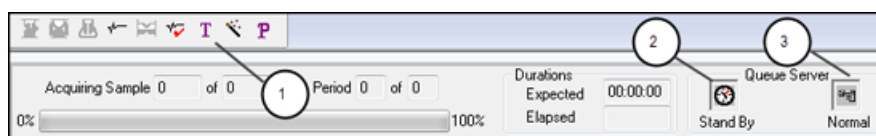Figure 12-4 Queue Manager



Item	Description
1	The Reserve Instrument for Tuning icon should not be pressed in.
2	Queue status should be in Ready mode.
3	Queue Server should be in Normal. Refer to Queue States on page 145 .

5. Click **Acquire > Start Sample**.

Set Sample Locations in the Batch Editor

If an autosampler is used in the acquisition method, then the vial positions of the samples must be defined in the acquisition batch. Define the location in the **Sample** tab or in the **Locations** tab. For more information about creating batches, refer to [Add Sets and Samples to a Batch on page 135](#).

1. In the **Sample** tab, from the **Set** list, select the set.
2. For each sample in the set, do the following if applicable:
 - In the **Rack Code** column, select the rack type.
 - In the **Rack Position** column, select the position of the rack in the autosampler.
 - In the **Plate Code** column, select the plate type.
 - In the **Plate Position** column, select the position of the plate on the rack.
 - In the **Vial Position** column, type the position of the vial in the plate or tray.
3. Save the file.

Select Vial Positions Using the Locations Tab (Optional)

1. Click the **Locations** tab in the **Batch Editor**.
2. Select the set from the **Set** list.
3. Select the autosampler from the **Autosampler** list.
4. In the space associated with the rack, right-click and then select the rack type.

The plates or trays are shown in the rack.

5. Double-click in the white space labeled rack type. A visual sample rack layout is shown.


The appropriate number of rack spaces for the autosampler is shown in the graphic rack view.

6. Double-click one of the rectangles.

The circles depicting the wells or vials for the plate or tray are shown.

Tip! To see the corresponding vial number in the graphical representation, move the cursor over the sample position. Use this information to confirm that the vial positions in the software match the vial positions in the autosampler.

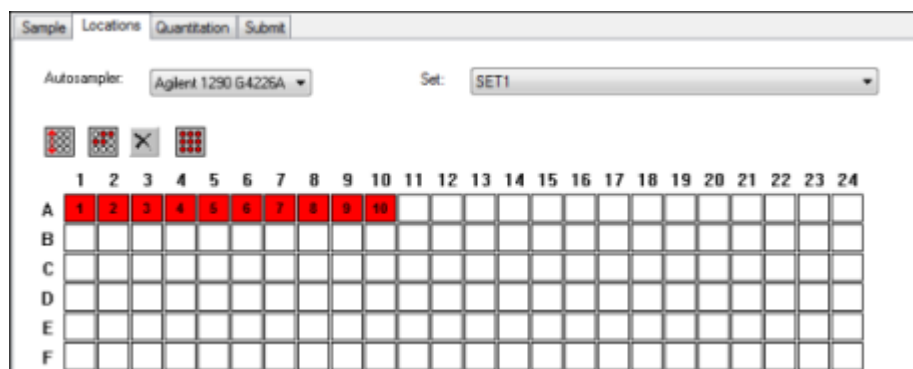Figure 12-5 Locations Tab



Note: Depending on the autosampler being used, it might not be necessary to type details in additional columns.

7. To select whether samples are marked by row or column, click the **Row/Column selection** selector button.

If the button shows a red horizontal line, then the **Batch Editor** marks the samples by row. If the button shows a red vertical line, then the **Batch Editor** marks the samples by column.

- Click the sample wells or vials in the order to be analyzed.

Tip! Click a selected well or vial again to clear it.

Tip! To fill in the samples automatically, press the **Shift** key while clicking the first and last vial within a set. To perform multiple injections from the same vial, press the **Ctrl** key while clicking the vial location. The red circle changes to a green circle.

Set Quantitation Details in the Batch Editor (Optional)

If a quantitation method is used with a batch and if the user does not want to select quantitation details after acquisition, then the quantitation details (sample type, sample concentration) must be defined before the batch is submitted.

The appropriate Internal Standard and Standard columns are shown in the Quantitation tab according to the quantitation method selected in the Sample tab.

- With a batch file open in the **Batch Editor** window, click the **Quantitation** tab.
- Select the set containing the samples.
- Select a **Quant Type** for all the samples from the list in the cell.
- If applicable, type the analyte concentration in the **Analyte** column.
- If applicable, type the internal standard concentration in the **Internal Standard** column.
- Repeat this procedure for each set in the batch.

Stop Sample Acquisition

When a sample acquisition is stopped, the current scan finishes before the acquisition is stopped.

- In the **Queue Manager**, click the sample in the queue after the point where acquisition should stop.
- On the Navigation bar, click **Acquire**.
- Click **Acquire > Stop Sample**.

The queue stops after the current scan in the selected sample is complete. The sample status in the **Queue Manager (Local)** window changes to **Terminated**, and all other samples following in the queue are **Waiting**.

- To continue processing the batch, click **Acquire > Start Sample**.

Import Batch Files

Users can import a text file containing batch information instead of creating a batch in the **Batch Editor**. If all the sample details are in a spreadsheet, then it is faster to rearrange and import the data in the spreadsheet than to manually type the data into the **Batch Editor**.

Before importing batch information from a text file, make sure that the data in the file is organized and formatted correctly. In particular, the column headings in the spreadsheet must match the **Batch Editor** column headings.

Build a Batch as a Text File

Prerequisites

Make sure that the active hardware profile includes all of the devices to be used to acquire the samples.

To make sure that the text file includes the proper headings, create a batch using the Batch Editor, export it as a text file, type the appropriate values in a spreadsheet editor, and then import the file back into the Batch Editor. Users can export a batch only if it contains at least one set with at least one sample. The saved text file can be used again later as a template.

1. In the **Batch Editor**, create a single-set, single-sample batch.
2. Click **File > Export**.

The **Save As** dialog opens.

3. Type a name for the text file in the **File name** field and then click **Save**.
4. Open the text file in a spreadsheet program such as Microsoft Excel.
5. Type, or copy and paste, the details for the samples: one sample per row, with the details under the appropriate headings.

Note: Do not delete any of the columns. The columns in the spreadsheet must match the columns in the Batch Editor.

6. Save the modified text file as a .txt or .csv file and then close the spreadsheet program.

The text file is now ready to be imported into the Batch Editor.

Import a Batch from a Text File

1. In the **Batch Editor**, in the **Sample** tab, right-click, and then click **Import From > File**.

The **Open** dialog opens.

2. Click the required text file and then click **Open**.

If an autosampler is used, then the **Select Autosampler** dialog opens.

Note: If the saved text file is not visible in the **Files of type** list, then select **Microsoft Text Driver (*.txt; *.csv)**. Files with the extension .txt are shown in the field.

3. In the autosampler list, select the autosampler and then click **OK**.


The sample table fills with the details from the text file.

4. Submit the batch.

Batch Editor Right-click Menu

Right-click in the **Batch Editor** table to access the options.

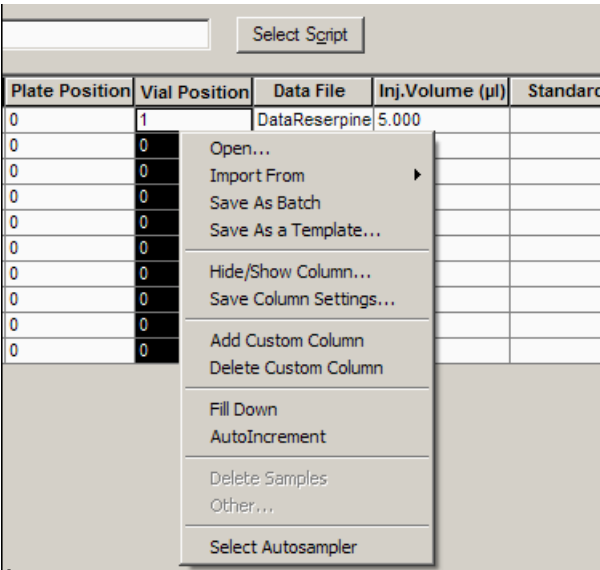Figure 12-6 Batch Right-Click Menu



Menu	Function
Open	Opens a batch file.
Import From	Imports a file.
Save As Batch	Saves the batch with a different name.
Save As a Template	Saves the batch as a template.
Hide/Show Column	Hides or shows a column.

Menu	Function
Save Column Settings	Saves the batch column settings.
Add Custom Column	Adds a custom column.
Delete Custom Column	Deletes a custom column.
Fill Down	Copies the same data into the selected cells.
AutoIncrement	Automatically increments data into the selected cells.
Delete Samples	Deletes the selected row.
Select Autosampler	Selects an autosampler.

Queue States and Device Status

The **Queue Manager** shows queue, batch, and sample status. Detailed information about a particular sample in the queue can also be viewed.

Tip! Click  to view the queue.

Queue States

The current state of the queue is indicated in the **Queue Server**.

Figure 12-7 Queue Server Indicator Showing Normal Mode


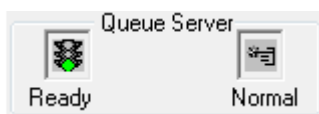


Figure 12-8 Queue Server Indicator Showing Tune Mode

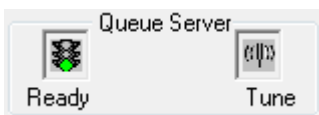

The first icon indicates the queue state. The second icon indicates whether the queue is in **Tune** mode (for tuning) or **Normal** mode (for running samples). [Table 12-1](#) describes the icons and queue states.

Table 12-1 Queue States

Icons	State	Definition
	Not Ready	The hardware profile is deactivated and the queue is not accepting any sample submissions.
	Stand By	The hardware profile has been activated, but all devices are idle. Pumps are not running and gases are turned off.
	Warming Up	The mass spectrometer and devices are equilibrating, columns are being conditioned, the autosampler needle is being washed, and column ovens are reaching temperature. The duration of equilibration is selected by the operator. From this state, the system can go to the Ready state.
	Ready	The system is ready to start running samples and the devices have been equilibrated and are ready to run. In this state, the queue can receive samples and will run after samples are submitted.
	Waiting	The system will automatically begin acquisition when the next sample is submitted.
	PreRun	The method is being downloaded to each device and device equilibration is occurring. This state occurs before the acquisition of each sample in a batch.
	Acquiring	The method is running and data acquisition is occurring.
	Paused	The system has been paused during acquisition.

View Instrument and Device Status Icons

Icons representing the mass spectrometer and each device in the active hardware configuration are shown on the status bar in the bottom right corner of the window. The user can view the detailed status of an LC pump to









determine whether the LC pump pressure is appropriate or view the detailed status of the mass spectrometer to confirm the temperature of the ion source.

Note: For each status, the background color can be red. A red background indicates that the device encountered an error while in that state.

- On the status bar, double-click the icon for the device or mass spectrometer.

The Instrument Status dialog opens.


Table 12-2 Instrument and Device Status Icons

Status	Icon	Background Color	Description
Idle		Green or yellow	The device is not running. If the background color is yellow, then the device should be equilibrated before it is ready to run. If the background color is green, the device is ready to run.
Equilibrating		Green or yellow	The device is equilibrating.
Waiting		Green	The device is waiting for a command from the software or another device, or for some action by the operator.
Running		Green	The device is running a batch.
Aborting		Green	The device is aborting a run.
Downloading		Green	A method is being transferred to the device.
Ready		Green	The device is not running, but is ready to run.
Error		Red	The device has encountered an error that should be investigated.

Queue Right-click Menu

Right-click in the Queue table to access the options.

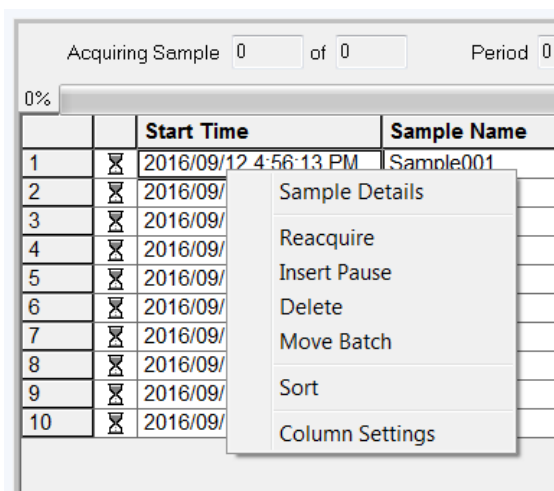Figure 12-9 Queue Manager Right-Click Menu



Menu	Function
Sample Details	Opens the Sample Details dialog.
Reacquire	Acquires a sample again.
Insert Pause	Inserts a pause, in seconds, between two samples.
Delete	Deletes either the batch or the selected samples.
Move Batch	Moves the batch within the queue.
Sort	Sorts on the preselected column.
Column Settings	Changes the column settings.

Service and Maintenance Information

13

Regularly clean and maintain the system for optimal performance.



WARNING! Electrical Shock Hazard. Do not remove the covers. Removing the covers might cause injury or malfunctioning of the system. The covers need not be removed for routine maintenance, inspection, or adjustment. Contact a SCIEX Field Service Employee (FSE) for repairs that require the covers to be removed.



WARNING! Radiation Hazard, Biohazard, or Toxic Chemical Hazard. Determine whether decontamination is required prior to cleaning or maintenance. The customer must decontaminate the system prior to cleaning or maintenance if radioactive materials, biological agents, or toxic chemicals have been used with the system.

Recommended Maintenance Schedule

[Table 13-1](#) provides a recommended schedule for cleaning and maintaining the system.

Tip! Perform maintenance tasks regularly to make sure that the mass spectrometer is performing optimally.

For information on maintaining the ion source, refer to the ion source *Operator Guide*.

Contact a Qualified Maintenance Person (QMP) to order consumable parts. Contact a SCIEX Field Service Employee (FSE) for maintenance service and support.

Table 13-1 Maintenance Tasks

Component	Frequency	Task	For more information...
System			
Tubing	Daily	Inspect the tubing and fittings to make sure that they are securely connected, and that there are no leaks.	Refer to Chemical Precautions on page 11 .

Service and Maintenance Information

Table 13-1 Maintenance Tasks (continued)

Component	Frequency	Task	For more information...
Mass Spectrometer			
Curtain plate	Daily	Clean	Refer to Clean the Curtain Plate on page 155 .
Orifice plate (front)	Daily	Clean	Refer to Clean the Front of the Orifice Plate on page 156 .
Roughing pump oil	Weekly	Inspect the level	Refer to Inspect the Roughing Pump Oil Level on page 157 .
Mass spectrometer air filter	Every 6 months	Replace	Contact the local QMP or FSE.
Roughing pump oil	Annually	Replace	Contact the local QMP or FSE.
Instrument surfaces	As needed	Clean	Refer to Clean the Surfaces on page 151 .
Source exhaust drain bottle	As needed	Empty	Refer to Empty the Source Exhaust Drain Bottle on page 157 .
Orifice plate (front and rear)	As needed	Clean	Contact the local QMP or FSE.
QJet [®] ion guide and IQ0 lens	As needed	Clean	Contact the local QMP or FSE.
Q0 rod set and IQ1 lens	As needed	Clean	Contact the local QMP or FSE.
Roughing pump oil	As needed	Refill	Contact the local QMP or FSE.
Interface heater	As needed	Replace	Contact the local QMP or FSE.
Ion Source			
TurbolonSpray [®] and APCI electrodes	As needed	Inspect and replace	Refer to the ion source <i>Operator Guide</i> .
Corona discharge needle	As needed	Replace	Refer to the ion source <i>Operator Guide</i> .
Turbo heater	As needed	Replace	Contact the local QMP or FSE.
Sample tubing	As needed	Replace	Refer to the ion source <i>Operator Guide</i> .

For “As needed” tasks, follow these guidelines:

- Clean the mass spectrometer surfaces after a spill or when they become dirty.
- Empty the drain bottle before it becomes full.
- Clean the orifice plate, QJet[®] ion guide, and Q0 region if system sensitivity degrades.

Tip! Clean the Q0 region regularly to minimize the impact of charging (a significant loss of sensitivity of the ions of interest over a short period of time) on the quadrupoles and lenses. Contact a QMP or FSE.

- Clean the QJet[®] ion guide and Q0 region if system sensitivity degrades.

Tip! Clean the Q0 region regularly to minimize the impact of charging (a significant loss of sensitivity of the ions of interest over a short period of time) on the quadrupoles and lenses. Contact a QMP or FSE.

- Refill the roughing pump oil when it falls below the minimum oil level.

Clean the Surfaces

Clean the external surfaces of the mass spectrometer after a spill or when they become dirty.

CAUTION: Potential System Damage. Use only the recommended cleaning methods and materials to avoid damaging the equipment.

1. Wipe the external surfaces with a soft cloth dampened with warm, soapy water.
2. Wipe the external surfaces with a soft cloth moistened with water to remove any soap residue.

Clean the Front-End

The following warning applies to all of the procedures in this section:



WARNING! Hot Surface Hazard. Let the ion source cool for at least 30 minutes before starting any maintenance procedures. Surfaces of the ion source and the vacuum interface components become hot during operation.

Clean the mass spectrometer front-end using the routine cleaning method, to:

- Minimize unscheduled system downtime.
- Maintain optimum sensitivity.
- Avoid more extensive cleaning that requires a service visit.

Service and Maintenance Information

When contamination occurs, perform an initial routine cleaning. Clean up to and including the front of the orifice plate. If routine cleaning does not resolve issues with sensitivity, then a full cleaning might be necessary. Contact the local QMP or FSE.

This section provides instructions for performing routine cleaning without breaking vacuum.

Note: Follow all of the applicable local regulations. For health and safety guidelines, refer to [Chemical Precautions on page 11](#).

Symptoms of Contamination

The system might be contaminated if any one of the following is observed:

- Significant loss in sensitivity
- Increased background noise
- Additional peaks that are not part of the sample are shown in full scan or survey scan methods

If any of these issues are observed, then clean the mass spectrometer front-end.

Required Materials

Note: U.S. customers can call 877-740-2129 for ordering information and inquiries. International customers can visit sciex.com/contact-us.

- Powder-free gloves (nitrile or neoprene recommended)
- Safety glasses
- Laboratory coat
- Fresh, high-quality (pure) water (at least 18 MΩ de-ionized [DI] water or ultra-pure HPLC-grade water). Old water can contain contaminants that can further contaminate the mass spectrometer.
- MS-grade methanol, isopropanol (2-propanol), or acetonitrile
- Cleaning solution. Use one of:
 - 100% methanol
 - 100% isopropanol
 - 1:1 acetonitrile:water solution (freshly prepared)
 - 1:1 acetonitrile:water with 0.1% acetic acid solution (freshly prepared)
- Clean 1 L or 500 mL glass beaker to prepare cleaning solutions
- 1 L beaker to catch used solvent
- Organic waste container

- Lint-free wipes. Refer to [Tools and Supplies Available from the Manufacturer on page 153](#).
- (Optional) Polyester (poly) swabs

Tools and Supplies Available from the Manufacturer

Description	Part Number
Small poly swab (thermally bonded). Also available in the Cleaning kit.	1017396
Lint-free wipe (11 cm x 21 cm, 4.3 inches x 8.3 inches). Also available in the Cleaning kit.	018027
Cleaning kit. Contains the small poly swab, lint-free wipes, Q0 cleaning tool, straight QJet [®] ion guide cleaning brush, and Alconox packets.	5020761

Cleaning Best Practices



WARNING! Hot Surface Hazard. Let the ion source cool for at least 30 minutes before starting any maintenance procedures. Surfaces of the ion source and the vacuum interface components become hot during operation.



WARNING! Toxic Chemical Hazard. Refer to the chemical product *Safety Data Sheets* and follow all of the recommended safety procedures when handling, storing, and disposing of chemicals. For health and safety precautions, refer to the *System User Guide*.



WARNING! Radiation Hazard, Biohazard, or Toxic Chemical Hazard. Determine whether decontamination is required prior to cleaning or maintenance. The customer must decontaminate the system prior to cleaning or maintenance if radioactive materials, biological agents, or toxic chemicals have been used with the system.



WARNING! Environmental Hazard. Do not dispose of system components in municipal waste. Follow local regulations when disposing of components.

- Allow the ion source to cool before removing it.
- Always wear clean, powder-free gloves (nitrile or neoprene recommended) for the cleaning procedures.
- After cleaning the mass spectrometer components, and before reassembling them, put on a new, clean pair of gloves.

Service and Maintenance Information

- Do not use cleaning supplies other than those specified in this procedure.
- If possible, prepare cleaning solutions just before cleaning.
- Prepare and store all of the organic solutions and organic-containing solutions in very clean glassware only. Never use plastic bottles. Contaminants can leach from these bottles and further contaminate the mass spectrometer.
- To avoid contaminating the cleaning solution, pour the solution on the wipe or swab.
- Allow only the center area of the wipe to contact the mass spectrometer surface. Cut edges can leave fibers behind.

Tip! Wrap the wipe around a thermally-bonded poly swab.

Figure 13-1 Example: Folding the Wipe



- To avoid cross-contamination, discard the wipe or swab after it has touched the surface once.
- Larger parts of the vacuum interface, such as the curtain plate, might require several cleanings, using multiple wipes.
- Only dampen the wipe or swab slightly when applying water or cleaning solution. Water, more often than organic solvents, might cause the wipe to deteriorate, leaving residue on the mass spectrometer.
- Do not rub the wipe across the aperture. Wipe around the aperture to prevent fibers from the wipes from entering the mass spectrometer.
- Do not insert the brush in the aperture on the curtain plate or orifice plate.

Prepare the Mass Spectrometer


Note: Mass spectrometers with a NanoSpray[®] ion source might require a full cleaning for best results. Contact a local QMP or an FSE.



WARNING! Hot Surface Hazard. Let the ion source cool for at least 30 minutes before starting any maintenance procedures. Surfaces of the ion source and the vacuum interface components become hot during operation.

CAUTION: Potential System Damage. Do not drop anything into the source drain when the ion source is removed.

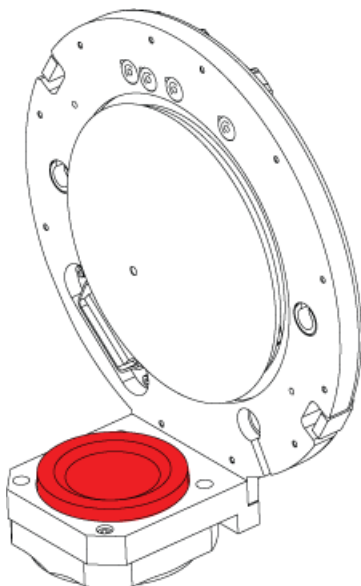Figure 13-2 Source Drain on the Vacuum Interface



1. Deactivate the hardware profile.
2. Remove the ion source. Refer to the ion source *Operator Guide*.

When the ion source is not in use, store it to protect it from damage and to maintain operating integrity.

Clean the Curtain Plate

CAUTION: Potential System Damage. Do not rest the curtain plate or orifice plate on the aperture tip. Make sure that the conical side of the curtain plate faces up.

CAUTION: Potential System Damage. Do not insert a wire or metal brush into the aperture on the curtain plate, orifice plate, or interface heater to avoid damaging the aperture.

1. Pull the curtain plate off of the vacuum interface and then put it, conical side up, on a clean, stable surface.

The curtain plate is held in place by three retaining ball catches mounted on the orifice plate.

Tip! If the curtain plate does not immediately separate from the orifice plate, then turn the curtain plate slightly (less than 90 degrees) to release the ball spring latches.

Service and Maintenance Information

2. Dampen a lint-free wipe with pure water and then clean both sides of the curtain plate.

Note: Use multiple wipes, as required.

3. Repeat step 2 using the cleaning solution.
4. Using a dampened wipe or small poly swab, clean the aperture.
5. Wait until the curtain plate is dry.
6. Inspect the curtain plate for solvent stains or lint, removing any residue with a clean, slightly damp, lint-free wipe.

Note: Persistent spotting or filming is an indicator of contaminated solvent.

Clean the Front of the Orifice Plate

CAUTION: Potential System Damage. When cleaning the surface of the orifice plate, do not remove the interface heater. Frequent removal of the interface heater can result in damage to the interface heater. Surface cleaning of the interface heater is adequate for routine cleaning.

CAUTION: Potential System Damage. Do not insert a wire or metal brush into the aperture on the curtain plate, orifice plate, or interface heater to avoid damaging the aperture.

1. Dampen a lint-free wipe with water and then wipe the front of the orifice plate, including the interface heater.
2. Repeat step 1 using the cleaning solution.
3. Wait until the orifice plate is dry.
4. Inspect the orifice plate for solvent stains or lint, removing any residue with a clean, slightly damp, lint-free wipe.

Note: Persistent spotting or filming is an indicator of contaminated solvent.

Put the Mass Spectrometer Back in Service


1. Install the curtain plate on the mass spectrometer.
2. Install the ion source on the mass spectrometer. Refer to the ion source *Operator Guide*.
3. Activate the hardware profile. Refer to the *System User Guide*.

Inspect the Roughing Pump Oil Level

- Inspect the sight glass on the roughing pump to verify that the oil is above the minimum level.

If the oil is below the minimum level, then contact the qualified maintenance person (QMP) or SCIEX field service employee (FSE).

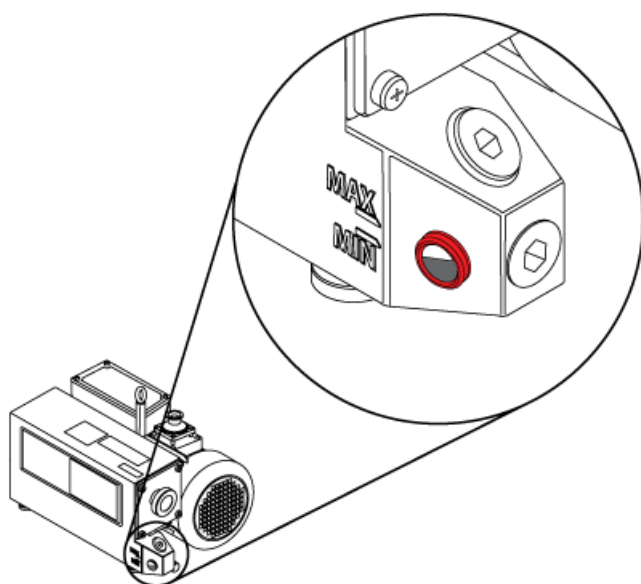Figure 13-3 Sight Glass



Empty the Source Exhaust Drain Bottle



WARNING! Radiation Hazard, Biohazard, or Toxic Chemical Hazard. Deposit hazardous materials in appropriately labeled waste containers and dispose of them according to local regulations.




WARNING! Radiation Hazard, Biohazard, or Toxic Chemical Hazard. Take care to vent exhaust gases to a dedicated laboratory fume hood or exhaust system and make sure that the ventilation tubing is secured with clamps. Make sure that the laboratory has appropriate air exchange for the work performed.

Service and Maintenance Information

Inspect the source exhaust drain bottle regularly, and empty it before it becomes full. Also inspect the bottle and its fitting for leaks, and tighten connections or replace components, if required. Follow the steps in this procedure to empty the bottle.

1. Remove the ion source. Refer to the ion source *Operator Guide*.
2. Loosen the clamps connecting the hoses to the cap of the source exhaust drain bottle.
3. Disconnect the hoses from the cap.
4. If applicable, lift the drain bottle out of the holder.
5. Remove the cap from the drain bottle.
6. Empty the drain bottle and then dispose of the waste according to laboratory procedures and local waste regulations.
7. Install the cap on the bottle and then put the bottle in the holder.
8. Attach the hoses to the cap and then secure them tightly with clamps.

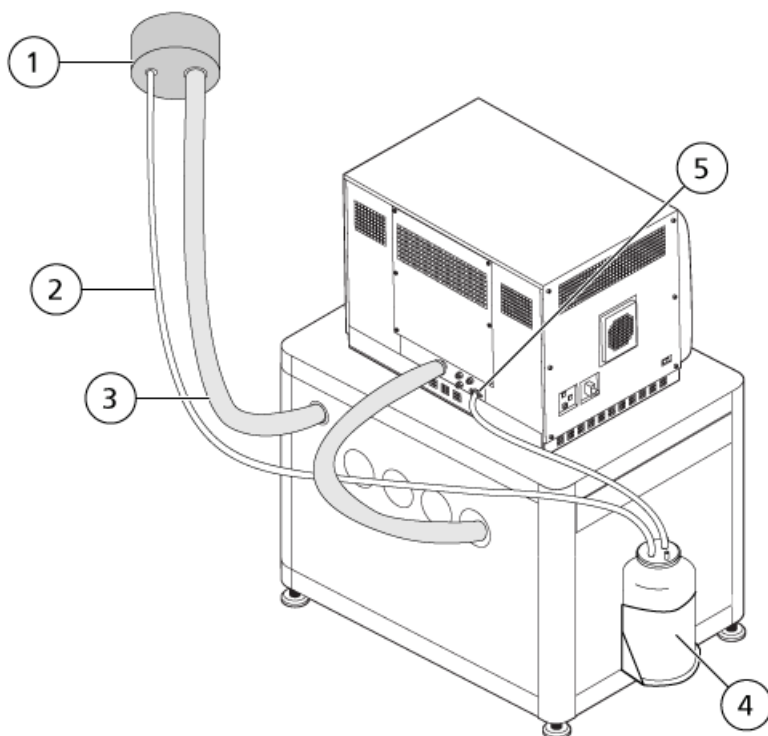Figure 13-4 Source Exhaust Drain Bottle



Item	Description
1	Connection to vent
2	Source exhaust drain tubing: 2.5 cm (1.0 inch) inside diameter (i.d.)
3	Roughing pump exhaust hose: 3.2 cm (1.25 inch) i.d.
4	Source exhaust drain bottle. Make sure that the bottle is secured to prevent spillage.
5	Connection to the mass spectrometer: 1.6 cm (0.625 inch) i.d.

Note: Source exhaust hose connections at the drain bottle, mass spectrometer, and the lab vent are secured with hose clamps.

Storage and Handling



WARNING! Environmental Hazard. Do not dispose of system components in municipal waste. Follow local regulations when disposing of components.

If the mass spectrometer must be stored for a long time or prepared for shipping, then contact a SCIEX FSE for decommissioning information. To disconnect power from the mass spectrometer, remove the mains supply connector from the AC mains supply.

Note: The ion source and mass spectrometer must be transported and stored between $-30\text{ }^{\circ}\text{C}$ to $+60\text{ }^{\circ}\text{C}$ ($-22\text{ }^{\circ}\text{F}$ to $140\text{ }^{\circ}\text{F}$). Store the system at an altitude not exceeding 2 000 m (6 562 feet) above sea level.

Mass Spectrometer Troubleshooting

14

This section contains information for troubleshooting basic system issues. Certain activities can only be carried out by a SCIEX -trained Qualified Maintenance Person (QMP) in the laboratory. For advanced troubleshooting, contact a SCIEX Field Service Employee (FSE).

Table 14-1 System Issues

Symptom	Possible Cause	Corrective Action
The QJet [®] ion guide is extremely dirty or frequently becomes dirty.	The Curtain Gas [™] flow rate is too low.	Verify the setting for the CUR parameter, and increase it, if applicable.
A system fault has occurred because the vacuum pressure is too high.	<ol style="list-style-type: none">1. The roughing pump oil level is too low.2. There is a leak.3. The wrong orifice plate is installed.	<ol style="list-style-type: none">1. Inspect the roughing pump oil level, and then contact the local QMP or an FSE to add oil.2. Inspect and repair leaks.3. Install the correct orifice plate.
A system fault has occurred because the QPS exciter module temperature is too high.	<ol style="list-style-type: none">1. The mass spectrometer air filter is blocked.2. Ambient temperature is too high.3. The coil box is not tuned.	Contact the local QMP or FSE.
The Analyst [®] software reports that the mass spectrometer is in Fault status because of the ion source.	<ol style="list-style-type: none">1. The probe is not installed.2. The probe is not connected securely.	<ol style="list-style-type: none">1. Confirm the fault in the Status panel of the device details page.2. Install the probe. Refer to the ion source <i>Operator Guide</i>.3. Remove and replace the probe. Tighten the retaining ring securely. Refer to the ion source <i>Operator Guide</i>.
The Analyst [®] software indicates that the APCI probe is in use, but the TurbolonSpray [®] probe is installed.	The F3 fuse is blown.	Contact an FSE.

Table 14-1 System Issues (continued)

Symptom	Possible Cause	Corrective Action
The spray is not uniform.	The electrode is blocked.	Clean or replace the electrode. Refer to the ion source <i>Operator Guide</i> .
Sensitivity is reduced.	<ol style="list-style-type: none"> 1. The ion source parameters are not optimized. 2. The mass spectrometer is not optimized. 3. The curtain plate is dirty. 4. The orifice plate is dirty. 5. The QJet[®] ion guide or IQ0 lens is dirty. 6. The Q0 region is dirty. 7. The syringe or sample line is leaking. 8. The sample has degraded or has a low concentration. 9. The probe is not installed properly. 10. The ion source is not installed properly, or it is faulty. 11. One or more of the O-rings on the vacuum interface is missing. 12. There is an issue with the LC system or connections. 	<ol style="list-style-type: none"> 1. Optimize the ion source parameters. Refer to the Analyst[®] software Help system. 2. Refer to Clean the Curtain Plate on page 155. 3. Refer to Clean the Front of the Orifice Plate on page 156 or contact the local QMP or FSE. 4. Clean the Q0 region. Contact the QMP or FSE. 5. Inspect the syringe or sample line for leaks, and repair any leaks found. Make sure that all fittings are the correct type and size. 6. Verify the sample concentration. Use a fresh sample. 7. Remove and install the probe. 8. Remove and install the ion source, making sure that the latches are properly secured. If this does not resolve the issue, then install and optimize an alternate ion source. 9. If the O-rings are on the ion source, then install them on the vacuum interface. If they are missing, then contact an FSE. 10. Troubleshoot the LC system.

Mass Spectrometer Troubleshooting

Table 14-1 System Issues (continued)

Symptom	Possible Cause	Corrective Action
The mass spectrometer performance has degraded.	<ol style="list-style-type: none">1. The probe is not optimized.2. The sample was not prepared correctly or the sample has degraded.3. The sample inlet fittings are leaking.	<ol style="list-style-type: none">1. Optimize the probe. Refer to the ion source <i>Operator Guide</i>.2. Confirm that the sample was prepared correctly.3. Verify that the fittings are the right size and type and make sure that they are tight. Do not overtighten the fittings. Replace the fittings if leaks continue.4. Install and optimize an alternate ion source.5. Contact an FSE if the issue persists.
Arcing or sparks occur.	The position of the corona discharge needle is incorrect.	Turn the corona discharge needle toward the curtain plate, and away from the stream of heater gas. Refer to the ion source <i>Operator Guide</i> .

For sales, technical assistance, or service, contact an FSE or visit the SCIEX Web site at sciex.com for contact information.

Operating Instructions — Manual Compound Optimization

A

The user must control the autosampler and injection valve manually, as these devices cannot be controlled through the system while it is in Tune and Calibrate mode.

Prerequisites

- The mass spectrometer is tuned and calibrated.
- Conditions for an LC separation are known.
- All of the required peripheral devices, including the syringe pump, if needed, and LC components are in the hardware profile.

Required Materials

To tune instrument parameters for particular compounds the following solutions are recommended. The mixture of four compounds is used for illustrating the steps of the procedure.

- Mobile phase: 1:1 acetonitrile:water + 2 mM ammonium acetate + 0.1% formic acid.
- LC pump and autosampler.
- Autosampler vials.
- Four-compound mix (50 ng/mL), consisting of reserpine, minoxidil, tolbutamide, and rescinnamine. Use a solution that is 49.9% acetonitrile, with 50% deionized water and 0.1% formic acid as diluent. Other compounds can be substituted as long as their molecular weight is known and the compound is reasonable to be ionized by a API (Atmospheric Pressure Ionization) source.

Table A-1 Compounds and Molecular Weights

Compound	<i>m/z</i>
Minoxidil	210.2
Tolbutamide	271.1
Reserpine	609.3
Rescinnamine	635.3

About Manual Compound Optimization

Manual compound optimization is used to optimize compound- and ion source-dependent parameters for an analyte. When the user manually optimizes for an analyte, an MS acquisition method is created in Tune and Calibrate mode. Depending on the method of sample introduction selected, add an LC method to the acquisition method so that infusion or LC can be used.

Optimizing to give the highest signal does not always give the highest signal-to-noise ratios. Noise can scale with signal for some parameters and should be verified during optimization if the goal is to obtain maximum signal-to-noise ratio.

When optimizing ion source-dependent parameters, introduce the sample at flow rates that will be used during sample analysis, using either flow injection analysis (FIA) or Tee-infusion as the method of sample introduction. The CAD gas is the only compound-dependent parameter that is shown on the Source/Gas tab and can be easily optimized while infusing the analyte.

Optimize the position of the ion source before optimizing the ion source-dependent parameters. Refer to the ion source *Operator Guide*.

About Scan Types

For this example, use the Q1 MS, Q1 MI, Product Ion, and MRM scan types. The Q1 MS scan type is used to confirm the presence of compounds of interest. The Q1 MI scan is used to optimize MS or pre-collision cell voltages. The Product Ion scan type is used to determine the product ions of each compound. The MRM scan type is used to optimize the collision energy (CE) and collision cell exit potential (CXP) for each product ion or fragment. Use the methods created in this section for quantitative or qualitative analysis.

Manually Optimize an Analyte

After the acquisition method is created, optimize compound-dependent parameters using the **Edit Ramp** function or by manually editing the parameters in the Tune Method Editor. Ion source-dependent parameters can be optimized only by manually adjusting the parameters in the Tune Method Editor. Depending on the scan type used, different parameters are available to optimize.

Follow the procedures in the order given:

1. [Confirm the Presence of Compounds on page 165](#)
2. [Optimize MS-Specific Parameters on page 166](#)
3. [Determine the Product Ions for Optimization on page 168](#)
4. [Determine the Product Ions for Optimization on page 168](#)
5. [Optimize Collision Cell Exit Potential for each Product Ion on page 169](#)

Confirm the Presence of Compounds

1. Create a project.
2. Activate the hardware profile.
3. Prepare the sample:
 - a. Aspirate the compound solution into a syringe and remove the air from the syringe.
 - b. Use the tubing with the special fitting to connect the syringe to the mass spectrometer.
 - c. Install the syringe in the integrated syringe pump.
4. Infuse the compound in solution at a rate of 5 $\mu\text{L}/\text{min}$ to 10 $\mu\text{L}/\text{min}$.
5. On the Navigation bar, under **Tune and Calibrate**, double-click **Manual Tuning**.
6. On the Syringe Pump Method Properties tab, type the parameter values shown in [Table A-2](#).

Table A-2 Syringe Pump Method Properties Tab

Parameter	Value
Syringe Diameter	Syringe dependent. 4.610 mm if a 1.0 mL syringe is used
Flow Rate	10
Unit	$\mu\text{L}/\text{min}$

Figure A-1 Syringe Pump Method Properties Tab

7. Click **Start Syringe Pump**.
8. Click **MS Method** from the method list.
9. Click **Start**.
10. Wait until an even TIC is shown on the left and peaks are shown on the right and then click **Stop**.

11. Select the **MCA** check box.
12. Type **10** in the **Cycles** field.
13. Click **Start**.
14. When ten scans are complete, the masses of the four compounds are shown as ions.

Tip! If the previous or expected ion chromatogram pane is hidden, then select the window to be shown from the **Window** menu.

Note: Ion intensities of the compounds can show large variances. To facilitate moving to a solution of higher or lower concentration as needed during the optimization, have several concentration levels prepared before beginning optimization.

15. Right-click the bottom-right spectral pane and then click **Open File**.
16. Find the compounds of interest and then write down the m/z values for the highest peaks. These values should be within 0.1 Da to 0.2 Da of the expected m/z . In the next procedure, use the m/z values.

Optimize MS-Specific Parameters

Declustering potential (DP) and entrance potential (EP) are pre-collision cell voltages. Optimization of these values involves gradually changing the voltage range while monitoring the signal intensity of the compound.

The DP is the difference between the orifice and the ground. The higher the potential difference, the greater the amount of declustering.

The DP parameter has a significant effect on analyte signal. Typical DP values range from 20 V to 150 V. A DP value that is too low will result in lower ion intensity and potential interferences from clusters. A DP value that is too high can cause fragmentation of the analyte in the source. Generally, the DP should be set to the value that provides the highest intensity.

The EP parameter controls the entrance potential, which guides and focuses the ions through the high-pressure Q0 region. It is typically set at 10 V for positive ions or –10 V for negative ions. The EP has a minor effect on compound optimization and can generally be left at default values without impact on analyte detection limits.

1. Return to the Tune Method Editor and change the method to the **Q1 Multiple Ions (Q1 MI)** scan type.
2. In the mass table, type the appropriate parameter values. Refer to [Table A-3](#).

Table A-3 Mass Table Parameters—Q1 Multiple Ions (Q1 MI)

Compound	Q1 Mass	Time
Reserpine	609.4	1
Minoxidil	210.2	1

Table A-3 Mass Table Parameters—Q1 Multiple Ions (Q1 MI) (continued)

Compound	Q1 Mass	Time
Tolbutamide	271.3	1
Rescinnamine	635.4	1

Start with reserpine for a simple case. Repeat the manual optimization process for the remaining compounds.

- Click **Edit Ramp**.
- Select **Declustering Potential (DP)** in the Ramp Parameter Settings dialog.

Note: Start with the DP parameter and then optimize the other parameters in the order that they are shown in the dialog. The parameters might not be optimized correctly if they are optimized out of order.

- Type the required **Start**, **Stop**, and **Step** values.

Tip! The existing values are good starting points. Use the **Edit Ramp** function to change these values to be more efficient.

- Click **OK**.
- Click **Start**.
- Right-click the bottom-right XIC pane and then click **Open File** to maximize the XIC view.
- Monitor the XICs.

Note: The value that gives the best signal per second for the ion of interest is the optimal value.

- Note the optimal value for the ion of interest.
- Move the cursor to the mass table, right-click, and then add the parameter just optimized.
A column is added to the table.
- Add the optimized value to the appropriate row.
- Repeat these steps for each mass in the acquisition method until there is a list of optimal values for all of the masses.
- Repeat these steps to optimize for the other MS-specific parameters.

Table A-4 MS-Specific Parameters

Parameter	Comment
DP	Set the DP to the value that provides the highest intensity.
EP	Rarely optimize this parameter because it has a smaller effect.

Determine the Product Ions for Optimization

The collision energy (CE) controls the amount of energy that the precursor ions receive as they are accelerated into the collision cell.

Perform this procedure, one compound at a time, using the MS-specific optimized values that were obtained previously. The product ions provide the Q3 mass of MRM transitions.

In this example, the compound reserpine is used.

1. In the Tune Method Editor, close the XIC panes.
2. Click **Product Ion (MS2)** in the **Scan type** field.
3. Click the **Compound** tab and then type the optimal value noted previously.
4. On the MS tab, in the **Product Of** field, type **609.4**. This value is the mass assignment for reserpine that was noted in [Confirm the Presence of Compounds on page 165](#).
5. Make sure that the **Center / Width** check box is not selected.
6. In the mass table, type the appropriate Start, Stop, and Time values. Refer to [Table A-5](#).

Table A-5 Mass Table Parameters (Product Ion Scan)

Field	Value
Start (Da)	100
Stop (Da)	650
Time (sec)	2

7. Click **Edit Ramp**.
8. In the Ramp Parameter Settings dialog, select **Collision Energy** and then type the required **Start**, **Stop**, and **Step** values.

Note: The existing values are good starting points. Use the **Edit Ramp** function to change these values to be more efficient.

9. Click **OK**.
10. Select the **MCA** check box.
11. Click **Start**.
12. Right-click the bottom-right XIC pane and then click **Open File**.
13. Select the product ions with the greatest intensity and note the product ion m/z value to the first decimal point, such as 195.1.

We recommend that two or three product ions be optimized for each compound. The additional transitions can be used for confirmation, or to avoid the need to re-optimize a compound in case an interference is found.

Note: Make sure that the highest peaks selected for optimization do not represent a common loss from the precursor ion, such as water or carbon dioxide. Also make sure that the product ion is not too low in mass or interferences might occur for that transition in real samples or clusters from the mobile phase when analyzing on column.

14. Repeat this procedure for the remaining compounds.

Optimize Collision Cell Exit Potential for each Product Ion

1. In the Tune Method Editor, close the XIC panes.
2. Open the previously saved method.
3. In the mass table, verify the Q1 and Q3 m/z values for the compound.
4. Click the **Compound** tab and then type the optimal DP and CE values noted previously.
5. Click **Edit Ramp**.
6. Select **Collision Cell Exit Potential (CXP)** in the Ramp Parameter Settings dialog.
7. Type the required **Start**, **Stop**, and **Step** values.

Tip! The existing values are good starting points. Use the **Edit Ramp** function to change these values to be more efficient.

8. Click **OK**.
9. Click **Start**.
10. Right-click the bottom-right XIC pane and then click **Open File**.
11. Note the optimal value for the ion of interest.

The value that gives the best signal is the optimal value.

12. In the mass table, right-click and then select the parameter just optimized.

This adds a column to the table.

13. Repeat if other product ions were monitored.
14. Add the optimized values to the appropriate row.
15. Save the method.
16. Repeat this procedure for any other compounds that were previously optimized.

Manually Optimize the Ion Source and Gas Parameters

The ion source and gas settings must be set correctly to prevent contamination of the mass spectrometer and to make sure that the compounds of interest are optimally put into the gas phase as ions.

Note: The optimal ion source and gas settings are related to the solvent composition and the flow rate.

The ion source and gas settings must be adjusted when the LC conditions change significantly.

To optimize the ion source and gas parameters, set up a syringe pump with the compounds of interest and connect the line with a tee to the LC device. Controlling the LC pump can be done manually or through the software.

Another way to manually optimize ion source and gas settings is to use the autosampler to manually inject the compound of interest while manually varying the parameters in manual tuning to find the optimal settings.

Prepare the Ion Source

1. Set the horizontal micrometer to 5.
2. Set the vertical micrometer on the ion source for the flow rate.

Use the parameters in [Table A-6](#).

Table A-6 Turbo V™ Ion Source Vertical Parameters

Flow Rate	Initial Vertical Parameters
1 µL/min to 20 µL/min	10
20 µL/min to 250 µL/min	5
250 µL/min to 500 µL/min	2
500 + µL/min	0

Depending on how the optimization is run, a hardware profile with the LC pumps might need to be configured.

Refer to the ion source *Operator Guide*.

Optimize the Ion Source Parameters

Ion source parameters are optimized for best signal-to-noise for the compound of interest. The Curtain Gas™ supply is optimized at the highest setting without losing sensitivity. Refer to the ion source *Operator Guide*.

Use the following procedure to optimize the Curtain Gas™ flow parameter. The main function of the Curtain Gas™ flow parameter is to prevent the contamination of the ion optics. The Curtain Gas™ flow parameter should always be maintained as high as possible without losing sensitivity. The value depends on the type of mass spectrometer and ion source.

Do not set the parameter below the starting value.

The minimum flow rate supported by the APCI probe is 200 µL/min.

Note: Depending on the mass spectrometer, all parameters might not be available.

1. On the Navigation bar, under **Tune and Calibrate**, double-click **Manual Tuning**.
2. Click **File > Open**.
3. In the **Files** list, click the acquisition method used to optimize the compound parameter and then click **OK**.

The method opens in the Tune Method Editor.

4. Click the **Source/Gas** tab.
5. Using the ion source and gas flow guide, set all of the ion source and gas parameters so that they are appropriate for the flow rate.
6. Set the run time long enough so that many parameters can be adjusted. A recommended starting time is 15 minutes.
7. Click **Start**.

Data is shown in panes below the Tune Method Editor.

8. Note the signal of the peak of interest.
9. In the **Curtain Gas (CUR)** field, increase the value by five.
10. Continue increasing the **Curtain Gas (CUR)** value until the highest value without losing sensitivity is found.

As with most Source/Gas parameters, if two values give the same result, then use the higher value.

11. Repeat this procedure for the other Source/Gas parameters.

When optimizing for these parameters, look for the value that gives the highest signal-to-noise value.

Advanced Parameters

The following parameters should only be optimized by an experienced operator.

Optimize AF2

The AF2 parameter controls the fragmentation of the second precursor ion in an MS3 scan. The amount of excitation energy used depends on the compound and the desired amount of fragmentation.

1. For both positive and negative ion modes, ramp AF2 from 0.07 V to 0.3 V, with a step value of 0.01 V.

The allowable range is 0 V to 1 V.

2. For both positive and negative ion modes, ramp AF2 from 0 mV to 200 mV with a step value of 5 mV.
3. Select the AF2 value that gives the best amount of fragmentation.

If the MS3 scan is selected when using the 6500 and 6500⁺ series of instruments, the **No Fragmentation** option is not available. However, users can achieve the effect of the **No Fragmentation** option by setting the **AF2** parameter to **0** on the **Compound** tab of the Tune Method Editor.

About Collision Energy Spread (CES)

Under the compound-dependent parameters in Enhanced Product Ion (EPI) and MS3 experiments, the collision energy spread (CES) field allows users to specify the difference in collision energies that will be applied in the experiment. For example, if CE value of 30 and a CES value of 15 are used, then collision energies of 15, 30, and 45 are used.

Parameters for SCIEX Triple Quad 3500 Instruments

B

The following table contains generic parameters for the SCIEX Triple Quad 3500 instrument. The first number under each scan type is the preset value. The range of numbers is the accessible range for each parameter.

Table B-1 System Parameters for Triple Quadrupole Scan Types

Access ID	Positive Ion Mode			Negative Ion Mode		
	Q1	Q3	MS/MS	Q1	Q3	MS/MS
CUR ⁽¹⁾⁽²⁾	10 10 to 55	10 10 to 55	10 10 to 55	10 10 to 55	10 10 to 55	10 10 to 55
CAD	0 Fixed	5 Fixed	9 0 to 12	0 Fixed	5 Fixed	9 0 to 12
IS ⁽¹⁾⁽²⁾	5500 0 to 5500	5500 0 to 5500	5500 0 to 5500	-4500 -4500 to 0	-4500 -4500 to 0	-4500 -4500 to 0
NC ⁽³⁾	3 0 to 5	3 0 to 5	3 0 to 5	-3 -5 to 0	-3 -5 to 0	-3 -5 to 0
TEM ⁽²⁾⁽³⁾	0 0 to 750	0 0 to 750	0 0 to 750	0 0 to 750	0 0 to 750	0 0 to 750
DP	130 0 to 300	130 0 to 300	120 0 to 300	-60 -300 to 0	-60 -300 to 0	-150 -300 to 0
EP	10 2 to 15	10 2 to 15	10 2 to 15	-10 -15 to -2	-10 -15 to -2	-10 -15 to -2
CEM	2000 0 to 3300	2000 0 to 3300	2000 0 to 3300	2000 0 to 3300	2000 0 to 3300	2000 0 to 3300
GS1	15 0 to 90	15 0 to 90	15 0 to 90	15 0 to 90	15 0 to 90	15 0 to 90
GS2	0 0 to 90	0 0 to 90	0 0 to 90	0 0 to 90	0 0 to 90	0 0 to 90

Parameters for SCIEX Triple Quad 3500 Instruments

Table B-1 System Parameters for Triple Quadrupole Scan Types (continued)

Access ID	Positive Ion Mode			Negative Ion Mode		
IQ1 (IQ1 = Q0 + offset)	Q0 + (–0.5)	Q0 + (–0.5)	Q0 + (–0.5)	Q0 + 0.5	Q0 + 0.5	Q0 + 0.5
ST (ST = Q0 + offset)	Q0 + (–8) –12 to –5	Q0 + (–8) –12 to –5	Q0 + (–8) –12 to –5	Q0 + 8 12 to 5	Q0 + 8 12 to 5	Q0 + 8 12 to 5
IE1 (IE1 = Q0 – RO1)	0.9 0 to 3	n/a	0.9 0 to 3	–1 –3 to 0	n/a	–1 –5 to 0
CE (CE = Q0 – RO2)	n/a	n/a	53 5 to 180	n/a	n/a	–40 –180 to –5
CXP (CXP = RO2 – ST3)	n/a	9 0 to 55	27 0 to 55	n/a	–17 –55 to 0	–12 –55 to 0
IE3	n/a	2 0 to 5	1.5 0 to 5	n/a	–2.500 –5 to 0	–1.2 –5 to 0
(1) Turbo V™ ion source (2) TurbolonSpray® probe (3) APCI probe						

Calibration Ions and Solutions

C

Table C-1 Tuning Frequency

Calibration			Resolution Optimization	
Scan Type	Frequency	Manual/Automated	Frequency	Manual/Automated
Q1 and Q3	3 months to 6 months	Both	3 months to 6 months	Both

Table C-2 Suggested Tuning Solutions for the SCIEX Triple Quad™ 3500 System

System	Positive	Negative
SCIEX SCIEX Triple Quad™ 3500 LC-MS/MS system	POS PPG, 1e-5 M	NEG PPG, 3e-4 M

Table C-3 Q1 and Q3 PPG Scans

Polarity	Masses							
Positive	59.05	175.13	500.38	616.46	906.67	1254.93	1545.13	1952.43
Negative	45.00	411.26	585.39	933.64	1223.85	1572.10	1863.31	1979.39

Toolbar Icons

D

For additional toolbar icons, refer to the *Advanced User Guide*.

Table D-1 Tool Bar Icons











Icon	Name	Description
	New Subproject	Creates a subproject. Subprojects can only be created later in the process if the project was originally created with subprojects.
	Copy Subproject	Copies a Subproject folder. Subprojects can be copied only from another project that has existing subprojects. If the same folders exist at both the project and subproject levels, then the software uses the project level folders.

Table D-2 Acquisition Method Editor Icons

Icon	Name	Description
	Mass Spec	Click to show the MS tab in the Acquisition Method editor.
	Period	Right-click to add an experiment, add an IDA Criteria Level , or delete the period.
	Autosampler	Click to open the Autosampler Properties tab.
	Syringe Pump	Click to open the Syringe Pump Properties tab.
	Column Oven	Click to open the Column Oven Properties tab.
	Valve	Click to open the Valve Properties tab.
	DAD	Click to open the DAD Method Editor. Refer to Generate DAD Data on page 122 .
	ADC	Click to open the ADC Properties tab. Refer to Show ADC Data on page 113 .

Toolbar Icons

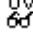Table D-3 Acquire Mode Icons

















Icon	Name	Description
	View Queue	Shows the sample queue.
	Instrument Queue	Shows a remote instrument station.
	Status for Remote Instrument	Shows the status of a remote instrument.
	Start Sample	Starts the sample in the queue.
	Stop Sample	Stops the sample in the queue.
	Abort Sample	Aborts a sample acquisition in the middle of the processing of that sample.
	Stop Queue	Stops the queue before it has completed processing all the samples.
	Pause Sample Now	Inserts a pause in the queue.
	Insert Pause before Selected Sample(s)	Inserts a pause before a specific sample.
	Continue Sample	Continues acquiring the sample.
	Next Period	Starts a new period.
	Extend Period	Extends the current period.
	Next Sample	Stops acquiring the current sample and starts acquiring the next sample.
	Equilibrate	Selects the method to be used to equilibrate the devices. This method should be the same as the method used with the first sample in the queue.
	Standby	Puts the instrument in Standby mode.
	Ready	Puts the instrument in Ready mode.

Table D-3 Acquire Mode Icons (continued)













Icon	Name	Description
	Reserve Instrument for Tuning	Reserves the mass spectrometer for tuning and calibrating.
	IDA Method Wizard	Starts the IDA Method Wizard.
	Purge Modifier	Starts the modifier purge from the modifier pump.

Table D-4 Tune and Calibrate Mode Icons

Icon	Name	Description
	Calibrate from spectrum	Opens the Mass Calibration Option dialog and uses the active spectrum to calibrate the mass spectrometer.
	Manual Tune	Opens the Manual Tune Editor.
	Compound Optimization	Optimizes for a compound using infusion by FIA.
	Instrument Optimization	Verifies the instrument performance, adjusts the mass calibration, or adjusts mass spectrometer settings.
	View Queue	Views the sample queue.
	Instrument Queue	Views a remote instrument.
	Status for Remote Instrument	Views the status of a remote instrument.
	Reserve Instrument for Tuning	Reserves the instrument for tuning and calibrating.
	IDA Method Wizard	Starts the IDA Method Wizard.

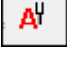Toolbar Icons

Table D-5 Explore Quick Reference: Chromatograms and Spectrum



















Icon	Name	Description
	Open Data File	Opens files.
	Show Next Sample	Goes to the next sample.
	Show Previous Sample	Goes to the previous sample.
	Go To Sample	Opens the Select Sample dialog.
	List Data	Views the data in tables.
	Show TIC	Generates a TIC from a spectrum.
	Extract Using Dialog	Extracts ions by selecting masses.
	Show Base Peak Chromatogram	Generates a BPC.
	Show Spectrum	Generates a spectrum from a TIC.
	Copy Graph to new Window	Copies the active graph to a new window.
	Baseline Subtract	Opens the Baseline Subtract dialog.
	Threshold	Adjusts the threshold.
	Noise Filter	Shows the Noise Filter Options dialog, which can be used define the minimum width of a peak. Signals below this minimum width are regarded as noise.
	Show ADC	Shows ADC data.
	Show File Info	Shows the experimental conditions used to collect the data.
	Add arrows	Adds arrows to the X-axis of the active graph.
	Remove all arrows	Removes arrows from the X-axis of the active graph.
	Offset Graph	Compensates for slight differences in the time during which the ADC data and the mass spectrometer data were recorded. This is useful when overlaying graphs for comparison.

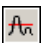Table D-5 Explore Quick Reference: Chromatograms and Spectrum (continued)

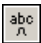













Icon	Name	Description
	Force Peak Labels	Labels all of the peaks.
	Expand Selection By	Sets the expansion factor for a portion of a graph to be viewed in greater detail.
	Clear ranges	Returns the expanded selection to normal view.
	Set Selection	Defines start and stop points for a selection. This feature provides more accurate selection than is possible by selecting the region using the cursor.
	Normalize To Max	Scales a graph to maximum size, so that the most intense peak is scaled to full scale, whether or not it is visible.
	Show History	Shows a summary of data processing operations performed on a particular file, such as smoothing, subtraction, calibration, and noise filtering.
	Open Compound Database	Opens the compound database.
	Set Threshold	Adjusts the threshold.
	Show Contour Plot	Shows selected data as either a spectrum graph or an XIC. Additionally, for data acquired by a DAD, a contour plot can show selected data as either a DAD spectrum or an XWC.
	Show DAD TWC	Generates a TWC of the DAD spectrum.
	Show DAD Spectrum	Generates a DAD spectrum.
	Extract Wavelength	Extracts up to three wavelength ranges from a DAD spectrum to view the XWC.

Table D-6 Explore Toolbar Quick Reference: Overlaying Graphs

Icon	Name	Description
	Home Graph	Click to return the graph to the original scale.
	Overlay	Click to overlay graphs.

Toolbar Icons

Table D-6 Explore Toolbar Quick Reference: Overlaying Graphs (continued)



Icon	Name	Description
	Cycle Overlays	Click to cycle between overlaid graphs.
	Sum Overlays	Click to add the graphs together.

Table D-7 Explore Toolbar Quick Reference: Fragment Interpretation Tool


Icon	Name	Description
	Show Fragment Interpretation Tool	Click to open Fragment Interpretation tool, which calculates the single, non-cyclic bond cleavage fragments from a .mol file.

Table D-8 Navigation Icons on the Explore Toolbar





























Icon	Name	Function
	Open File	Click to open files.
	Show Next Sample	Click to navigate to the next sample.
	Show Previous Sample	Click to navigate to the previous sample.
	GoTo Sample	Click to open the Select Sample dialog.
	List Data	Click to view the data in tables.
	Show TIC	Click to generate a TIC from a spectrum.
	Extract Using Dialog	Click to extract ions by selecting masses.
	Show Base Peak Chromatogram	Click to generate a BPC.
	Show Spectrum	Click to generate a spectrum from a TIC.
	Copy Graph to new Window	Click to copy the active graph to a new window.
	Baseline Subtract	Click to open the Baseline Subtract dialog.

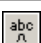Table D-8 Navigation Icons on the Explore Toolbar (continued)

Icon	Name	Function
	Threshold	Click to adjust the threshold.
	Noise Filter	Click to use the Noise Filter Options dialog to define the minimum width of a peak. Signals below this minimum width are regarded as noise.
	Show ADC	Click to view ADC data.
	Show File Info	Click to show the experimental conditions you used to collect your data.
	Add arrows	Click to add arrows to the x-axis of the active graph.
	Remove all arrows	Click to remove arrows from the x-axis of the active graph.
	Offset Graph	Click to compensate for slight differences in the time during which the ADC data and the mass spectrometer data were recorded. This is useful when overlaying graphs for comparison.
	Force Peak Labels	Click to label all the peaks.
	Expand Selection By	Click to set the expansion factor for a portion of a graph that you want to view in greater detail.
	Clear ranges	Click to return the expanded selection to normal view.
	Set Selection	Click to type start and stop points for a selection. This provides more accurate selection than is possible by highlighting the region using the cursor.
	Normalize to Max	Click to scale a graph to maximum, so that the most intense peak is scaled is to full scale, whether or not it is visible.
	Show History	Click to view a summary of data processing operations performed on a particular file, such as smoothing, subtraction, calibration, and noise filtering.
	Open Compound Database	Click to open the compound database.
	Set Threshold	Click to adjust the threshold.
	Show Contour Plot	Click to display selected data as either a spectrum graph or an XIC. Additionally, for data acquired by a DAD, a contour plot can display selected data as either a DAD spectrum or an XWC.
	Show DAD TWC	Click to generate a TWC of the DAD.

Toolbar Icons

Table D-8 Navigation Icons on the Explore Toolbar (continued)



Icon	Name	Function
	Show DAD Spectrum	Click to generate a DAD spectrum.
	Extract Wavelength	Click to extract up to three wavelength ranges from a DAD spectrum to view the XWC.

Table D-9 Integration Tab and Quantitation Wizard Icons








Icon	Name	Description
	Set parameters from Background Region	Uses the selected peak.
	Select Peak	Uses the selected background.
	Manual Integration Mode	Manually integrates peaks.
	Show or Hide Parameters	Toggles the peak-finding parameters between shown and hidden.
	Show Active Graph	Shows the analyte chromatogram only.
	Show Both Analyte and IS	Shows the analyte and its associated chromatogram (available only when an associated internal standard exists).
	Use Default View for Graph	Returns to the preset (view all data) view (if, for example, the user has zoomed in on a chromatogram).

Table D-10 Results Table Icons





Icon	Name	Description
	Sort Ascending by Selection	Sorts the selected column by ascending values.
	Sort Descending by selection	Sorts the selected column by descending values.
	Lock Or Unlock Column	Locks or unlocks the selected column. A locked column cannot be moved.
	Metric Plot By Selection	Creates a metric plot from the selected column.

Table D-10 Results Table Icons (continued)

















Icon	Name	Description
	Show all Samples	Shows all the samples in the Results Table.
	Delete Formula Column	Deletes formula columns.
	Report Generator	Opens the Reporter software.









Table D-11 Icon Quick Reference: Quantitate Mode










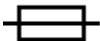
Icon	Name	Description
	Add/Remove Samples	Adds or removes samples from the Results Table.
	Export as Text	Saves the Results Table as a text file.
	Modify Method	Opens a wiff file.
	Peak Review - Pane	Opens peaks in a pane.
	Peak Review - Window	Opens peaks in a window.
	Calibration - Pane	Opens the calibration curve in a pane.
	Calibration - Window	Opens the calibration curve in a window.
	Show First Peak	Shows the first peak in the pane or window.
	Show Last Peak	Shows the last peak in the pane or window.
	Show Audit Trail	Shows the audit trail for the Results Table.
	Clear Audit Trail	Clears the audit trail for the Results Table. This functionality is not available.
	Statistics	Opens the Statistics window.
	Report Generator	Opens the Reporter software.

Glossary of Symbols













E



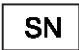









Note: Not all of the symbols in the following table are applicable to every instrument.

Symbol	Description
	Australian Regulatory Compliance Mark. Indicates the products complies with Australian communications Media Authority (ACMA) EMC Requirements.
	Alternating current
A	Amperes (current)
	Authorized representative in the European community
	Biohazard
	CE Marking of Conformity
	cCSAus mark. Indicates electrical safety certification for Canada and USA.
	Catalogue number
	Caution Note: In SCIEX documentation, this symbol identifies a personal injury hazard.



Symbol	Description
	China RoHS Caution Label. The electronic information product contains certain toxic or hazardous substances. The center number is the Environmentally Friendly Use Period (EFUP) date, and indicates the number of calendar years the product can be in operation. Upon the expiration of the EFUP, the product must be immediately recycled. The circling arrows indicate the product is recyclable. The date code on the label or product indicates the date of manufacture.
	China RoHS logo. The device does not contain toxic and hazardous substances or elements above the maximum concentration values, and it is an environmentally-friendly product that can be recycled and reused.
	Consult instructions for use.
	cTUVus mark for TUV Rheinland of North America.
	Data Matrix symbol that can be scanned by a barcode reader to obtain a unique device identifier (UDI).
	Ethernet connection
	Explosion Hazard
	Fire Hazard
	Fragile
	Fuse
Hz	Hertz

Glossary of Symbols

Symbol	Description
	High Voltage. Electrical Shock Hazard If the main cover must be removed, contact a SCIEX representative to prevent electric shock.
	Hot Surface Hazard
	In Vitro Diagnostic Device
	Ionizing Radiation Hazard
	Keep dry. Do not expose to rain. Relative humidity must not exceed 99%.
	Keep upright.
	Laser Radiation Hazard
	Lifting Hazard
	Manufacturer
	Moving Parts Hazard
	Pinch Hazard
	Pressurized Gas Hazard

Symbol	Description
	Protective Earth (ground)
	Puncture Hazard
	Serial number
	Toxic Chemical Hazard
	Transport and store the system within 66 kPa to 103 kPa.
	Transport and store the system within 75 kPa to 101 kPa.
	Transport and store the system within 10% to 90% relative humidity.
	Transport and store the system within -30 °C to +45 °C.
	Transport and store the system within -30 °C to +60 °C.
	USB 2.0 connection
	USB 3.0 connection
	Ultraviolet Radiation Hazard
VA	Volt Ampere (power)

Glossary of Symbols

Symbol	Description
V	Volts (voltage)
	WEEE. Do not dispose of equipment as unsorted municipal waste.
W	Watts
	<i>yyyy-mm-dd</i> Date of manufacture

Glossary of Warnings

F

Note: If any of the labels used to identify a component become detached, contact an FSE.

Label	Translation (if applicable)
EN61326—1, EN61326—2-6, CLASS A, GROUP 1, ISM EQUIPMENT	EN61326—1, EN61326—2-6, CLASS A, GROUP 1, ISM EQUIPMENT
FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.	FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.
IMPACT INDICATOR SENSITIVE PRODUCT WARNING	IMPACT INDICATOR SENSITIVE PRODUCT WARNING Note: If the indicator is tripped, then this container has been dropped or otherwise mishandled. Make a note on the Bill of Lading and then check for damage. Any claims for shock damage require a notation.
IMPORTANT! RECORD ANY VISIBLE CRATE DAMAGE INCLUDING TRIPPED "IMPACT INDICATOR" OR "TILT INDICATOR" ON THE WAYBILL BEFORE ACCEPTING SHIPMENT AND NOTIFY YOUR LOCAL AB SCIEX CUSTOMER SUPPORT ENGINEER IMMEDIATELY. DO NOT UNCRATE. CONTACT YOUR LOCAL CUSTOMER SUPPORT ENGINEER FOR UNCRATING AND INSTALLATION.	IMPORTANT! RECORD ANY VISIBLE CRATE DAMAGE INCLUDING TRIPPED "IMPACT INDICATOR" OR "TILT INDICATOR" ON THE WAYBILL BEFORE ACCEPTING SHIPMENT AND NOTIFY YOUR LOCAL AB SCIEX CUSTOMER SUPPORT ENGINEER IMMEDIATELY. DO NOT UNCRATE. CONTACT YOUR LOCAL CUSTOMER SUPPORT ENGINEER FOR UNCRATING AND INSTALLATION.
MINIMUM OF SIX PERSONS REQUIRED TO SAFELY LIFT THIS EQUIPMENT	MINIMUM OF SIX PERSONS REQUIRED TO SAFELY LIFT THIS EQUIPMENT
TIP & TELL	Tilt Indicator Note: Indicates whether the container was tipped or mishandled. Write on the Bill of Lading and inspect for damage. Any claims for tipping require a notation.

Glossary of Warnings

Label	Translation (if applicable)
TiltWatch PLUS ShockWatch	Tilt Indicator Note: Indicates whether the container was tipped or mishandled. Write on the Bill of Lading and inspect for damage. Any claims for tipping require a notation.
WARNING: DO NOT OPERATE WITHOUT FIRST ENSURING BOTTLE CAP IS SECURED.	WARNING: DO NOT OPERATE WITHOUT FIRST ENSURING BOTTLE CAP IS SECURED. Note: This warning is attached to the source exhaust waste bottle.
WARNING: NO USER SERVICEABLE PARTS INSIDE. REFER SERVICING TO QUALIFIED PERSONNEL.	WARNING: NO USER SERVICEABLE PARTS INSIDE. REFER SERVICING TO QUALIFIED PERSONNEL. Note: Consult instructions for use.